READ THIS MANUAL COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT

Keep this manual in a safe location for future reference

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Questions or comments?

Call our US-based customer service toll free at 1-800-536-0366.

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CONTENTS

1. Important safety notes	2
2. Intended use.....	4
3. Package contents.....	4
4. Device description.....	4
5. Set up	4
6. Use	4
7. Changing the infrared bulb	5
8. Care, maintenance and disposal.....	5
9. Troubleshooting Guide	5
10. Technical specifications	5
11. FCC Compliance Information.....	6
12. Warranty	7

Caution, hot surface

Do not stare at light source

Eye protection must be worn

Caution, infrared radiation

Protective insulation class 2

Protected against solid foreign objects 0.49 in (12.5 mm) in diameter and larger

DANGER

Electric shock:

- Do not place or store product where it can fall or be pulled into a bathtub, sink, water, or other liquid.
- Do not use this product while bathing, showering, washing dishes, or close to water sources of any kind.
- Do not place device in or drop into water or other liquid. Do not allow any water to be sprayed onto the unit. The unit must be operated only when it is completely dry.
- Open the unit only to change the light bulb. Do not attempt to repair the unit in case of malfunctions. If you need to have the unit repaired, please contact our customer service.
- Do not pinch or bend the power cord or pull it over sharp objects. This can cause an electric shock.
- Do not touch the plug with wet hands.
- Do not operate this device outdoors.
- Contact Customer Service or your dealer if the cord or housing are damaged.

Safety advice

- The power cord presents a strangulation hazard. NEVER LEAVE CHILDREN OR THOSE WHO REQUIRE CLOSE SUPERVISION UNATTENDED WITH THIS DEVICE.
- Do not touch the unit with wet hands if it is plugged in.
- Route the power cord so that it does not pose a trip hazard.
- Do not use on babies and small children.

Operation / proper use:

- Packaging materials are a deadly hazard for children and can cause suffocation. Remove all packaging materials immediately and keep them away from children at all times.

SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE

If the instruction manual is damaged or if you no longer have the instruction manual in your possession, please contact Beurer customer service. Please see warranty for service contact.

The manufacturer shall not be held liable for damage or injuries caused by improper or incorrect use.

1. IMPORTANT SAFETY NOTES

Signs and symbols

Read these instructions for use carefully. Observe the warnings and safety notes. Keep these instructions for use for future reference. Make the instructions for use accessible to other users. If the device is passed on, provide the instructions for use to the next user as well.

Whenever used, the following signs identify safety and property damage messages and designate a level of hazard or seriousness.

DANGER

Indicates a hazardous situation which, if not avoided, will result in death or serious injury.

WARNING

Indicates a hazardous situation which, if not avoided, could result in death or serious injury.

NOTICE

Addresses practices not related to personal injury, such as product and/or property damage.

Observe the instructions for use

- Remove all packaging materials before using the unit.
- Before use, ensure that there is no visible damage to the unit or accessories. Never operate the unit if it
 - has any damaged parts, e.g. housing,
 - has been exposed to extreme humidity or to any - liquids inside the housing,
 - has been dropped,
 - is not working properly.
- If the unit does not function correctly or if you notice any abnormalities, e.g. an unusual noise or smell, stop using the unit immediately. When in doubt, do not use the unit and contact Beurer customer service. (Please see warranty for service contact).
- This product should never be left unattended when in operation.
- Do not carry this appliance by supply cord or use cord as handle.
- The unit must be protected and stored away from sources of heat. (e.g. stoves, heating radiators)
- Connect the unit only to the voltage listed on the data plate.
- Please make sure the unit is not placed on flammable materials such as table cloths or carpets.
- Do not leave any easily flammable or meltable items close to the unit.
- The manufacturer is not liable for damage caused by improper or incorrect use.
- Using the unit is not a substitute for consultation with or treatment from a physician. Whenever you have any pain, illness or discomfort, always contact your physician first. If you have an emergency please call 911 immediately.
- If you start feeling any discomfort, pain, drowsiness or other change in your condition, stop using the unit immediately. For medical advice please contact a licensed physician. If you have an emergency please call 911 immediately.
- This device is not intended to be used by persons (including children) with restricted physical (e.g. activity limitations, disabilities), sensory (e.g. insensitivity to pain) or mental disabilities (e.g. mental disorders) or persons lacking the required experience or knowledge for safe operation of the device, unless closely supervised by a person responsible for their safety.

WARNING

Safety advice:

- Never use this product while in bed.
- Use the unit only on human beings and for its intended purpose and in the manner described in these instructions for use. Any improper use can be dangerous!
- Do not operate the unit in the presence of flammable gases, oxygen and aerosol spray products. Non-compliance may result in e.g. fire, burns, explosions, etc.
- Do not use this unit while operating a vehicle or a machine.
- The glass screen and unit housing become very hot during operation. You could be burned by touching it!
- Switch the unit off immediately if it becomes faulty or malfunctions.
- Do not pull on the power cord to remove the plug from the outlet.
- Do not look directly at the infrared light; close or cover the eyes. Infrared heat does not tan your skin!
- Certain products, medicines, cosmetics and food can lead to oversensitive or allergic skin reactions. Treatment should be discontinued immediately in these cases.
- Do not use the unit if there is the risk of you falling asleep!
- Excessive exposure can result in skin burns.
- Maintain a distance of at least 40 inches from other objects, both to protect the unit itself from overheating and to prevent the risk of fire.

Operation / proper use:

- This unit is not intended for commercial or clinical use, but only for private household use!
- This device is not intended for exposure to overhead use. Positioning it above the head can lead to safety hazards.
- Disconnection from the power supply is only guaranteed if the plug has been removed from the wall outlet.

NOTICE

- Use an easily accessible power supply so that you can disconnect the plug quickly in an emergency.
- Protect the unit or the parts from strong impacts such as dropping the unit on the floor.
- Protect the unit against dust, dirt and moisture. This may damage the device.
- Do not bump the unit while it is operating or cooling down or you could damage the lamp.

2. INTENDED USE

2.1 Intended purpose

This infrared radiant heater is only intended for radiating the human body (especially the muscles). Infrared light causes warmth to be transported to the human body. Blood supply to the radiated skin is improved and the metabolic rate is increased in the presence of temperature. The improved blood circulation stimulates and improves the healing processes in the radiated part of the body.

2.2 Indications for use

Infrared light can be used for muscle relaxation, relief of muscle spasm and temporary relief of minor pain. Infrared heat can also be used as a support tool in the temporary relief of pain and stiffness, also for stiffness associated with arthritis, because the heat temporarily stimulates the local blood circulation. You should first ask your GP whether its use is medically advisable in particular cases.

Contraindications

DO NOT USE INFRARED HEAT AFTER RECENT USE OF ACCUTANE OR IF YOU ARE PRESENTLY USING PHOTSENSITIZING DRUGS OR COSMETICS!

It is strongly recommended NOT to use the unit in the following cases in order to prevent health injuries:

- for persons insensitive to heat. Heat sensitivity can be reduced or increased in the following cases:
- in diabetic patients
- in people suffering from drowsiness, dementia or a lack of concentration
- in people with skin changes caused by illness
- in people with scarred skin in the treatment area
- in people with allergies
- in children and elderly people
- after taking medicine or consuming alcohol.

3. PACKAGE CONTENTS

- Infrared heat lamp
- Instruction manual
- Protective eyewear

4. DEVICE DESCRIPTION


- | | |
|---|-----------------|
| 1 | Housing |
| 2 | Protective grid |
| 3 | Device base |

5. SET UP

1. Place the unit on a solid and even surface. The unit must not be set up on a slope or laid down.

NOTE: The minimum distance between flammable items and the infrared lamp must not be less than 40 inches (1 m).

2. Unwind the power cord fully.
3. Insert the plug into a wall outlet.

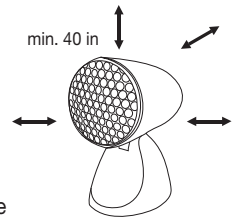6. USE

Preparing for treatment

1. Arrange yourself in front of the unit so you can relax and enjoy the treatment on the targeted body part.
2. Tilt the lamp so the unit is angled directly at the part of the body being treated.
3. The distance between the body part and the infrared lamp is generally determined based on individual sensitivity to heat and the type of treatment being performed.
4. Remove any metal items, e.g. belt buckles, necklaces, bras, jewellery or piercings, from the treatment area before switching on the device as they can become very hot.

Treatment duration

We recommend beginning with short treatment times. However, even when the unit is used correctly, the skin of some people can have an oversensitive or allergic reaction to the heat treatment (e.g. intense redness, blisters, itching, heavy sweating).

NOTE: Should you notice any discomfort or unusual reactions, stop using the device.

Starting the unit

Important

Always use the supplied protective eyewear during treatment.

The eye-protection cover may have an initial scent, which will dissipate after a few uses.

1. Insert the mains plug into a suitable socket. The infrared lamp begins to light up.
2. You are now ready to start treatment.
3. To end treatment, pull the mains plug out of the socket.
4. Remove the protective eyewear.

NOTE: Do not use the unit for more than 15

minutes at a time (risk of overheating). Allow the unit to cool down before you use it again.

After treatment

1. The infrared lamp becomes hot after use. Allow it to cool down sufficiently before touching.
2. Please note that the mains cable must not be wound around the infrared lamp, either during use or storage.
3. Unplug the device from the socket after use to ensure it is completely switched off.


NOTE: Do not cover, pack, or store the device while it is still warm.

7. CHANGING THE INFRARED BULB


Procedure

NOTE: Remember that a faulty bulb should only be replaced with a bulb of the same type.


1. Before changing the bulb, switch the device off, unplug it and allow it to cool down completely.
2. Using a suitable Phillips screwdriver, unscrew the screw on the left side of the device.

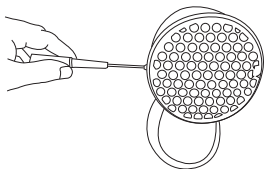3. Remove the protective grid by unscrewing the side screw with a screwdriver and carefully pulling the grid off the device.

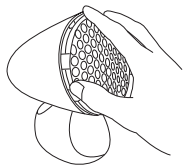4. Using a cloth or gloves, unscrew the old infrared lamp from the device..

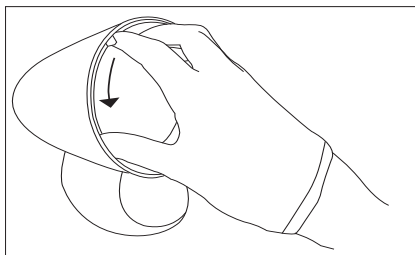NOTE: Do not use bare fingers, as the grease on your skin could damage the coating on the bulb.

5. Screw the new infrared lamp into the device.

6. The unit is closed by following the reverse sequence of steps used for opening.

NOTE: Make sure that no foreign objects are left in the body of the unit and check that the protective grid is firmly attached to the device.

Replacement bulb

Type: Halogen bulb 110V – 120V / 100 W E27 socket
Original Order no. 110.248

The service life of the bulb should be approximately 300 hours

8. CARE, MAINTENANCE AND DISPOSAL

Care and Maintenance

- Before cleaning the unit, unplug it and allow it to cool down completely before you touch it or pack it away.
- Never immerse the device in water or other liquids.
- The unit can be cleaned with a damp cloth. Do not use cleaning agents.
- Ensure that no water enters the unit.

Disposal

Dispose of the device in accordance with local regulations.

If you have any questions, please contact the local authorities responsible for waste disposal.

Unlawful disposal may cause environmental pollution.

9. TROUBLESHOOTING GUIDE

Problem	Possible cause	Solution
Infrared lamp does not illuminate	No power	Plug the power plug into the outlet correctly.
	Infrared lamp burn time exceeded	The infrared lamp is faulty. Contact Customer Service or your retailer
	After replacing the infrared bulb, the new infrared bulb was not fitted correctly	Check that the new infrared bulb has been fitted correctly


10. TECHNICAL SPECIFICATIONS

Voltage	AC 110V – 120V / 60 Hz
Power consumption	100 W
Risk group	RG2, Risk Group 2 (moderate risk)

Max. irradiation intensity according to EN 62471	< 3200 W/m ²
Dimensions (W x H x D)	7.1 x 4.8 x 8.1 in 180 x 122 x 205 mm
Weight	Approx. 0,77 lbs (349 g)
Operating conditions	Temperature: 41 to 86°F (+5°C to +30°C) Relative humidity: 15 to 93%
Storage conditions	Temperature: -13 to 158°F (-25°C to +70°C) Relative humidity: ≤ 93%
Emitted infrared wavelength range	IRA + IRB / 590 - 1700 nm
Protection class	II
Expected burn time of light source	300 hours
IP protection class	IP20
Protective eyewear	Corresponds to the requirements of ANSI Z87.1 and serves to protect the eyes.
Safety distance to the eyes	5.75 in (14.6 cm)

The serial number is located on the device.

Spectral irradiation intensity

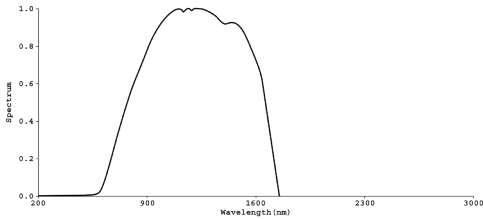We retain the right to make technical modifications to improve and further develop the product. Due to manufacturing processes, the color of the glass may change during operation but this does not affect the performance or quality of the unit.

NOTE: If the device is used outside of the specifications, proper function is no longer guaranteed!

This equipment complies with Standard IEC 60601-2-83 and is subject to specific safety provisions for electromagnetic compatibility. Please note, however, that portable and mobile HF communication equipment may influence this unit. Further details may be obtained from our Customer Service Department.

ELECTROMAGNETIC COMPATIBILITY

- The device complies with current specifications with regard to electromagnetic compatibility and is suitable for use in all premises, including those designated for private residential purposes. The radio frequency emissions of the device are extremely low and in all probability do not cause any interference with other devices in the proximity.
- It is recommended that you do not place the device on top of or close to other electronic devices. Should you notice any interference with other electrical devices, move the device or connect it to a different socket.
- Radio equipment and TV interference may affect the operation of this device.

11. FCC COMPLIANCE INFORMATION

Infrared heat lamp IL 11
Responsible Party – U.S. Contact Information:
Beurer North America LP
1 Oakwood Blvd, Suite 255
Hollywood, FL 33020
United States
1-800-536-0366
info@beurer.com

FCC Compliance Statement

This device complies with Part 18 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Consumer ISM equipment, pursuant to Part 18 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.

- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

If you're unsure whether your seller is authorized, please contact our support team.

Exclusions

This warranty does not cover damage caused by misuse or abuse; accident; the attachment of unauthorized accessories; alterations; improper installation; misapplication; lack of reasonable care; unauthorized repairs or modifications; improper electrical/power supply; worn batteries; normal wear and tear; power loss; dropped product; failure to follow instructions; lack of recommended maintenance; transit damage; theft; neglect; vandalism; or environmental conditions.

This warranty is void if the product is used in a commercial or business environment.

Beurer's maximum liability under this warranty is limited to the actual purchase price paid, as confirmed by proof of purchase. The warranty is effective only in the country where the product was purchased and used.

Legal Notice

This warranty gives you specific legal rights, and you may have other rights which vary from state to state. THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. ANY IMPLIED WARRANTIES, INCLUDING BUT NOT LIMITED TO MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO THE LENGTH OF THIS WRITTEN WARRANTY.

IN NO EVENT SHALL BEURER BE LIABLE FOR ANY INCIDENTAL OR CONSEQUENTIAL DAMAGES ARISING FROM USE OR INABILITY TO USE THE PRODUCT.

Beurer does not authorize any third party — including retailers or resellers — to modify the terms of this warranty.

This warranty does not apply to opened, used, repackaged, or resold products, including those sold through online auctions or unauthorized channels.

For Canada only:

Manufacturer:

**Beurer GmbH, Söflinger Str. 218
89077 Ulm, Germany**

Distributed by:

**Beurer North America LP
1 Oakwood Blvd, Suite 255
Hollywood, FL 33020, USA
www.shop-beurer.com**

Questions or comments?

**Call our US-based customer service toll free
at 1-800-536-0366 or contact info@beurer.com.**

Assembled in Hungary

Subject to error and change

12. WARRANTY

Limited Lifetime Warranty For Original Purchaser

Your Beurer Infrared heat lamp, Model IL 11, is warranted to be free from defects in materials and workmanship for the life of the product under normal conditions of intended use and service. This warranty extends only to the original retail purchaser and does not extend to retailers or subsequent owners.

We will, at our option, repair or replace the Beurer Infrared heat lamp, Model IL 11, without additional charge, for any part or parts covered by this written warranty. No refunds will be given. Repair or replacement is our only responsibility and your only remedy under this written warranty. If replacement parts for defective materials are not available, Beurer reserves the right to make product substitutions in lieu of repair or replacement.

For warranty service, contact our customer service department at 1-800-536-0366 or info@beurer.com with a description of the problem. If the issue is deemed to fall within the scope of the limited lifetime warranty, you will be asked to mail the product (at your cost) in its original package along with proof of purchase, your name, address, and phone number. If the issue is not deemed eligible under warranty, we will provide a quotation for repair or replacement, including return shipping fees.

Third-Party Repairs

This warranty complies with the federal Magnuson-Moss Warranty Act. You may use third-party parts or services for repairs without voiding this warranty, unless the repair itself causes damage to the product. We encourage you to retain your proof of purchase and contact our customer service team with any questions.

Unauthorized Sellers

This warranty applies only to products purchased from Beurer or authorized retailers and distributors. Products purchased from unauthorized sellers — including third-party platforms or online marketplaces (such as eBay or unverified Amazon sellers) — are not covered by this warranty. This policy is in place to protect customers from counterfeit, altered, or used products that may not meet Beurer's quality standards.

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