READ THIS MANUAL COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT

Keep this manual in a safe location for future reference

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Hollywood, FL 33020, USA

For Canada only/Solo para Canada/

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Manufacturer/ Fabricante/ Fabricant:

Beurer GmbH, Soflinger Str. 218

89077 Ulm Germany

www.shop-beurer.com

Assembled in China

Questions or comments?

Call our US-based customer service toll free at 1-800-536-0366.

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CONTENTS

1. Important safety notes	2
2. Intended use.....	4
3. Package contents.....	4
4. Device description.....	4
5. Set up	4
6. Use	4
7. Changing the infrared bulb.....	5
8. care, Maintenance and disposal	6
9. Troubleshooting Guide	6
10. Technical specifications	6
11. FCC Compliance Information.....	7
12. Warranty	7

SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE

If the instruction manual is damaged or if you no longer have the instruction manual in your possession, please contact Beurer customer service. Please see warranty for service contact.

The manufacturer shall not be held liable for damage or injuries caused by improper or incorrect use.

1. IMPORTANT SAFETY NOTES

Signs and symbols

Whenever used, the following signs identify safety and property damage messages and designate a level of hazard or seriousness.

READ AND UNDERSTAND THIS ENTIRE MANUAL, INCLUDING THE SAFETY SECTION AND ALL INSTRUCTIONS AND WARNINGS, COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT. FOLLOW ALL SAFETY INSTRUCTIONS AND WARNINGS TO AVOID HAZARDOUS SITUATIONS AND TO MAKE CORRECT USE OF THIS PRODUCT.

	DANGER	Indicates a hazardous situation which, if not avoided, will result in death or serious injury.
	WARNING	Indicates a hazardous situation which, if not avoided, could result in death or serious injury.
NOTICE		Addresses practices not related to personal injury, such as product and/or property damage.

	Distance
	Observe the instructions Read the instructions before starting work and/or operating devices or machines
	Caution, hot surface
	Do not stare at light source
	Eye protection must be worn
	Caution, infrared radiation
	Protective insulation class 2
	Protected against solid foreign objects 0.49 in (12.5 mm) in diameter and larger

DANGER

Electric shock:

- Do not place or store product where it can fall or be pulled into a bathtub, sink, water, or other liquid.
- Do not use this product while bathing, showering, washing dishes, or close to water sources of any kind.
- Do not place device in or drop into water or other liquid. Do not allow any water to be sprayed onto the unit. The unit must be operated only when it is completely dry.
- Open the unit only to change the light bulb. Do not attempt to repair the unit in case of malfunctions. If you need to have the unit repaired, please contact our customer service.
- Do not pinch or bend the power cord or pull it over sharp objects. This can cause an electric shock.
- Do not touch the plug with wet hands.
- Do not operate this device outdoors.
- Contact Customer Service or your dealer if the cord or housing are damaged.

Safety advice

- The power cord presents a strangulation hazard. **NEVER LEAVE CHILDREN OR THOSE WHO REQUIRE CLOSE SUPERVISION UNATTENDED WITH THIS DEVICE.**
- Do not touch the unit with wet hands if it is plugged in.

- Route the power cord so that it does not pose a trip hazard.
- Do not use on babies and small children.

Safety advice:

- Never use this product while in bed.
- Use the unit only on human beings and for its intended purpose and in the manner described in these instructions for use. Any improper use can be dangerous!
- Do not operate the unit in the presence of flammable gases, oxygen and aerosol spray products. Non-compliance may result in e.g. fire, burns, explosions, etc.
- Do not use this unit while operating a vehicle or a machine.
- The glass screen and unit housing become very hot during operation. You could be burned by touching it!
- Switch the unit off immediately if it becomes faulty or malfunctions.
- Do not pull on the power cord to remove the plug from the outlet.
- Do not look directly at the infrared light; close or cover the eyes. Infrared heat does not tan your skin!
- Certain products, medicines, cosmetics and food can lead to oversensitive or allergic skin reactions. Treatment should be discontinued immediately in these cases.
- Do not use the unit if there is the risk of you falling asleep!
- Excessive exposure can result in skin burns.
- Maintain a distance of at least 40 inches from other objects, both to protect the unit itself from overheating and to prevent the risk of fire.

Operation / proper use:

- Packaging materials are a deadly hazard for children and can cause suffocation. Remove all packaging materials immediately and keep them away from children at all times.
 - Remove all packaging materials before using the unit.
 - Before use, ensure that there is no visible damage to the unit or accessories. Never operate the unit if it
 - has any damaged parts, e.g. housing,
 - has been exposed to extreme humidity or to any - liquids inside the housing,
 - has been dropped,
 - is not working properly.
 - If the unit does not function correctly or if you notice any abnormalities, e.g. an unusual noise or smell, stop using the unit immediately. When in doubt, do not use the unit and contact Beurer customer service. (Please see warranty for service contact).
 - This product should never be left unattended when in operation.
 - Do not carry this appliance by supply cord or use cord as handle.
 - The unit must be protected and stored away from sources of heat. (e.g. stoves, heating radiators)
 - Connect the unit only to the voltage listed on the data plate.
 - Please make sure the unit is not placed on flammable materials such as table cloths or carpets.
 - Do not leave any easily flammable or meltable items close to the unit.
 - The manufacturer is not liable for damage caused by improper or incorrect use.
 - Using the unit is not a substitute for consultation with or treatment from a physician. Whenever you have any pain, illness or discomfort, always contact your physician first. If you have an emergency please call 911 immediately.
 - If you start feeling any discomfort, pain, drowsiness or other change in your condition, stop using the unit immediately. For medical advice please contact a licensed physician. If you have an emergency please call 911 immediately.
 - This device is not intended to be used by persons (including children) with restricted physical (e.g. activity limitations, disabilities), sensory (e.g. insensitivity to pain) or mental disabilities (e.g. mental disorders) or persons lacking the required experience or knowledge for safe operation of the device, unless closely supervised by a person responsible for their safety.
- This unit is not intended for commercial or clinical use, but only for private household use!
 - This device is not intended for exposure to overhead use. Positioning it above the head can lead to safety hazards.
 - Disconnection from the power supply is only guaranteed if the plug has been removed from the wall outlet.

NOTICE

- Use an easily accessible power supply so that you can disconnect the plug quickly in an emergency.
- Protect the unit or the parts from strong impacts such as dropping the unit on the floor.
- Protect the unit against dust, dirt and moisture. This may damage the device.
- Do not bump the unit while it is operating or cooling down or you could damage the lamp.

2. INTENDED USE

2.1 Intended purpose

This infrared radiant heater is only intended for radiating the human body (especially the muscles). Infrared light causes warmth to be transported to the human body. Blood supply to the radiated skin is improved and the metabolic rate is increased in the presence of temperature. The improved blood circulation stimulates and improves the healing processes in the radiated part of the body.

2.2 Indications for use

Infrared light can be used for muscle relaxation, relief of muscle spasm and temporary relief of minor pain. Infrared heat can also be used as a support tool in the temporary relief of pain and stiffness, also for stiffness associated with arthritis, because the heat temporarily stimulates the local blood circulation. You should first ask your GP whether its use is medically advisable in particular cases.

Contraindications

DO NOT USE INFRARED HEAT AFTER RECENT USE OF ACCUTANE OR IF YOU ARE PRESENTLY USING PHOTSENSITIZING DRUGS OR COSMETICS!

It is strongly recommended NOT to use the unit in the following cases in order to prevent health injuries:

- for persons insensitive to heat. Heat sensitivity can be reduced or increased in the following cases:
- in diabetic patients
- in people suffering from drowsiness, dementia or a lack of concentration
- in people with skin changes caused by illness
- in people with scarred skin in the treatment area
- in people with allergies
- in children and elderly people
- after taking medicine or consuming alcohol.

3. PACKAGE CONTENTS

- Infrared heat lamp
- Instruction manual
- Protective eyewear

4. DEVICE DESCRIPTION


- 1 Housing
- 2 Glass pane
- 3 Time selector button
- 4 Remaining time display
- 5 ON/OFF button
- 6 Device base

5. SET UP

1. Place the unit on a solid and even surface. The unit must not be set up on a slope or laid down.

NOTE: The minimum distance between flammable items and the infrared lamp must not be less than 40 inches (1 m).

2. Unwind the power cord fully.
3. Insert the plug into a wall outlet.

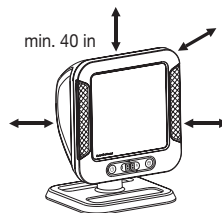6. USE

Preparing for treatment

1. Arrange yourself in front of the unit so you can relax and enjoy the treatment on the targeted body part.
2. Tilt the lamp so the unit is angled directly at the part of the body being treated.
3. Make sure the distance between the infrared heat lamp and the part of the body being treated is not less than 12 inches (30 cm) and not more than 16 inches (40 cm).
NOTE: The distance to the infrared heat lamp depends largely on personal heat sensitivity and the treatment concerned.
4. Remove any metal items, e.g. belt buckles, necklaces, bras, jewellery or piercings, from the treatment area before switching on the device as they can become very hot.

Treatment duration

We recommend beginning with short treatment times. However, even when the unit is used correctly, the skin of some people can have an oversensitive or allergic reaction to the heat treatment (e.g. intense redness, blisters, itching, heavy sweating).

NOTE: Should you notice any discomfort or unusual reactions, stop using the device.

Starting the unit

Important

Always use the supplied protective eyewear during treatment.

The eye-protection cover may have an initial scent, which will dissipate after a few uses.

1. Press the On/Off button. The unit will turn on automatically after a short time.
2. The remaining time display will show 15 minutes (recommended treatment time).

3. You can press the time selector button to reduce the treatment time. Each push reduces the time by 1 minute. To reset the timer to 15 minutes, turn the device off and then on again.

NOTE: Examine the part of the body treated frequently and check for skin reactions. Cease exposure immediately if there are any signs of oversensitivity or an allergic reaction, and consult a doctor.

4. The infrared heat lamp switches off automatically after the set time.

NOTE: Do not use the unit for more than 15 minutes at a time (risk of overheating). Allow the unit to cool down before you use it again.

5. To end the treatment before the timer runs out, press the On/Off button.

After treatment

1. The fan will continue running for 5 minutes after you have finished using the device in order to cool it down. Please do not unplug the device during this period.

NOTE: Do not move, cover up or pack the unit away while it is still warm.


2. Once the cooling period has finished, the fan will switch itself off. Remove the plug from the socket.
3. If possible, store the infrared heat lamp in its original packaging in a dry place.

7. CHANGING THE INFRARED BULB


Procedure

NOTE: Remember that a faulty bulb should only be replaced with a bulb of the same type.


1. Before changing the bulb, switch the device off, unplug it and allow it to cool down completely.
2. Place the cool unit with the glass ceramic plate down on a flat surface.
3. Using a suitable Phillips screwdriver, unscrew the 4 screws in the back panel.

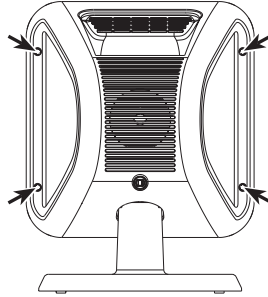4. Leave the glass ceramic plate on the surface and set the housing back on its base.

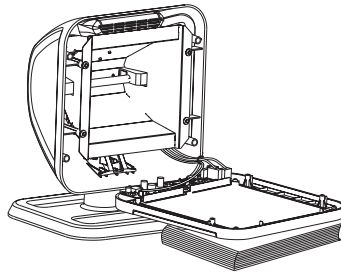5. Using a cloth or gloves, push the bulb to one side, then slide the opposite end forward out of the socket.

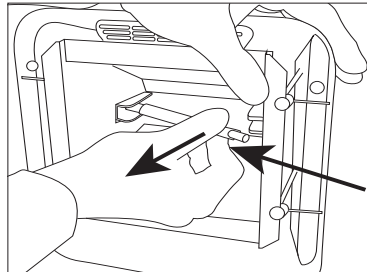NOTE: Do not use bare fingers, as the grease on your skin could damage the coating on the bulb.

6. To insert the bulb, place one end of the bulb onto either pin inside the housing. Push the bulb to the inserted side and slide the opposite end of the bulb onto the other pin. Ensure it is securely fitted.

7. The unit is closed by following the reverse sequence of steps used for opening.

NOTE: Make sure that no foreign objects are left in the body of the unit and that no cables become trapped.

Replacement bulb

Type: Halogen bulb 110V – 120V / 300 W

118 mm long for R7S socket

Original Order no. 110.150

The service life of the bulb should be approximately 1500 hours

8. CARE, MAINTENANCE AND DISPOSAL

Care and Maintenance

- Before cleaning the unit, switch it off, unplug it and allow it to cool down completely before you touch it or pack it away.
- Never immerse the device in water or other liquids.
- The unit can be cleaned with a damp cloth. Do not use cleaning agents.
- Ensure that no water enters the unit.

Disposal

Dispose of the device in accordance with local regulations.

If you have any questions, please contact the local authorities responsible for waste disposal.

Unlawful disposal may cause environmental pollution.

9. TROUBLESHOOTING GUIDE


Problem	Possible cause	Solution
Infrared lamp does not illuminate	ON/OFF button was not pressed	Press the ON/OFF button
	No power	Plug the power plug into the outlet correctly.
	Infrared lamp burn time exceeded	The infrared lamp is faulty. Contact Customer Service or your retailer
	After replacing the infrared bulb, the new infrared bulb was not fitted correctly	Check that the new infrared bulb has been fitted correctly
Infrared lamp shuts off during operation although timer has not yet elapsed	Unit automatically shut down to protect from overheating.	Allow the device to cool down and restart it

10. TECHNICAL SPECIFICATIONS

Voltage	AC 110V – 120V / 60 Hz
Power consumption	300 W
Risk group	RG2, Risk Group 2 (moderate risk)
Timer	Adjustable from 1 to 15 minutes in one-minute steps.
Max. irradiation intensity according to EN 62471	< 3200 W/m ²
Dimensions (W x H x D)	10.6 x 11.2 x 7.7 in 225 x 288 x 195 mm
Weight	Approx. 2.9 lbs (1.3 kg)
Operating conditions	Temperature: 41 – 95°F (+5°C to +35°C) Relative humidity: 15-93%
Storage and transport	Temperature: -13 – 158°F (-25°C to +70°C) Relative humidity: ≤ 93%
Emitted infrared wavelength range	IRA + IRB / 780 – 3000 nm
Protection class	II
Expected burn time of light source	1500 hours
Fuse	Type H 4AT/250V
Protective eyewear	Corresponds to the requirements of ANSI Z87.1 and serves to protect the eyes.

The serial number is located on the device.

Spectral irradiation intensity

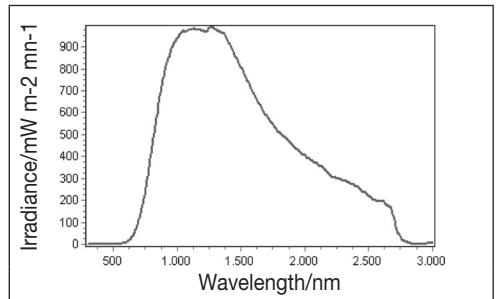We retain the right to make technical modifications to improve and further develop the product. Due to manufacturing processes, the color of the glass may

change during operation but this does not affect the performance or quality of the unit.

NOTE: If the device is used outside of the specifications, proper function is no longer guaranteed!

This equipment complies with Standard IEC 60601-2-83 and is subject to specific safety provisions for electromagnetic compatibility. Please note, however, that portable and mobile HF communication equipment may influence this unit. Further details may be obtained from our Customer Service Department.

ELECTROMAGNETIC COMPATIBILITY

- The device complies with current specifications with regard to electromagnetic compatibility and is suitable for use in all premises, including those designated for private residential purposes. The radio frequency emissions of the device are extremely low and in all probability do not cause any interference with other devices in the proximity.
- It is recommended that you do not place the device on top of or close to other electronic devices. Should you notice any interference with other electrical devices, move the device or connect it to a different socket.
- Radio equipment and TV interference may affect the operation of this device.

11. FCC COMPLIANCE INFORMATION

Infrared heat lamp IL 60

Responsible Party – U.S. Contact Information:

Beurer North America LP
1 Oakwood Blvd, Suite 255
Hollywood, FL 33020
United States
1-800-536-0366
info@beurer.com

FCC Compliance Statement

This device complies with Part 18 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Consumer ISM equipment, pursuant to Part 18 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can

radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

12. WARRANTY

Your Beurer Infrared heat lamp, Model IL 60, excluding all enclosed accessories as these require routine maintenance, is warranted to be free from defects in materials and workmanship for the life of the product under normal conditions of intended use and service. This warranty extends only to the original retail purchaser and does not extend to retailers or subsequent owners.

We will, at our option, repair or replace the Beurer Product, Model IL 60, without charge, for any part or parts covered by these written warranties. No refunds will be given. Repair or replacement is our only responsibility and your only remedy under this written warranty. If replacement parts for defective materials are not available, Beurer reserves the right to make product substitutions in lieu of repair or replacement.

For warranty service contact our customer service department at 1-800-536-0366 or at info@beurer.com to provide a description of the problem. If the problem is deemed to be within the scope of the limited lifetime warranty, you will be asked to mail the product at your costs in its original package with proof of purchase, your name, address and phone number. If the problem is not deemed to be within the scope of the limited lifetime warranty, we will provide quotation for repair, respectively, replacement and return shipping fee.

This warranty does not cover damage caused by misuse or abuse; accident; failures of the Infrared heat lamp, tubing, or accessories; parts or cleaning solutions; the attachment of unauthorized accessories; alteration to the product; improper installation; misapplication; lack of reasonable care with respect to the product; unauthorized repairs or modifications; improper use of electrical/power supply; normal wear; loss of power; dropped product; malfunction

or damage of an operating part as a result of failure to comply with instructions for use or to provide manufacturer's recommended maintenance; transit damage; theft; neglect; vandalism; or environmental conditions; loss of use during the period the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions whatsoever that are beyond the control of Beurer. This warranty is void if the product is ever used in a commercial or business environment. The maximum liability of Beurer under this warranty is limited to the purchase price actually paid by the customer for the product covered by the warranty, as confirmed by proof of purchase, regardless of the amount of any other direct or indirect damage suffered by the customer.

This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. A product that requires modifications or adaptation to enable it to operate in any other country than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under this warranty.

THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO THE DURATION OF THIS APPLICABLE WRITTEN WARRANTY. Some states do not allow limitations on how long an implied warranty lasts, so the above limitations may not apply to you.

IN NO EVENT SHALL BEURER BE LIABLE FOR ANY SPECIAL, INCIDENTAL, INDIRECT OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS, IMPLIED OR ANY OTHER THEORY OF LIABILITY, WHATSOEVER. Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.

Beurer does not authorize anyone, including, but not limited to, retailers, the subsequent consumer purchaser of the product from a retailer or remote purchaser, to obligate Beurer in any way beyond the terms set forth herein.

This warranty does not extend to the purchase of opened, used, repaired, repackaged and/or resealed products, including but not limited to sale of such products on Internet auction sites and/or products by surplus or bulk resellers. Any and all warranties or guarantees shall immediately cease and terminate in connection with any products or parts thereof which

are repaired, replaced, altered, or modified, without the prior explicitly written consent of Beurer.

This warranty gives you specific legal rights, and you may also have other rights which may vary from state to state.

For more information regarding our product line in the USA, please visit: www.shop-beurer.com

Distributed by:
Beurer North America LP
1 Oakwood Blvd, Suite 255
Hollywood, FL 33020
www.shop-beurer.com

For Canada only:
Manufacturer:
Beurer GmbH
Söflinger Str. 218
89077 Ulm, Germany

Questions or comments?
Call our US-based customer service toll free
at 1-800-536-0366 or contact info@beurer.com.

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