**READ THIS MANUAL COMPLETELY AND CAREFULLY BEFORE
USING THIS PRODUCT**

Keep this manual in a safe location for future reference

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 **Questions or comments? Call our US-based customer service toll free at 1-800-536-0366.**

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CONTENT

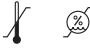
1. IMPORTANT SAFETY INSTRUCTIONS	4
DO NOT DESTROY	2
2. Package Contents	4
3. Parts and Controls.....	4
4. Intended Use	4
5. To Set Up	4
6. To Use	4
7. Care, Maintenance and Disposal.....	5
8. Technical Specifications	6
9. Troubleshooting Guide	6
10. FCC Compliance Information	6
11. Warranty	7

1. IMPORTANT SAFETY INSTRUCTIONS DO NOT DESTROY

READ THIS ENTIRE MANUAL, THE SAFETY SECTION AND ALL INSTRUCTIONS AND WARNINGS COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT. FOLLOW ALL SAFETY INSTRUCTIONS AND WARNINGS TO AVOID HAZARDOUS SITUATIONS AND TO MAKE CORRECT USE OF THIS PRODUCT.

Signs and symbols

Whenever used, the following signs identify safety and property damage messages and designate a level of hazard or seriousness.

	This is the safety alert symbol. It alerts you to potential personal injury hazards. Obey all safety messages that follow this symbol to avoid possible injury or death.
 WARNING	Indicates a hazardous situation which, if not avoided, could result in death or serious injury.
 NOTICE	Addresses practices not related to personal injury, such as product and/or property damage.
Storage 	Permissible storage and transport temperature and humidity
Operating 	Permissible operating temperature and humidity
	Observe the instructions for use
	Do not use the device with a pacemaker

WARNING

To reduce the risk of burns, electric shock and fire, this product must be used in accordance with the following instructions:

- Read all instructions before using the product.
- Do not use on infants.
- Never leave heating pad unattended.
- Do not use pins or other metallic means to fasten this pad in place.
- Never pull this pad by the power cord. Do not use the cord as a handle. Unplug when not in use.
- Do not use a heating pad when symptoms of appendicitis are present. Consult your physician.
- Inspect this product before using it in order to make sure that there are no broken parts, missing parts and/or damaged areas. Do not use heating pad if there are broken parts, missing parts or damaged areas.
- This appliance has a polarized plug (one blade is wider than the other). As a safety feature, this plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not attempt to defeat this safety feature.
- Do not use an electrical outlet that has become loose or does not engage the power plug tightly.
- Consult with a physician for use with controlled substances.
- Do not tuck in the wired area of this heating pad.
- Do not route electrical cords between the mattress and boxspring. Do not allow the cords to be pinched.
- Keep the controller away from drafty areas that may become damp or wet such as an open window.
- Do not use this heating pad if you or the product is wet.
- Be sure that this heating pad is used on an alternating current supply circuit of the proper voltage.
- Turn off and unplug when not in use.
- Indoor use only.
- The electronic and magnetic fields emitted by this heating pad may interfere with the function of a pacemaker. However, they are still well below the limits for electrical field strength (5000 V/m max), magnetic field strength (80 A/m max), and magnetic flux density (0.1 millitesla max). Please consult your doctor and the manufacturer of your pacemaker before using this heating pad.
- This model contains filling materials. If filling material is released due to product damage, discontinue use of product. Do not ingest filling material. If filling material is ingested, seek medical attention.
- Do not use while sleeping.
- This pad should not be used on or by an invalid, sleeping or unconscious person, or a person with poor blood circulation unless carefully attended.
- Do not use on areas of insensitive skin.
- Burns can occur regardless of control setting. Check skin under pad frequently.
- Do not use pad with liniment, salve or ointment preparations that contain heat-producing ingredients. This can result in skin burns.
- Do not tuck the heating pad or supply cord in a sofa or other furniture. All portions of the heating pad or controller should be visible at all times.
- Do not trap or pinch the cord.
- Do not use this heating pad on animals. Keep all pets away.
- Do not crush pad. Avoid sharp folds.
- This heating pad is intended for use on top of the body. Do not sit or lie on top of the heating pad. Never place pad between yourself and a chair, sofa, bed or pillow.
- Loop cord loosely when storing. Tight wrapping may damage cord and internal parts.
- This product is intended solely for household use and not for medical or commercial purposes. Not for hotel, motel or other institutional use.

SAVE THESE INSTRUCTIONS

i NOTICE

- Only use the heating pad together with the control unit specified on the label. Otherwise the heating pad could be damaged.
- With the heating pad switched on, never place any sharp object or any heating source such as a hot-water bottle, heating pad or similar items on it.
- The electronic components inside the control unit will heat up whilst operating the heating pad. Therefore never cover the control unit or place it onto the heating pad during operation.
- Protect the heating pad from knocks, damp, dust, chemicals, marked temperature fluctuations and nearby sources of heat (ovens, heaters).

2. Package Contents

- 1x Heating Pad UHP 27
- 1x Controller with Power Cord
- 1x Instruction Manual

3. Parts and Controls

4. Intended Use

This heating pad is intended for applying heat to the human body. It is not suitable for use in hospitals or for any other commercial use. In particular, it should not be used for babies or infants, persons who are insensitive to heat, helpless persons, or animals. Any use not described above, or modifications of the heating pad are not permissible and may lead to injury and /or damage to the heating pad. The manufacturer assumes no liability for damage due to improper use.

5. To Set Up


NOTICE

1. To begin, first connect the control connector to the heating pad connector.
2. Plug the power plug into a 120V AC electrical outlet.

6. To Use

Switching On

Set the controller to the highest temperature setting (setting 6). This will cause the heating pad to warm up quickly. The heat indicator light illuminates after the heating pad is switched on. When using the heating pad for several hours, set it to the lowest temperature setting (setting 1).

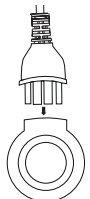Switching Off

Set the controller to 0 to switch off the heating pad. The heat indicator light will switch off.

Auto-off feature

The heating pad switches itself off automatically after about 120 minutes. The heat indicator light then starts to flash. To switch on the heating blanket again, select the setting 0 and after about 5 seconds, reset the heating pad to the desired temperature setting.

If the heating pad is no longer used after the timer expires:

- Switch off the heating pad (setting 0),
- Unplug from electrical outlet,
- Disconnect the controller from the heating pad.

Safety System

The heating pad is equipped with a safety system. The electronic sensor technology prevents the entire surface of the heating pad from overheating by automatically switching off in the event of a fault. If the safety system switches the heating pad off automatically due to a fault, then the heating indicator light on the controller is no longer illuminated, even if the heating pad is switched on. Please note that, for safety reasons, the heating pad should not be used after a fault and must be returned to the service address indicated. Never connect the defective heating pad to another controller of the same type. This would also result in a permanent shutdown by the safety system in the controller.

7. Care, Maintenance and Disposal

7.1 Washing Instructions

Heating Pad

If the heating pad is slightly soiled, use a dampen cloth and some mild detergent, if necessary. Follow these washing instructions to preserve the softness of the fabric and protect the electrical components: Wash by hand or as follows in automatic washer using slow agitation and normal spin speed: Fill washer with lukewarm water. Add minimum amount of all-purpose detergent, agitate to dissolve, measure heating pad before washing. Add heating pad and soak five minutes. Agitate two minutes. Spin. Fill for cool rinse, agitate one minute then spin. Dry in automatic dryer on high heat, five minutes only, or line dry. Remove damp heating pad from dryer, block to original size, and drape to finish drying. Do not use clothes pins. Do not iron your heating pad. Do not use wringer. Do not dry clean, use dry-cleaning fluids, cleaning solvents or bleach as it may have a deteriorating effect on the insulation of the heating element.

WARNING

- Disconnect the mains plug. Do NOT pull by the cord.
- Do not reconnect the controller to the heating pad until the control connector and the heating pad are completely dry. Otherwise there is a risk of electric shock.
- Never switch the heating pad on to help it dry. Otherwise there is a risk of electric shock

Controller

Use a dry, lint-free cloth for controller cleaning and care.

7.2 Storage

- Let the heating pad cool before putting it away. The heating pad may otherwise be damaged.
- To prevent the heating pad from becoming sharply creased, do not store other objects on top of the heating pad. The heating pad may otherwise be damaged.
- If you do not intend to use the heating pad for any length of time, store it in the original packaging in a dry environment and without weighing it down.
- Disconnect the mains connection and the controller from the heating pad to do so.

7.3 Disposal

Dispose of the device in accordance with local regulations. If you have any questions, please contact the local authorities responsible for waste disposal.

8. Technical Specifications

Model:	UHP 27
Heating Pad Size:	12 x 24 in (30.5 x 61 cm)
Power supply:	120 V / 60 Hz
Wattage:	50 W
Pad Material:	100% polyester
Timer:	120 minutes
Temperature levels:	6 levels Level 0: OFF Level 1: minimum heat Level 2-5: medium heat Level 6: maximum heat
Operating Conditions:	32°F to 95°F (0°C - 35 °C) < 80% humidity
Storage Conditions:	23°F to 104°F (-5°C to 40°C) < 80% humidity

9. Troubleshooting Guide

Problem	Cause	Solution
The heat indicator light is not illuminated while <ul style="list-style-type: none">• The controller is properly connected to the heating pad.• The mains plug is connected to the working electrical outlet.• The sliding switch is set to 1,2,3,4,5 or 6.	The safety system has switched the heating pad off permanently.	Send the heating pad and the controller for servicing

10. FCC Compliance Information

Heating Pad UHP 27

Responsible Party – U.S. Contact Information:

Beurer North America LP

1 Oakwood Boulevard, Suite 255

Hollywood, FL 33020

United States

1-800-536-0366

info@beurer.com

FCC Compliance Statement

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

1. This device may not cause harmful interference, and
2. This device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

11. Warranty

5-Year Limited Warranty for Original Purchaser

Your Beurer Heating Pad, Model UHP 27, is warranted to be free from defects in materials and workmanship for a period of five (5) years from the date of purchase under normal conditions of intended use and service. This warranty extends only to the original retail purchaser and does not extend to retailers or subsequent owners. We will, at our option, repair or replace the Beurer Heating Pad, Model UHP 27, without additional charge, for any part or parts covered by this written warranty. No refunds will be given. Repair or replacement is our only responsibility and your only remedy under this written warranty. If replacement parts for defective materials are not available, Beurer reserves the right to make product substitutions in lieu of repair or replacement. For warranty service, contact our customer service department at 1-800-536-0366 or info@beurer.com with a description of the problem. If the issue is deemed to fall within the scope of the 5-Year Limited Warranty, you will be asked to mail the product (at your cost) in its original package along with proof of purchase, your name, address, and phone number. If the issue is not deemed eligible under warranty, we will provide a quotation for repair or replacement, including return shipping fees.

Third-Party Repairs

This warranty complies with the federal Magnuson–Moss Warranty Act. You may use third-party parts or services for repairs without voiding this warranty, unless the repair itself causes damage to the product. We encourage you to retain your proof of purchase and contact our customer service team with any questions.

Unauthorized Sellers

This warranty applies only to products purchased from Beurer or authorized retailers and distributors. Products purchased from unauthorized sellers — including third-party platforms or online marketplaces (such as eBay or unverified Amazon sellers) — are not covered by this warranty. This policy is in place to protect customers from counterfeit, altered, or used products that may not meet Beurer's quality standards. If you're unsure whether your seller is authorized, please contact our support team.

Exclusions

This warranty does not cover damage caused by misuse or abuse; accident; the attachment of unauthorized accessories; alterations; improper installation; misapplication; lack of reasonable care; unauthorized repairs or modifications; improper electrical/power supply; worn batteries; normal wear and tear; power loss; dropped product; failure to follow instructions; lack of recommended maintenance; transit damage; theft; neglect; vandalism; or environmental conditions.

This warranty is void if the product is used in a commercial or business environment.

Beurer's maximum liability under this warranty is limited to the actual purchase price paid, as confirmed by proof of purchase. The warranty is effective only in the country where the product was purchased and used.

Legal Notice

This warranty gives you specific legal rights, and you may have other rights which vary from state to state. THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. ANY IMPLIED WARRANTIES, INCLUDING BUT NOT LIMITED TO MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO THE LENGTH OF THIS WRITTEN WARRANTY. IN NO EVENT SHALL BEURER BE LIABLE FOR ANY INCIDENTAL OR CONSEQUENTIAL DAMAGES ARISING FROM USE OR INABILITY TO USE THE PRODUCT.

Beurer does not authorize any third party — including retailers or resellers — to modify the terms of this warranty.

This warranty does not apply to opened, used, repackaged, or resold products, including those sold through online auctions or unauthorized channels.

For more information regarding our product line in the USA, please visit: www.shop-beurer.com

Distributed by:

Beurer North America LP

1 Oakwood Boulevard, Suite 255

Hollywood, FL 33020, USA

www.shop-beurer.com

Questions or comments? Call our US-based customer service toll free at 1-800-536-0366.

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