**READ THIS MANUAL COMPLETELY AND CAREFULLY BEFORE
USING THIS PRODUCT**

Keep this manual in a safe location for future reference

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 **Questions or comments? Call our US-based customer service toll free at 1-800-536-0366.**

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CONTENT

1. IMPORTANT SAFETY INSTRUCTIONS	2	6. Care, Maintenance and Disposal	7
2. Package Contents	4	7. Technical Specifications	7
3. Parts and Controls.....	5	8. Troubleshooting Guide	8
4. Intended Use	5	9. FCC Compliance Information	8
5. To Use	6	10. Warranty	9

SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE

If the instruction manual is damaged or if you no longer have the instruction manual in your possession, please contact Beurer customer service. See warranty for service contact.

Dear Valued Customer,

Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for application in the areas of heating, weight management, blood pressure, body temperature, pulse, gentle therapy, massage, and beauty.

Register your product today to receive further benefits:

- added protection for your purchase
- health and wellness tips
- new product information
- continual outstanding customer service and support

Visit www.registerbeurer.com or contact customer service at 1-800-536-0366.

Please read these instructions for use carefully and keep them for later use.

Make them accessible to other users and observe the information they contain.

With kind regards,

Your Beurer team

1. IMPORTANT SAFETY INSTRUCTIONS - DO NOT DESTROY

READ THIS ENTIRE MANUAL, THE SAFETY SECTION AND ALL INSTRUCTIONS AND WARNINGS COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT. FOLLOW ALL SAFETY INSTRUCTIONS AND WARNINGS TO AVOID HAZARDOUS SITUATIONS AND TO MAKE CORRECT USE OF THIS PRODUCT.

Signs and symbols

Whenever used, the following signs identify safety and property damage messages and designate a level of hazard or seriousness.

	This is the safety alert symbol. It alerts you to potential personal injury hazards. Obey all safety messages that follow this symbol to avoid possible injury or death.
 WARNING	WARNING indicates a hazardous situation which, if not avoided, could result in death or serious injury.
 CAUTION	CAUTION indicates a hazardous situation which, if not avoided, may result in minor or moderate injury.
 NOTICE	NOTICE addresses practices not related to personal injury, such as product and/or property damage.
<p>Storage</p> 	Permissible storage and transport temperature and humidity.
<p>Operating</p> 	Permissible operating temperature and humidity.
	Observe the instructions for use.
	Do not use the device with a pacemaker.

WARNING – To reduce the risk of burns, electric shock and fire, this product must be used in accordance with the following instructions:

- Read all instructions before using the product.
- Do not use on infants.
- Keep children away from packaging materials (risk of suffocation).
- Never leave Microplush Heated Pillow unattended.
- Do not use pins or other metallic means to fasten this in place.
- Never pull this Microplush Heated Pillow by the power cord. Do not use the cord as a handle. Unplug when not in use.
- Do not use a Microplush Heated Pillow when symptoms of appendicitis are present. Consult your physician.
- Inspect this product before using it in order to make sure that there are no broken parts, missing parts and/or damaged areas. Do not use Microplush Heated Pillow if there are broken parts, missing parts or damaged areas.
- This appliance has a polarized plug (one blade is wider than the other). As a safety feature, this plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not attempt to defeat this safety feature.
- Do not use an electrical outlet that has become loose or does not engage the power plug tightly.
- Do not immerse in water or use chemicals or solvents during cleaning, except as directed in the owner's manual.
- Do not remove product labels that contain warnings or safety instructions.

- Consult with a physician for use with controlled substances.
- Do not tuck in the wired area of this Microplush Heated Pillow.
- Do not route electrical cords between the mattress and boxspring. Do not allow the cords to be pinched.
- Keep the controller away from drafty areas that may become damp or wet such as an open window.
- Do not use this Microplush Heated Pillow if you or the product is wet.
- Be sure that this Microplush Heated Pillow is used on an alternating current supply circuit of the proper voltage.
- Turn off and unplug when not in use.
- Indoor use only.
- With the Microplush Heated Pillow switched on, never place any sharp object or any heating source such as a hot-water bottle, Microplush Heated Pillow or similar items on it.
- The electronic components inside the control unit will heat up whilst operating the Microplush Heated Pillow. Therefore never cover the control unit or place it onto the Microplush Heated Pillow during operation.
- The electronic and magnetic fields emitted by this Microplush Heated Pillow may interfere with the function of a pacemaker. However, they are still well below the limits for electrical field strength (5000 V/m max), magnetic field strength (80 A/m max), and magnetic flux density (0.1 millitesla max). Please consult your doctor and the manufacturer of your pacemaker before using this Microplush Heated Pillow.
- Household plug use only.

CAUTION – To reduce the risk of personal injury or product/property damage:

- Do not use while sleeping.
- This Microplush Heated Pillow should not be used on or by persons who are unconscious, persons with poor blood circulation, or persons with diabetes, unless they are carefully monitored.
- Do not use on areas of insensitive skin.
- Burns can occur regardless of control setting. Check skin under pillow frequently.
- Do not use Microplush Heated Pillow with liniment, salve or ointment preparations that contain heat-producing ingredients. This can result in skin burns.
- Do not use with alcohol.
- Do not tuck the Microplush Heated Pillow or supply cord in a sofa or other furniture. All portions of the Microplush Heated Pillow or controller should be visible at all times.
- Do not trap or pinch the cord.
- Protect the Microplush Heated Pillow from knocks, damp, dust, chemicals, marked temperature fluctuations and nearby sources of heat (ovens, heaters).
- Do not use this Microplush Heated Pillow on animals. Keep all pets away.

NOTICE


- Do not crush Microplush Heated Pillow. Avoid sharp folds.
- Loop cord loosely when storing. Tight wrapping may damage cord and internal parts.
- The Microplush Heated Pillow is not suitable for use in hospitals or for any other commercial use.
- Not for Hotel, Motel or other institutional use.
- Not for medical use in hospitals.
- Only use the Microplush Heated Pillow together with the control unit specified on the label. Otherwise the Microplush Heated Pillow could be damaged.

2. Package Contents

- 1x Microplush Heated Pillow UHP 48
- 1x Controller with Power Cord
- 1x Instruction Manual

3. Parts and Controls

1. Controller
2. Microplush Heated Pillow (cover and inner pad)
3. Power plug
4. Cord
5. Control Connector
6. Power button for ON/OFF and temperature settings
7. Heat indicator light
8. Microplush Heated Pillow Connector (on the inside of the cover)

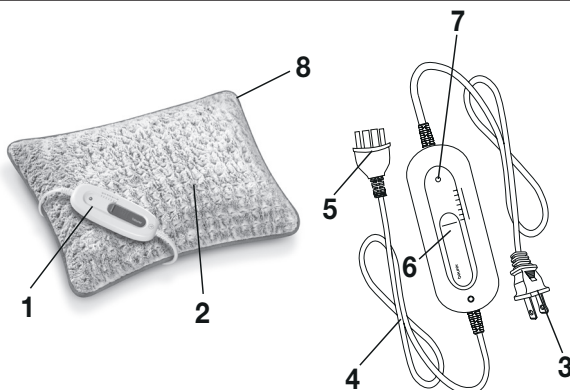4. Intended Use

This Microplush Heated Pillow is intended for applying heat to the human body. IHeat therapy is a natural way to relieve joint pain and pain associated with muscle tension and stress. It has been clinically proven that applying heat to a vulnerable area helps soothe sore muscles, relieve cramps and spasms. Any use not described above, or modifications of the Microplush Heated Pillow are not permissible and may lead to injury and /or damage to the Microplush Heated Pillow. The manufacturer assumes no liability for damage due to improper use.

Safety System

The Microplush Heated Pillow is equipped with a safety system. The electronic sensor technology prevents the entire surface of the Microplush Heated Pillow from overheating by automatically switching off in the event of a fault. If the safety system switches the Microplush Heated Pillow off automatically due to a fault, then the heating indicator light on the controller is no longer illuminated, even if the Microplush Heated Pillow is switched on.

NOTICE


- Please note that, for safety reasons, the Microplush Heated Pillow must not be used after a fault and must be returned to the service address indicated.
- Never connect the defective Microplush Heated Pillow to another controller of the same type. This would also result in a permanent shutdown by the safety system in the controller.

5. To Use

NOTICE

- The Microplush Heated Pillow may smell of plastic the first time you use it, but this will disappear after a short period of time.

To begin, first connect the control connector to the Microplush Heated Pillow connector (on the inside of the cover).

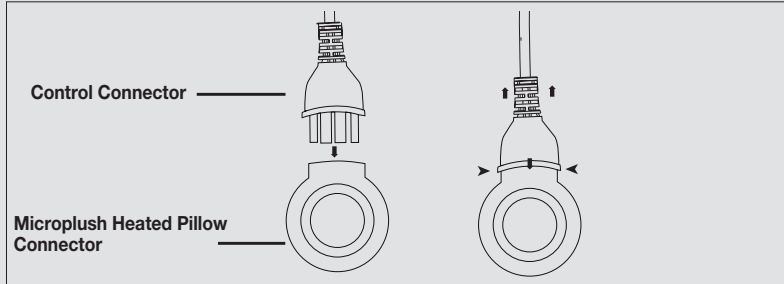The control connector must be connected to the Microplush Heated Pillow connector before plugging the unit into an electrical outlet. After connecting the Microplush Heated Pillow connector to the control connector, plug the unit into a 120V AC electrical outlet.

Switching On

Set the controller to 1, 2, 3, 4, 5 or 6 to switch on the Microplush Heated Pillow.

NOTICE

- The heat indicator light illuminates after the Microplush Heated Pillow is switched on.
- This Microplush Heated Pillow features rapid heating, which allows it to heat up within the first 10 minutes of operation.

Switching Off

Set the controller to 0 to switch off the Microplush Heated Pillow.

NOTICE

- The heat indicator light will switch off.

Selecting Temperature

Just before using the Microplush Heated Pillow, set it to the highest temperature setting (setting 6). This will cause the Microplush Heated Pillow to warm up quickly. When using the Microplush Heated Pillow for several hours, set it to the lowest temperature setting (setting 1).

Automatic shut-off

The Microplush Heated Pillow switches itself off automatically after about 120 minutes. The heat indicator light then starts to flash. Then select the setting 0 and after about 5 seconds, reset the Microplush Heated Pillow to the desired temperature setting to switch it on again.

If the Microplush Heated Pillow is no longer used after the timer expires:

- Switch off the Microplush Heated Pillow (setting 0),
- Unplug from electrical outlet,
- Disconnect the controller from the Microplush Heated Pillow.

6. Care, Maintenance and Disposal

Washing Instructions

DO NOT DRY CLEAN! DO NOT USE BLEACH! DO NOT USE WRINGER! DO NOT IRON YOUR HEATED PILLOW!

If the Microplush Heated Pillow is slightly soiled, use a dampen cloth and some mild detergent, if necessary. Use a dry, lint-free cloth for controller cleaning and care.

NOTE: Do not wash inner pad!

WARNING

1. Switch off the heated pillow (setting 0).
2. Disconnect power cord. **DO NOT** pull by the cord to unplug.
3. Disengage pillow/controller connector. Always remove the inner pad before cleaning. To do this undo the zip fastening. The inner pad is not washable.
4. Machine wash cover in warm water on gentle cycle; tumble dry on medium heat. Remove promptly from the dryer. Or, hand wash in warm water and line dry.
5. Reshape the still damp heated pillow to the original size immediately after washing, and dry flat, spread on a clothes drying rack
6. **DO NOT** reconnect the controller to the heated pillow until the control connector and the heated pillow are completely dry.

NOTICE

- Never immerse the controller in water or other liquids. It may otherwise be damaged.
- Do **NOT** use chemical cleaning or scouring agents to clean the Microplush Heated Pillow and controller. The Microplush Heated Pillow and controller may otherwise be damaged.
- Do not hang the Microplush Heated Pillow up to dry with pegs or similar items. The Microplush Heated Pillow may otherwise be damaged.

Storage

NOTICE

- Let the Microplush Heated Pillow cool before storing. The Microplush Heated Pillow may otherwise be damaged.
- Do not store other objects on top of the Microplush Heated Pillow. The Microplush Heated Pillow may otherwise be damaged.
- If you do not intend to use the Microplush Heated Pillow for any length of time, store it in the original packaging in a dry environment and without weighing it down.
- Disconnect the controller from the Microplush Heated Pillow.

Disposal

Observe the local regulations for material disposal. Dispose of the device in accordance with local regulations. If you have any questions, please contact the local authorities responsible for waste disposal.

7. Technical Specifications

Model:	Microplush Heated Pillow UHP 48
Type:	U-SP16RC
Microplush Heated Pillow Size:	12 x 16 in (40 cm x 30 cm)
Power supply:	120 V / 60 Hz
Wattage:	50 W
Material:	100% polyester

Timer:	120 minutes
Temperature levels:	6 levels Level 0: OFF Level 1: minimum heat Level 2-5: medium heat Level 6: maximum heat
Operating Conditions:	32°F to 95°F (0°C - 35 °C) < 80% humidity
Storage Conditions:	23°F to 104°F (-5°C to 40°C) < 80% humidity

8. Troubleshooting Guide

Problem	Cause	Solution
The heat indicator light is not illuminated while <ul style="list-style-type: none"> • The controller is properly connected to the Microplush Heated Pillow. • The power plug is connected to the working electrical outlet. • The sliding switch is set to 1,2,3,4,5 or 6. 	The safety system has switched the Microplush Heated Pillow off permanently.	Send the Microplush Heated Pillow and the controller for servicing

9. FCC Compliance Information

Microplush Heated Pillow UHP 48

Responsible Party – U.S. Contact Information:

Beurer North America LP

1 Oakwood Boulevard, Suite 255

Hollywood, FL 33020

United States

1-800-536-0366

info@beurer.com

FCC Compliance Statement

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

1. This device may not cause harmful interference, and
2. This device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.

- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

10. Warranty

Limited Lifetime Warranty For Original Purchaser

Your Beurer Microplush Heated Pillow, model UHP 48, is warranted to be free from defects in materials and workmanship for the life of the product under normal conditions of intended use and service. This warranty extends only to the original retail purchaser and does not extend to retailers or subsequent owners.

We will, at our option, repair or replace the Beurer Microplush Heated Pillow, model UHP 48, without additional charge, for any part or parts covered by these written warranties. No refunds will be given. Repair or replacement is our only responsibility and your only remedy under this written warranty. If replacement parts for defective materials are not available, Beurer reserves the right to make product substitutions in lieu of repair or replacement.

For warranty service contact our customer service department at 1-800-536-0366 or at info@beurer.com to provide a description of the problem. If the problem is deemed to be within the scope of the limited lifetime warranty, you will be asked to mail the product at your costs in its original package with proof of purchase, your name, address and phone number. If the problem is not deemed to be within the scope of the limited lifetime warranty, we will provide a quotation for repair respectively replacement and return shipping fee.

This warranty does not cover damage caused by misuse or abuse; accident; the attachment of unauthorized accessory; alteration to the product; improper installation; misapplication; lack of reasonable care with respect to the product; unauthorized repairs or modifications; improper use of electrical/power supply; old worn batteries; normal wear; loss of power; dropped product; malfunction or damage of an operating part as a result of failure to comply with instructions for use or to provide manufacturer's recommended maintenance; transit damage; theft; neglect; vandalism; or environmental conditions; loss of use during the period the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions whatsoever that are beyond the control of Beurer. This warranty is void if the product is ever used in a commercial or business environment. The maximum liability of Beurer under this warranty is limited to the purchase price actually paid by the customer for the product covered by the warranty, as confirmed by proof of purchase, regardless of the amount of any other direct or indirect damage suffered by the customer.

This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. A product that requires modifications or adaptation to enable it to operate in any other country than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under this warranty.

THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO THE DURATION OF THIS APPLICABLE WRITTEN WARRANTY. Some states do not allow limitations on how long an implied warranty lasts, so the above limitations may not apply to you.

IN NO EVENT SHALL BEURER BE LIABLE FOR ANY SPECIAL, INCIDENTAL, INDIRECT OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS, IMPLIED OR ANY OTHER THEORY OF LIABILITY, WHATSOEVER. Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.

Beurer does not authorize anyone, including, but not limited to, retailers, the subsequent consumer purchaser of the product from a retailer or remote purchaser, to obligate Beurer in any way beyond the terms set forth herein.

This warranty does not extend to the purchase of opened, used, repaired, repackaged and/or resealed products, including but not limited to sale of such products on Internet auction sites and/or products by surplus or bulk resellers. Any and all warranties or guarantees shall immediately cease and terminate in connection with any products or parts thereof which are repaired, replaced, altered, or modified, without the prior explicitly written consent of Beurer.

This warranty gives you specific legal rights, and you may also have other rights which may vary from state to state.

For more information regarding our product line in the USA, please visit: www.beurer.com

Distributed by:

Beurer North America LP

1 Oakwood Boulevard, Suite 255

Hollywood, FL 33020, USA

www.beurer.com

Questions or comments? Call our US-based customer service toll free at 1-800-536-0366.

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