beurer

BR 60

EN Insect bite healer	
Instructions for use	3
ES Curador de picaduras de insectos	
Instrucciones de uso	26
FR Appareil soulageant les piqûres d'insectes	
Mode d'emploi	50

For Canada only/ Solo para Canadá/
Pour le Canada seulement :
Manufactured by / Fabricado por / Fabriqué par :
Beurer GmbH, Söflinger Str. 218,
89077 Ulm, Germany

Distributed by/ Distribuido por/ Distribué par :
Beurer North America LP,
1 Oakwood Boulevard, Suite 255,
Hollywood, FL 33020, USA
Assembled in Indonesia / Ensamblado en
Indonesia / Assemblé en Indonésie
www.shop-beurer.com

**Questions or Comments? Call our
US-based customer service toll free at
1-800-536-0366.**

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**READ THIS MANUAL COMPLETELY AND CAREFULLY
BEFORE USING THIS PRODUCT**

Keep this manual in a safe location for future reference

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ENGLISH

Contents

1. IMPORTANT SAFETY NOTES.....	4	7. Care, Maintenance and Disposal	16
2. What to know about your device	13	8. Troubleshooting Guide.....	17
3. Package Contents.....	13	9. Technical Specifications.....	19
4. Device overview	14	10. FCC Compliance Information.....	20
5. Setting up before use	15	11. Warranty	22
6. Using the insect bite healer.....	15		

Read these instructions for use carefully. Observe the warnings and safety notes. Keep these instructions for use for future reference. Make the instructions for use accessible to other users. If the device is passed on, provide the instructions for use to the next user as well.

1. IMPORTANT SAFETY NOTES

Signs and symbols

Whenever used, the following signs identify safety and property damage messages and designate a level of hazard or seriousness.

	This is the safety alert symbol. It alerts you to potential personal injury hazards. Obey all safety messages that follow this symbol to avoid possible injury or death.
	Product information Note on important information
 WARNING	Indicates a potentially hazardous situation, which, if not avoided, could result in death or serious injury.
 CAUTION	Indicates a potentially hazardous situation, which, if not avoided, may result in minor/moderate injury.
NOTICE	Addresses practices not related to personal injury, such as product and/or property damage.
	Follow instructions for use

	Type BF applied part
I	On
O	Off
SN	Serial number
IP42	IP class Device protected against solid foreign objects ≥ 1 mm and against water dripping at an angle
Storage/Transport 	Permissible storage and transport temperature and humidity
Operating 	Permissible operating temperature and humidity

WARNING

Before using the device, read the safety notes carefully and follow the instructions exactly whenever using the device. Non-observance of the following instructions can result in injury or material damage. Store the instructions for use for future reference and make them accessible to other users. Any person using it on themselves must be over 12 years old. The product should only be used on children under 12 years old if supervised or assisted by an adult.

- This device may be used by persons with limited physical, sensory or mental capabilities, provided they are supervised by a responsible person or they receive detailed instructions on how to use the device.
- Do not use the device if you suffer from diabetes. Sensitivity to pain could, in certain circumstances, be reduced and diabetics may burn themselves.
- Do not use the device on people with disabilities, children under three years old or people with reduced sensitivity to heat. This includes people with skin alterations due to illness or scarred tissue in the application area, after taking pain relief medication, drinking alcohol, or taking drugs.
- Do not use the device on acutely or chronically diseased (injured or irritated) skin (e.g. inflamed skin, whether painful or not, reddened skin, rashes, (e.g. allergies),

burns, frostbite, bruises, swellings, both open and healing wounds, and post-operative scars where the healing process could be affected).

- Do not use the device if you are suffering from a fever.
- Do not use the device on sensitive areas of skin.
- Do not use the device if you have any sensory impairment that reduces the feeling of pain (e.g. metabolic disorders).
- Do not use the device if you are using lotions, creams, and gels locally at the same time.
- Do not use the device if you suffer from persistently irritated skin due to long-term heat application on the same area of skin.
- Before using the device, consult your doctor if you have any skin conditions, serious illness (in particular if you have a propensity to thromboembolic conditions or recurrent malignant growths), or unexplained chronic pain in any part of the body.
- Seek medical advice immediately if the symptoms persist or worsen despite using the insect bite healer. Stop using the device immediately.
- Seek medical advice immediately for a tick bite. The bite may result in transmission of pathogens such as tick-borne encephalitis (TBE) or borreliosis.

- Seek medical advice immediately at the first signs of increased itching, skin redness, swelling, dizziness, nausea, shortness of breath, formation of welts, or a marked change in blood pressure.
- Seek medical advice immediately if the symptoms worsen whilst or after using the insect bite healer and stop using it.
- Do not use near the eyes, on the eyelids, on the mouth, or on mucous membranes.
- Do not allow children to play with the device.
- Store the device out of the reach of children and pets.
- Keep packaging material away from children (risk of suffocation).
- The device is not a substitute for medical consultation and treatment.
 Always consult your doctor if you are experiencing any pain or are suffering from an illness.
- If the device does not work properly, or if you feel unwell or experience pain, stop using it immediately.

CAUTION

- Do not touch the device with wet hands. Do not allow water to spray onto the device. Only operate the device if it is completely dry.

- The device should only be used by one person; use by multiple people is not recommended.
- Depending on the treated person's skin sensitivity, the heat from the device may feel unpleasant and cause slight skin redness and, in rare cases, irritated skin. If the heat function feels too hot, stop treatment immediately
- Protect the device from impacts, moisture, dust, chemicals, marked temperature fluctuations and nearby sources of heat (ovens, heaters).
- Ensure that no water enters the device. If this should happen, remove the batteries and only use the device again when it is completely dry.
- If you store the device, no heavy objects should be placed on top of it.
- Do not leave the device unattended while in use. If the device malfunctions, switch it off and stop using it.

NOTICE

- Do not use any detergents or solvents to clean and maintain the device. If liquid penetrates the device this could cause damage to the electronics or other device parts and lead to a malfunction.
- Protect the device against direct sunlight and store in a dry place.

- Before use, ensure that there is no visible damage to the device or accessories. In case of any doubt, do not use the device, and contact your retailer or the specified Customer Service address.
- Store the device in a location protected against climatic influences. The device must be stored in the environmental conditions specified.
- No user-serviceable parts inside; do not repair or adjust the device yourself. Opening the device case will void the warranty.

Battery Handling Safety Precautions

- Use only the size and type of batteries specified.
- If your skin or eyes come into contact with battery fluid, rinse the affected areas with water and seek medical assistance.
- Be sure to follow the correct polarity when installing the batteries. Reversed batteries may cause damage to the device.
- Do not mix different types of batteries together (e.g. Alkaline and Carbon-zinc or rechargeable batteries) or old batteries with fresh ones. Always replace batteries as a simultaneous set.

- If the batteries in the device are depleted or the device will not be used for a long period of time, remove the batteries to prevent damage or injury from possible battery leakage.
- Do not try to recharge batteries not intended to be recharged; they can overheat and rupture (follow battery manufacturer's directions.)
- Do not dispose of batteries in fire, batteries may explode or leak.
- Clean the battery contacts and also those of the device prior to battery installation.
- Remove discharged batteries from the product and dispose/recycle in compliance with all applicable laws.
- Keep batteries away from children and pets. Batteries may be harmful if swallowed. Should a child or pet swallow a battery, seek medical assistance immediately.

Notes on electromagnetic compatibility

- The device is suitable for use in all environments listed in these instructions for use, including domestic environments.
- The use of the device may be limited in the presence of electromagnetic disturbances. This could result in issues such as error messages or the failure of the display/device.
- Avoid using this device directly next to other devices or stacked on top of other devices, as this could lead to faulty operation. If, however, it is necessary to use the device in the manner stated, this device as well as the other devices must be monitored to ensure they are working properly.
- The use of accessories other than those specified or provided by the manufacturer of this device can lead to an increase in electromagnetic emissions or a decrease in the device's electromagnetic immunity; this can result in faulty operation.
- Failure to comply with the above can impair the performance of the device.

2. What to know about your device

Intended Use

For household use only. The device uses heat to relieve itching and swelling from insect bites and stings. Do not use this product for any other purpose.


How it works

Many insect bites and stings result in painful and swollen bite / sting marks which can often cause severe itching. Remove any insect stinger before use. For best result, treat affected area immediately. Later treatment may still help reduce itching and speed healing but will be less effective. The insect bite healer applies heat to the affected area which relieves itching and discomfort. It is often noted that the treated area heals more comfortably.

3. Package Contents

- 1 x insect bite healer device BR 60
- 2 x batteries, type AAA (LR03 1.5 V)
- 1 x instruction manual

4. Device overview

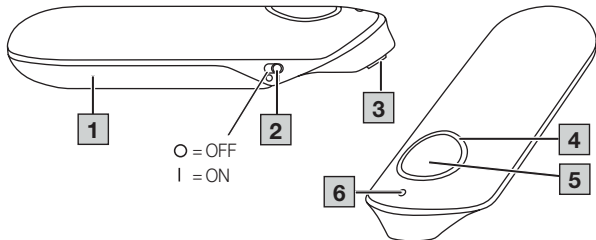1 Battery Compartment Cover

2 Power Switch

3 Heating Plate

4 Ring Light

5 Activation Button

6 Status LED

5. Setting up before use

1. Slide the Battery Compartment open and place the two “AAA” batteries (included) according to the polarity markings then slide the battery compartment cover closed until it clicks.
2. Make sure the device is clean, dry, and not damaged before use.

6. Using the insect bite healer

1. Slide the power switch to the “I” (ON) position to turn on the device. The Status LED will illuminate blue.
2. Place the white circular heating plate on the affected area.
3. Quick press the activation button. When treatment is complete the device will beep, and both the heating plate and ring light will shut off automatically.
 - Quick press once for sensitive skin (3-second treatment, only half of the ring light will illuminate)
 - Quick press twice for normal skin (6-second treatment, the full ring light will illuminate)

NOTE: Holding the button down, turns off the device. If this happens use the power switch to turn it off and on again.

4. Always slide the power switch to the “**○**” (OFF) position to turn off the device after use to save battery life.

NOTE: The device automatically enters sleeping mode after 2 minutes to prevent accidental use. To restart, use the power switch to turn it off and on again.

7. Care, Maintenance and Disposal

Battery Levels

When the battery is low, the Status LED will flash and the device will beep three times. Replace both batteries as instructed in “Setting up before use”.

Cleaning and Storage

- Allow the heating plate to cool completely before storing..
- Slide the power button to “**○**” (OFF) position to turn off the device and remove the batteries before cleaning.
- Wipe the heating plate with an alcohol-based cleaner after each use and between uses on multiple areas on the same person.

- Keep the device dry and away from bathtubs, showers, sinks, or other water containers.
- Avoid using substances for cleaning that could potentially be toxic for the skin and mucous membranes.
- Do not use abrasive cleaners or mechanical cleaning devices such as brushes.
- For long term storage, clean the device and keep it in a dry location out of direct sunlight.

Disposal

Follow local regulations for disposal. If you have any questions please contact the local authorities responsible for waste disposal.

8. Troubleshooting Guide

Problem	Solution
If the device does not turn on	Ensure the power switch is in the “I” (ON) position.
	Make sure the batteries are not depleted.

Problem	Solution
If the device does not emit a sufficient heat pulse	Battery level is low; replace batteries.
Device beeps three times and Status LED flashes	Battery level is low; replace batteries.
Device is on and batteries replaced, but no heat pulse is emitted during treatment	Device has turned off automatically. Slide the power switch to “●” (OFF) position to turn it off and then to the “I” (ON) position to turn it back on again.
Device beeps five times and Ring Light flashes	Device may be damaged or faulty; Turn device off and wait 10 minutes before restarting. If the device does not begin working normally upon restart, contact Customer Service.

9. Technical Specifications

Type	BR 60
Max. temperature	122 °F $\pm$ 3.6 °F (50 °C $\pm$ 2 °C)
Operating Duration per Battery Set	Approx. 300 applications
Power supply	2 x 1.5 V AAA LR03 batteries
Dimensions (L x W x H)	3.85 x 1.4 x 0.8 in (98 x 36 x 20 mm)
Weight (with batteries)	1.65 oz (47 g)
Functions	2 application times (3s / 6s)
Operating Ambient Temperature Range	50 - 104° F (10 - 40° C)
Storage Temperature Range	14 - 122° F (-10 - 50° C)
Operating and Storage Humidity Range	30 - 85%
Automatic device switch-off	2 min after last application

The serial number is located on the device or in the battery compartment.
Subject to technical changes.

10. FCC Compliance Information

Insect Bite Healer BR 60

Responsible Party U.S. Contact Information:

Beurer North America LP

1 Oakwood Boulevard, Suite 255

Hollywood, FL 33020

United States

1-800-536-0366

info@beurer.com

FCC Compliance Statement

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

11. Warranty

5-Year Limited Warranty for Original Purchaser

Your Beurer Insect Bite Healer, Model BR 60, is warranted to be free from defects in materials and workmanship for a period of five (5) years from the date of purchase under normal conditions of intended use and service. This warranty extends only to the original retail purchaser and does not extend to retailers or subsequent owners.

We will, at our option, repair or replace the Beurer Insect Bite Healer, Model BR 60, without additional charge, for any part or parts covered by this written warranty. No refunds will be given. Repair or replacement is our only responsibility and your only remedy under this written warranty. If replacement parts for defective materials are not available, Beurer reserves the right to make product substitutions in lieu of repair or replacement.

For warranty service, contact our customer service department at 1-800-536-0366 or info@beurer.com with a description of the problem. If the issue is deemed to fall within the scope of the 5-Year Limited Warranty, you will be asked to mail the product (at your cost) in its original package along with proof of purchase, your name, address, and phone number. If the issue is not deemed eligible under warranty, we will provide a quotation for repair or replacement, including return shipping fees.

Third-Party Repairs

This warranty complies with the federal Magnuson–Moss Warranty Act. You may use third-party parts or services for repairs without voiding this warranty, unless the repair itself causes damage to the product. We encourage you to retain your proof of purchase and contact our customer service team with any questions.

Unauthorized Sellers

This warranty applies only to products purchased from Beurer or authorized retailers and distributors. Products purchased from unauthorized sellers — including third-party platforms or online marketplaces (such as eBay or unverified Amazon sellers) — are not covered by this warranty. This policy is in place to protect customers from counterfeit, altered, or used products that may not meet Beurer's quality standards.

If you're unsure whether your seller is authorized, please contact our support team.

Exclusions

This warranty does not cover damage caused by misuse or abuse; accident; the attachment of unauthorized accessories; alterations; improper installation; misapplication; lack of reasonable care; unauthorized repairs or modifications; improper electrical/power supply; worn batteries; normal wear and tear; power loss; dropped

product; failure to follow instructions; lack of recommended maintenance; transit damage; theft; neglect; vandalism; or environmental conditions.

This warranty is void if the product is used in a commercial or business environment.

Beurer's maximum liability under this warranty is limited to the actual purchase price paid, as confirmed by proof of purchase. The warranty is effective only in the country where the product was purchased and used.

Legal Notice

This warranty gives you specific legal rights, and you may have other rights which vary from state to state.

THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. ANY IMPLIED WARRANTIES, INCLUDING BUT NOT LIMITED TO MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO THE LENGTH OF THIS WRITTEN WARRANTY.

IN NO EVENT SHALL BEURER BE LIABLE FOR ANY INCIDENTAL OR CONSEQUENTIAL DAMAGES ARISING FROM USE OR INABILITY TO USE THE PRODUCT.

Beurer does not authorize any third party — including retailers or resellers — to modify the terms of this warranty.

This warranty does not apply to opened, used, repackaged, or resold products, including those sold through online auctions or unauthorized channels.

For more information regarding our product line in the USA, please visit:
www.shop-beurer.com

For Canada only:
Manufacturer by:
Beurer GmbH
Söflinger Str. 218
89077 Ulm, Germany

Distributed by:
Beurer North America LP
1 Oakwood Boulevard, Suite 255
Hollywood, FL 33020, USA
www.shop-beurer.com

Questions or comments? Call our US-based customer service toll free at
1-800-536-0366.

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