# beurer

## MG 153

|                 |                                           |
|-----------------|-------------------------------------------|
| <b>ENGLISH</b>  | <b>4D Neck Massager</b>                   |
|                 | Instruction for Use..... 2                |
| <b>ESPAÑOL</b>  | <b>Masajeador de cuello 4D</b>            |
|                 | Instrucciones de uso..... 9               |
| <b>FRANÇAIS</b> | <b>Appareil de massage de la nuque 4D</b> |
|                 | Mode d'emploi ..... 17                    |

**READ THIS MANUAL COMPLETELY AND CAREFULLY  
BEFORE USING THIS PRODUCT**

Keep this manual in a safe location for future reference

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Questions or Comments? Call our US-based customer service  
toll free at 1-800-536-0366.

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## Contents

|                                      |   |                                         |   |
|--------------------------------------|---|-----------------------------------------|---|
| 1. IMPORTANT SAFETY NOTES .....      | 2 | 6. Care, Maintenance and Disposal ..... | 6 |
| 2. Getting to know your device ..... | 4 | 7. Troubleshooting Guide .....          | 6 |
| 3. Package Contents .....            | 4 | 8. Technical Specifications .....       | 6 |
| 4. Parts and Controls .....          | 5 | 9. FCC Compliance Information .....     | 6 |
| 5. To Use .....                      | 5 | 10. Warranty .....                      | 7 |

## SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE

If the instruction manual is damaged or if you no longer have the instruction manual in your possession, please contact customer service. See warranty for service contact.

### Dear Valued Customer

Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for application in the areas of heating, weight management, blood pressure, body temperature, pulse, gentle therapy, massage, and beauty.

Register your product today to receive further benefits:

- added protection for your purchase
- health and wellness tips
- new product information
- continual outstanding customer service and support.

Visit [www.registerbeurer.com](http://www.registerbeurer.com) or contact customer service at 1-800-536-0366.

Please read these instructions for use carefully and keep them for later use.

Make them accessible to other users and observe the information they contain.

With kind regards,  
Your Beurer team

## 1. IMPORTANT SAFETY NOTES

### Signs and symbols

Whenever used, the following signs identify safety and property damage messages and designate a level of hazard or seriousness.

**READ THIS ENTIRE MANUAL, THE SAFETY SECTION AND ALL INSTRUCTIONS AND WARNINGS COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT. FOLLOW ALL SAFETY INSTRUCTIONS AND WARNINGS TO AVOID HAZARDOUS SITUATIONS AND TO MAKE CORRECT USE OF THIS PRODUCT.**

|                                                                                                   |                                                                                                                                                                          |
|---------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
|                 | This is the safety alert symbol. It alerts you to potential personal injury hazards. Obey all safety messages that follow this symbol to avoid possible injury or death. |
|  <b>WARNING</b> | Indicates a potentially hazardous situation, which, if not avoided, could result in death or serious injury.                                                             |
|  <b>CAUTION</b> | Indicates a potentially hazardous situation, which, if not avoided, may result in minor/moderate injury.                                                                 |
|                | Direct current                                                                                                                                                           |

|                                                                                   |                                                                          |
|-----------------------------------------------------------------------------------|--------------------------------------------------------------------------|
|   | The unit is double-protected and thus corresponds to protection class 2. |
|  | Only use indoors.                                                        |
|  | Polarity                                                                 |
|  | Energy efficiency level 6                                                |
|  | The cover is machine-washable.                                           |
|  | Do not bleach                                                            |
|  | Do not tumble dry                                                        |
|  | Do not iron                                                              |
|  | Do not dry clean                                                         |

**WARNING – To reduce the risk of fire, electric shock, or serious personal injury:**

This device is designed for massaging the neck, back, shoulders and thighs of the human body. It does not replace medical treatment. Do not use this massager if any of the following warnings apply to you. If you are unsure whether the massager is suitable for you, consult your doctor.

- Packaging materials are a deadly hazard for children and can cause suffocation. Remove all packaging materials immediately and keep them away from children at all times.
- Do not use this device:
  - On animals.
  - If you suffer from a medical abnormality or injury in the area to be massaged (e.g. herniated disc, open wounds).
  - To massage the heart region of persons fitted with a pacemaker. Always consult a doctor before using the massage device on other parts of the body of a person fitted with a pacemaker.
  - On swollen, burned, inflamed, or damaged skin or body parts.
  - On lacerated wounds, capillaries, varicose veins, acne, herpes, or other skin diseases.
  - During pregnancy.
  - On the face, eyes, throat or other particularly sensitive parts of the body.
  - For foot massage.
  - While you are asleep.
  - After taking medication or consuming alcohol (reduced levels of awareness).
  - In the car.
- Consult your doctor before using this device if you:
  - Suffer from a serious illness or are recovering from an operation on the upper body.
  - Are fitted with a cardiac pacemaker, implant or other aid.
  - Have thromboses.
  - Are diabetic.
  - Have pains of unclarified origin.

- To help prevent the risk of electric shock:
  - Use only with the supplied power supply adapter and with the appropriate electrical voltage.
  - Do not use if damage to the device itself or any of its accessories is evident.
  - Never use the appliance during a thunderstorm.
- In the event of a malfunction, switch the device off immediately and disconnect it from the power supply. Never pull on the power cord or on the device itself to disconnect the power supply from the electrical outlet. Never hold or carry the device by the power cord. Keep the cord away from hot surfaces. Do not pinch, kink or twist the cables. Do not stick in any pins, needles or sharp objects.
- Ensure the massage device and power supply do not come into contact with water or other liquids. Therefore, use the device only in dry interior rooms (e.g. never in the bath or sauna).
- Never try to retrieve the device if it falls into water. Disconnect the plug immediately from the electrical outlet.
- Do not use the device if damage to the device itself or any of its accessories is evident.
- Do not subject the device to any impact and do not drop it.
- Repairs to electrical appliances must only be carried out by qualified persons. Incorrect repairs can lead to considerable danger for the user. In the event of repairs, please contact customer service.
- Do not open the zipper on the device. It is included only for ease of manufacture.
- To help prevent the risk of fire, do not use the device:
  - Under covers such as blankets or pillows
  - Near flammable substances such as gasoline
- This device is for private use only and is not intended for medical or commercial purposes.

**CAUTION - To reduce the risk of personal injury or product/property damage:**

- Shut the device off and unplug the power cord after every use and before cleaning.
- Never insert any objects into openings on the device, and never insert anything into the rotating parts. Make sure that all moving parts can move freely at all times.
- Never sit or stand on the moving parts of the device with all your weight and do not place objects on the device.
- Keep the device away from high temperatures.
- Repairs are to be performed only by an authorized service center. Contact Customer Service with any questions regarding operation, usage, or repair of this product.
- The zip fastening on the massager must not be opened. It is only used for technical production reasons.

## 2. Getting to know your device

The rollers and pinch nodes in this neck massager imitate the pressure and kneading movements of a massage. Using this electric massager provides a soothing, effective massage with no help required. The device offers a powerful, intensive massage for the neck, back, shoulders and thighs.


In addition, this massager features:

- Simple operation
- Heat function
- Two speed settings
- Removable machine washable cover

## 3. Package Contents

1 x 4D Neck Massager  
 1 x Power Adapter  
 1 x Instruction Manual


## 4. Parts and Controls

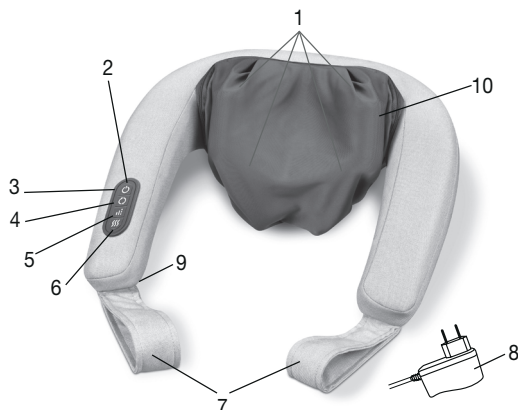1. Massage Heads (4 massage heads, 2 heated pinch nodes)
2. Control Panel
3. On/Off Button 
4. Rotation Selector Button 
5. Intensity Selector Button 
6. Heat Button 
7. Handles
8. Power Adapter
9. Power Cord
10. Removable Washable Cover

## 5. To Use

1. Remove all parts from the packaging and inspect for damage.
2. Connect the Power Adapter to the Power Cord on the massager and route it so that it will not be tripped over or pulled.
3. Press the body area to be treated carefully against the massager and check that the massage position is comfortable for you. The massage should feel pleasant and relaxing at all times. If the massage is painful or uncomfortable, stop the massage or change the pillow's position or the contact pressure. Do not trap or squeeze any part of your body, especially feet, hands, and fingers, between the rotating massage heads and their mounts.
4. Press the On/Off Button to switch the massager on; the massage heads will start to rotate and the pinch nodes will move up and down. Press the Rotation Selector Button again to change the direction of rotation.
5. Press the Heat Button to activate the heat function; press the Heat Button again to cancel.
6. Press the On/Off Button again to shut the massager off.
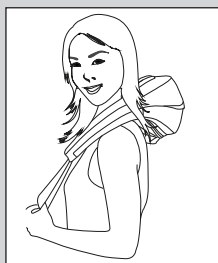
### Neck massage:

- Place massager around your neck.
- Adjust pressure using the handles

### Lower Back Massage

- Place massager around hips
- Adjust pressure using the handles

The above examples show possible ways to use the massager. Depending on your personal needs, the massager can also be used on the upper back area, thighs, and calves.

#### **CAUTION:**

- Do not use the massage device for longer than 15 minutes. Longer usage can over-stimulate the muscles and cause tension instead of relaxation. This device features an auto shutoff function, which will shut the device off after 15 minutes.

## **6. Care, Maintenance and Disposal**

### **WARNING**

- Before cleaning, unplug the power adapter from the electrical outlet and disconnect it from the device.
- Only clean this device as described. Keep water and other liquids out of the device and accessories.
- Small marks can be removed using a cloth or damp sponge and mild detergent. Do not use solvent-based cleaning products. Let dry completely before reuse.
- The removable cover can be machine-washed in warm water. Please note the symbols for care and cleaning on the label sewn into the cover. Set the washing machine to a delicate cycle at 86°F/ 30°C.
- Do not use the device again until it is completely dry.
- If not using the device for a long period, store it in the original packaging in a dry location. Do not place any heavy objects on it while in storage.

### **Disposal**

Observe the local regulations for material disposal. Dispose of the device in accordance with local regulations. If you have any questions, please contact the local authorities responsible for waste disposal.

## **7. Troubleshooting Guide**

| <b>Problem</b>                                  | <b>Cause</b>                                | <b>Solution</b>                                                                                           |
|-------------------------------------------------|---------------------------------------------|-----------------------------------------------------------------------------------------------------------|
| <b>Massage heads rotating at reduced speed.</b> | The load on the massage heads is too great. | Reduce pressure on the massage heads.                                                                     |
| <b>Massage heads not moving.</b>                | Device not plugged in.                      | Plug the power adapter to an electrical outlet after ensuring it is connected to the device's power cord. |

## **8. Technical Specifications**

|              |                                          |
|--------------|------------------------------------------|
| Type         | MG 153                                   |
| Power Supply |                                          |
| Input:       | 100-240 V~, 50/60 Hz, 0.7 A              |
| Output:      | 12.0 V === 2.0 A, 24.0 W                 |
| Dimensions   | 21.65 x 15.35 x 7.5 in (55 x 39 x 19 cm) |
| Weight       | About 4 Lbs (1.8 Kg)                     |

## **9. FCC Compliance Information**

4D Neck Massager MG 153

Responsible Party U.S. Contact Information:

Beurer North America LP

1 Oakwood Boulevard, Suite 255

Hollywood, FL 33020, United States

1-800-536-0366, info@beurer.com

## **FCC Compliance Statement**

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

## **10. Warranty**

### **Limited Lifetime Warranty For Original Purchaser**

Your Beurer 4D Neck Massager, model MG 153, is warranted to be free from defects in materials and workmanship for the life of the product under normal conditions of intended use and service. This warranty extends only to the original retail purchaser and does not extend to retailers or subsequent owners.

We will, at our option, repair or replace the Beurer 4D Neck Massager, model MG 153, without additional charge, for any part or parts covered by these written warranties. No refunds will be given. Repair or replacement is our only responsibility and your only remedy under this written warranty. If replacement parts for defective materials are not available, Beurer reserves the right to make product substitutions in lieu of repair or replacement.

For warranty service contact our customer service department at 1-800-536-0366 or at [info@beurer.com](mailto:info@beurer.com) to provide a description of the problem. If the problem is deemed to be within the scope of the limited lifetime warranty, you will be asked to mail the product at your costs in its original package with proof of purchase, your name, address and phone number. If the problem is not deemed to be within the scope of the limited lifetime warranty, we will provide a quotation for repair respectively replacement and return shipping fee.

This warranty does not cover damage caused by misuse or abuse; accident; the attachment of unauthorized accessory; alteration to the product; improper installation; misapplication; lack of reasonable care with respect to the product; unauthorized repairs or modifications; improper use of electrical/power supply; old worn batteries; normal wear; loss of power; dropped product; malfunction or damage of an operating part as a result of failure to comply with instructions for use or to provide manufacturer's recommended maintenance; transit damage; theft; neglect; vandalism; or environmental conditions; loss of use during the period the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions whatsoever that are beyond the control of Beurer. This warranty is void if the product is ever used in a commercial or business environment. The maximum liability of Beurer under this warranty is limited to the purchase price actually paid by the customer for the product covered by the warranty, as confirmed by proof of purchase, regardless of the amount of any other direct or indirect damage suffered by the customer.

This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. A product that requires modifications or adaptation to enable it to operate in any other country than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under this warranty.

**THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO THE DURATION OF THIS APPLICABLE WRITTEN WARRANTY.** Some states do not allow limitations on how long an implied warranty lasts, so the above limitations may not apply to you.

**IN NO EVENT SHALL BEURER BE LIABLE FOR ANY SPECIAL, INCIDENTAL, INDIRECT OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS, IMPLIED OR ANY OTHER THEORY OF LIABILITY, WHATSOEVER.** Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.

Beurer does not authorize anyone, including, but not limited to, retailers, the subsequent consumer purchaser of the product from a retailer or remote purchaser, to obligate Beurer in any way beyond the terms set forth herein.

This warranty does not extend to the purchase of opened, used, repaired, repackaged and/or resealed products, including but not limited to sale of such products on Internet auction sites and/or products by surplus or bulk resellers. Any and all warranties or guarantees shall immediately cease and terminate in connection with any products or parts thereof which are repaired, replaced, altered, or modified, without the prior explicitly written consent of Beurer.

This warranty gives you specific legal rights, and you may also have other rights which may vary from state to state.

For more information regarding our product line in the USA, please visit: [www.beurer.com](http://www.beurer.com)

Distributed by:  
Beurer North America LP  
1 Oakwood Boulevard, Suite 255  
Hollywood, FL 33020, USA  
[www.beurer.com](http://www.beurer.com)

Questions or comments? Call our US-based customer service toll free at 1-800-536-0366.

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