**READ THIS MANUAL COMPLETELY AND CAREFULLY
BEFORE USING THIS PRODUCT**

Keep this manual in a safe location for future reference

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Contents

1. Important safety notes 2	8. Additional information 7
2. Important information 3	9. Care and Maintenance 7
3. Package Contents 4	10. Technical Specifications 8
4. Introduction 4	11. Troubleshooting Guide 8
5. Parts and Controls 5	12. FCC Compliance Information 9
6. Before First Use 6	13. IC Compliance Statement 10
7. To Use 6	14. Warranty 10

SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE

If the instruction manual is damaged or if you no longer have the instruction manual in your possession, please contact Beurer customer service. Please see warranty for service contact.

Dear Valued Customer,

Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for application in the areas of heating, weight management, blood pressure, body temperature, pulse, gentle therapy, massage, and beauty.

Register your product today to receive further benefits:

- added protection for your purchase
- health and wellness tips
- new product information
- continual outstanding customer service and support.

Visit www.registerbeurer.com or contact customer service at **1-800-536-0366**.

Please read these instructions for use carefully and keep them for later use. Make them accessible to other users and observe the information they contain.

With kind regards,
Your Beurer Team

1. Important safety notes

Signs and symbols

Whenever used, the following signs identify safety and property damage messages and designate a level of hazard or seriousness.

READ AND UNDERSTAND THIS ENTIRE MANUAL INCLUDING THE SAFETY SECTION AND ALL INSTRUCTIONS AND WARNINGS COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT. FOLLOW ALL SAFETY INSTRUCTIONS AND WARNINGS TO AVOID HAZARDOUS SITUATIONS AND TO MAKE CORRECT USE OF THIS PRODUCT.

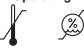
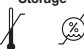
	This is the safety alert symbol. It alerts you to potential personal injury hazards. Obey all safety messages that follow this symbol to avoid possible injury or death.
 WARNING	Indicates a hazardous situation, which, if not avoided, could result in death or serious injury.
 CAUTION	Indicates a hazardous situation, which, if not avoided, may result in minor/moderate injury.

NOTICE	Addresses practices not related to personal injury, such as product and/or property damage.
---------------	---

2. Important information

Signs and symbols

The following symbols are used in these instructions for use, on the packaging and on the type plate for the device and accessories:

	Note Note on important information
Operating 	Permissible operating temperature and humidity
Storage 	Permissible storage temperature and humidity

WARNING

Safety notes

- Keep packaging material away from children (risk of suffocation).
- Swallowing batteries can be extremely dangerous. Keep the batteries and scale out of the reach of small children. Should a battery be swallowed, seek medical assistance immediately.
- Keep scale away from moisture, dust, chemicals, toiletries, liquid cosmetics, extreme hot and cold temperatures, and sources of heat (e.g., open flame, radiators).
- Prevent impact with hard objects.
- Do not press the button violently or with pointed objects.

IMPORTANT

- Do not place any object on the scale when it is not in use.
- Strong electromagnetic fields (e.g. cell phones) may impair the accuracy of the scale.
- Remove any transport locks

General notes

- The scale has a capacity of 44 lb (20 kg/705 oz). Measurements are displayed in 0.01 lb (5 g/0.18 oz) divisions.
- In case of differing measurements (between the scale and app), only use the measured values displayed on the scale.

CAUTION

- Cleaning: You can clean the scale with a damp cloth and a little washing up liquid, if required.
- Never immerse the scale in water or rinse it under running water.
- The scale is not intended to be used for commercial purposes.
- Repairs may only be performed by an authorized service center. Contact Customer Service with any questions regarding operation, usage, or repair of this product.

Battery Handling Safety Precautions

WARNING – To reduce the risk of serious personal injury:

- If your skin or eyes come into contact with battery fluid, rinse the affected area with water and seek medical assistance.
-  **Choking hazard!** Small children may swallow and choke on batteries. Store batteries out of the reach of small children.
- Keep batteries away from children. Only adults should handle batteries.
- If a battery has leaked, put on protective gloves and clean the battery compartment with a dry cloth.
- Protect batteries from excessive heat.
-  **Risk of explosion!** Never throw batteries into a fire. Dispose or recycle spent/discharged batteries in compliance with all applicable laws.
- Follow the battery manufacturer's safety and usage instructions.
- Do not charge or short-circuit batteries.

CAUTION

- Always use the size and type of battery indicated.
- Do not mix a used battery with a new battery, or mix alkaline, standard (carbon-zinc), or rechargeable (nickel-cadmium) batteries.
- Insert the batteries observing the proper polarity (+/-) as indicated.
- Always replace all batteries at the same time.
- Do not use rechargeable batteries.
- Do not disassemble, open or crush batteries.

NOTICE

- Clean the battery contacts and those of the device prior to installation.
- Remove spent or discharged batteries from the product.
- Remove the batteries when not using the product for an extended period.

3. Package Contents

- Baby Scale BY 90
- 3x 1.5 V AAA batteries (LR03)
- These operating instructions

4. Introduction

A baby's physical development requires a great deal of energy, which is why regular weight control is important; it puts parents' minds at ease to know their baby is getting sufficient nutrition and is developing properly. This scale ensures you can always keep an eye on your child's weight to ensure their healthy development.

This baby scale has the following features:

- Curved weighing surface.
- Large, easy-to-read LCD display.
- Automatic and manual hold function for recording baby's weight.
- Integrated tape measure (~59 in or 150 cm).
- Tare weighing.
- 44 lb (20 kg/705 oz) weight capacity.
- Automatic shutoff.
- Adjustable between weight units (kg/lb/oz).
- Storage of the last 10 measurements.
- Transfer of measured values to the „beurer BabyCare“ app.

- The product uses *Bluetooth*[®] low energy technology, frequency band = 2400 MHz - 2480 MHz. Max. 8.0 dBm transmission power. Compatible with *Bluetooth*[®] ≥ 4.0 smartphones/tablets.

System requirements


iOS version 10.0 or later; Android™ version 5.0 or later
Bluetooth[®] 4.0 or later

For a list of compatible devices, scan this QR Code with your mobile device:


5. Parts and Controls

Parts

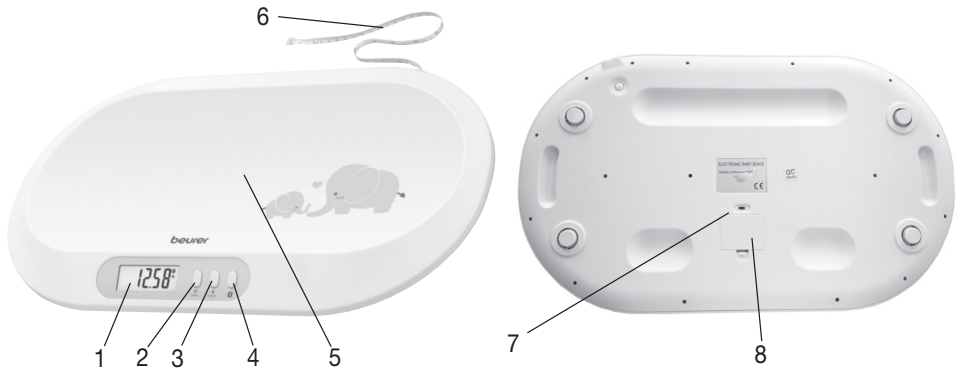1. LCD Display	5. Weighing surface
2. Power/Tare Button ⓘ	6. Tape measure
3. HOLD Button 🔒	7. Weight Selector Switch
4. <i>Bluetooth</i> [®] button ⓘ	8. Battery Compartment

Display

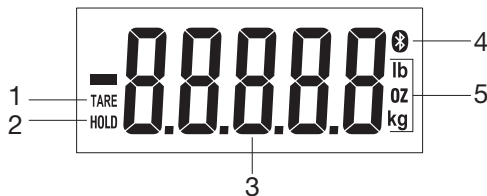1. „TARE“ (Tare function active)	4. <i>Bluetooth</i> [®] Pairing Icon
2. „HOLD“ (Hold function active)	5. Weight Unit
3. Weight Display	

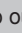

6. Before First Use

Remove the Battery Compartment Cover located underneath the scale and insert three „AAA“ size alkaline batteries (included) according to the polarity markings inside the compartment. Replace Battery Compartment Cover. When „L3“ appears in the LCD Display, the batteries need to be replaced. Note that when replacing batteries, you will need to pair the scale with your mobile device again.

After batteries are inserted properly, „APP“ will flash in the LCD Display. At this point, you have 120 seconds to pair the scale with your mobile device via *Bluetooth*®. If you do not pair the scale with your mobile device, note that „app“ will flash for 10 seconds whenever you first power on the scale. The indicator is no longer displayed once synchronisation with the app is successful.

7. To Use

You can transfer the measured weight values to your paired mobile device via *Bluetooth*® using the „beurer BabyCare“ app. Keep your mobile device close to the scale to maintain an active *Bluetooth*® connection during initial use. Press the Power/Tare Button to activate the scale. To use the scale with your mobile device, start with step 1; otherwise, press any button to cancel app use and skip to step 6.

1. Download the free „beurer BabyCare“ app from the Apple App Store (iOS) or from Google Play (Android™).
2. Activate *Bluetooth*® in your mobile device's settings.
NOTE: *Bluetooth*® devices, such as this scale, may not be visible in the general settings of your device's *Bluetooth*® device list.
3. Start the „beurer BabyCare“ app and follow the instructions.
4. Enter user data. Settings for baby's date of birth (Day, Month, and Year) must be entered into the „beurer BabyCare“ App.
5. The scale is successfully connected to the app when the *Bluetooth*® symbol  appears in the LCD Display.
6. When „0.00 lb“ appears in the LCD Display, the scale is ready for use.
7. Carefully place the baby onto the weighing surface, taking particular care with the baby's head.
8. If you do not want to place the baby directly on the weighing surface and prefer to use a towel or blanket, place it on the scale first and then press the Tare Button  to „zero out“ the LCD Display; „Tare“ appears in the LCD Display and the weight resets to zero. Now place the baby onto the scale.
9. This scale has both automatic and manual „HOLD“ functions. The auto hold function displays the baby's weight when it has stabilized enough to detect an exact value; it remains displayed for 120 seconds.
NOTE: If the weight cannot be determined automatically after 8 seconds, the baby scale averages the weight over 4 seconds (flashing display of the weight on the baby scale display).
10. If the baby is moving around too much, press the HOLD Button when the baby is as calm as he/she will become. Within six seconds, the word „HOLD“ stops flashing and the averaged weight is displayed.
11. To shut the auto hold function off completely, press and hold the HOLD Button for three seconds. „HOLD“ and „Off“ will appear in the LCD Display. To reactivate the auto hold function, press and hold the HOLD Button again; „HOLD“ and „On“ will appear in the LCD Display.
12. To transfer the measured weight to the app (if you are using the *Bluetooth*® function), open the app and press the *Bluetooth*® Button on the front panel of the scale. The scale will shut off when data transfer is complete.
13. If you decide not to transfer the measured weight to the app or are not using the *Bluetooth*® function, the scale will shut off after 120 seconds. To shut off the scale manually, press and hold the Power/Tare Button  for three seconds.

NOTE: Measurements saved on the baby scale before it is paired with your mobile device cannot be transferred.

Subsequent synchronization with the app

If you have not synchronized the baby scale with the „beurer BabyCare“ app during initial use, you can also do this later.

1. Press the Power/Tare Button  to power the scale on.
2. „App“ flashes in the LCD Display for 10 seconds. If the app is open, it will attempt to connect to the scale for 10 seconds.

Alternatively:

Press the *Bluetooth*[®] button whilst the baby scale is in operation. If the app is open, the baby scale can be connected for 120 seconds.

3. The baby scale is successfully connected to the app when the *Bluetooth*[®] symbol  appears in the LCD Display.

8. Additional information

Saving the measured values on the baby scale and transferring them to the app

A maximum of 10 measurements can be saved in the baby scale. The saved measurements are transferred automatically to the app if you open the app within the *Bluetooth*[®] range and *Bluetooth*[®] is switched on on the scale. It is not possible to automatically transfer data when the baby scale is switched off.

Deleting Measurements

To delete all measurements stored in the scale, press and hold the *Bluetooth*[®] Button for five seconds.

Using the tape measure

Pull out the tape measure on the top of the baby scale to measure the baby's length and other features. To retract, press the Button next to the tape measure on the bottom of the scale. Enter the length and head circumference of your baby into the „beurer BabyCare“ app.

Changing the batteries

Your scale is equipped with a low battery indicator. If you operate the scale with flat batteries, „L3“ will appear on the display and the scale will automatically switch off. In this case, the batteries must be replaced (3 x 1.5 V AAA). All saved measurements remain saved.

NOTE:

- You have to re-synchronise the baby scale with the app (see „Initial use with the app“ section) after replacing the batteries.
- When replacing the batteries, use batteries of the same type, make and capacity.
- Do not use rechargeable batteries.

9. Care and Maintenance

- The baby scale should be cleaned from time to time.
- Clean the scale with a damp cloth and a mild detergent. Do not immerse the scale in water or rinse it under running water.
- Store the scale in a cool, dry place when not in use; ensure no objects are placed on the scale during storage. If storing the scale for an extended period, remove the batteries.

IMPORTANT

- Never use harsh solvents or cleaning products!
- Do not submerge the baby scale in water.
- Do not rinse the baby scale under running water.

- Do not clean the baby scale in the dish washer.

10. Technical Specifications

Battery type:	3x 1.5 V AAA batteries (LR03)
Measuring range:	0.11 lb to 44 lb (0.05 kg to 20 kg/1.76 oz to 705 oz)
Weight unit:	0.01 lb (5 g/0.18 oz)
Temperature operating range:	50°F to 104°F (10°C to 40°C)
Storage temperature range:	14°F to 122°F (-10°C to 50°C)
Scale size:	22 in x 13 in x 1.85 in (560 x 330 x 47 mm)
Net weight of scale:	approx. 3 lb (1.36 kg)
Subject to technical modifications.	

Notes on electromagnetic compatibility

- The device is suitable for use in all environments listed in these instructions for use, including domestic environments.
- The use of the device may be limited in the presence of electromagnetic disturbances. This could result in issues such as error messages or the failure of the display/device.
- Avoid using this device directly next to other devices or stacked on top of other devices, as this could lead to faulty operation. If, however, it is necessary to use the device in the manner stated, this device as well as the other devices must be monitored to ensure they are working properly.
- The use of accessories other than those specified or provided by the manufacturer of this device can lead to an increase in electromagnetic emissions or a decrease in the device's electromagnetic immunity; this can result in faulty operation.
- Failure to comply with the above can impair the performance of the device.

11. Troubleshooting Guide

Display	Cause	Solution
Err	No flat, stable surface.	Place the scale on a flat, non-carpeted surface. Re-tare the scale.
	The maximum weight capacity of 44 lb (20 kg/705 oz) has been exceeded.	Only weigh the maximum permissible weight.
Incorrect weight displayed	No flat, stable surface.	Place the scale on a flat, non-carpeted surface. Re-tare the scale.
	The baby is wriggling too much.	Try to calm the baby so that it stops wriggling.
	Incorrect zero point of scale.	Switch the baby scale off and on again. Wait until the scale displays „0.00 lb“ and repeat the measurement.

No Bluetooth® connection 	Device outside range.	<i>Bluetooth</i> ® range is about 50 feet (~15m). Walls and ceilings reduce this range significantly. Radio waves may also cause interference, so do not place the scale near devices such as Wi-Fi routers and microwaves.
	<i>Bluetooth</i> ® is switched off on the scale.	Press the <i>Bluetooth</i> ® Button until the icon appears in the LCD Display.
	No connection to the app.	Completely close the app (including in the background). Switch <i>Bluetooth</i> ® off and back on again on the mobile device. Power your mobile device off and back on again. Remove the batteries from the scale for a short period and reinsert them. See FAQs at www.beurer.com .
Lo	Batteries in the scale are depleted.	Replace the batteries.
APP	Scale has not been paired with the app. Any measurements completed before pairing cannot be transferred.	Pair the scale with the app.

12. FCC Compliance Information

Baby Scale BY 90

Responsible Party – U.S. Contact Information:

Beurer North America LP

1 Oakwood Boulevard, Suite 255

Hollywood, FL 33020

United States

1-800-536-0366

info@beurer.com

FCC Compliance Statement

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

1. This device may not cause harmful interference, and
2. This device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.

—Consult the dealer or an experienced radio/TV technician for help.

13. IC Compliance Statement

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

1. This device may not cause interference.
2. This device must accept any interference, including interference that may cause undesired operation of the device.

14. Warranty

Limited Lifetime Warranty For Original Purchaser

Your Beurer Baby Scale, model BY 90, is warranted to be free from defects in materials and workmanship for the life of the product under normal conditions of intended use and service. This warranty extends only to the original retail purchaser and does not extend to retailers or subsequent owners. We will, at our option, repair or replace the Beurer Baby Scale, model BY 90, without additional charge, for any part or parts covered by these written warranties. No refunds will be given. Repair or replacement is our only responsibility and your only remedy under this written warranty. If replacement parts for defective materials are not available, Beurer reserves the right to make product substitutions in lieu of repair or replacement. For warranty service contact our customer service department at 1-800-536-0366 or at info@beurer.com to provide a description of the problem. If the problem is deemed to be within the scope of the limited lifetime warranty, you will be asked to mail the product at your costs in its original package with proof of purchase, your name, address and phone number. If the problem is not deemed to be within the scope of the limited lifetime warranty, we will provide a quotation for repair respectively replacement and return shipping fee. This warranty does not cover damage caused by misuse or abuse; accident; the attachment of unauthorized accessory; alteration to the product; improper installation; misapplication; lack of reasonable care with respect to the product; unauthorized repairs or modifications; improper use of electrical/power supply; old worn batteries; normal wear; loss of power; dropped product; malfunction or damage of an operating part as a result of failure to comply with instructions for use or to provide manufacturer's recommended maintenance; transit damage; theft; neglect; vandalism; or environmental conditions; loss of use during the period the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions whatsoever that are beyond the control of Beurer. This warranty is void if the product is ever used in a commercial or business environment. The maximum liability of Beurer under this warranty is limited to the purchase price actually paid by the customer for the product covered by the warranty, as confirmed by proof of purchase, regardless of the amount of any other direct or indirect damage suffered by the customer. This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. A product that requires modifications or adaptation to enable it to operate in any other country than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under this warranty.

THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO THE DURATION OF THIS APPLICABLE WRITTEN WARRANTY. Some states do not allow limitations on how long an implied warranty lasts, so the above limitations may not apply to you. **IN NO EVENT SHALL BEURER BE LIABLE FOR ANY SPECIAL, INCIDENTAL, INDIRECT OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS, IMPLIED OR ANY OTHER THEORY OF LIABILITY, WHATSOEVER.** Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you. Beurer does not authorize anyone, including, but not limited to, retailers, the subsequent consumer purchaser of the product from a retailer or remote purchaser, to obligate Beurer in any way beyond the terms set forth herein. This warranty does not extend to the purchase of opened, used, repaired, repackaged and/or resealed products, including but not limited to sale of such products on Internet auction sites and/or products by surplus or bulk resellers. Any and all warranties or guarantees shall immediately cease and terminate in connection with any products or parts thereof which are repaired, replaced, altered, or modified, without the prior explicitly written consent of Beurer.

This warranty gives you specific legal rights, and you may also have other rights which may vary from state to state.

For more information regarding our product line in the USA, please visit: www.beurer.com

Distributed by:

Beurer North America LP

1 Oakwood Boulevard, Suite 255

Hollywood, FL 33020

United States

www.beurer.com

Questions or Comments?

Call our US-based customer service toll free at 1-800-536-0366.

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Fabriqué en Chine.

The *Bluetooth*[®] word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by Beurer GmbH is under license. Other trademarks and trade names are those of their respective owners.

Apple and the Apple logo are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc., registered in the U.S. and other countries.

Google Play and the Google Play logo are trademarks of Google LLC.

Android is a trademark of Google LLC.