PO 13 Kids Pulse Oximeter

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READ THIS MANUAL COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT
Keep this manual in a safe location for future reference

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1. IMPORTANT SAFETY NOTES

Signs and symbols

Whenever used, the following signs identify safety and property damage messages and designate a level of hazard or seriousness.
READ THIS ENTIRE MANUAL, THE SAFETY SECTION AND ALL INSTRUCTIONS AND WARNINGS COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT. FOLLOW ALL SAFETY INSTRUCTIONS AND WARNINGS TO AVOID HAZARDOUS SITUATIONS AND TO MAKE CORRECT USE OF THIS PRODUCT.

⚠	This is the safety alert symbol. It alerts you to potential personal injury hazards. Obey all safety messages that follow this symbol to avoid possible injury or death.	**WARNING**	Indicates a potentially hazardous situation, which, if not avoided, could result in death or serious injury.
CAUTION	Indicates a potentially hazardous situation, which, if not avoided, may result in minor/moderate injury.	**NOTICE**	Addresses practices not related to personal injury, such as product and/or property damage.
ⓘ	**Note** Note on important information		No alarm function available
	Observe the instructions for use	⎓	Direct current
%SpO_2	Arterial oxygen saturation of haemoglobin (in percent)		Application part, type BF
PR bpm	Pulse rate (beats per minute)	SN	Serial number
Storage	Storage and transport conditions	**IP 22**	Device protected against foreign objects ≥ 12.5 mm and against falling drops of water
Operating	Operating conditions		

⚠ WARNING

Read these instructions for use carefully. Non-observance of the following information may result in personal injury or product/property damage. Store these instructions for use and make them accessible to other users. Make sure you include these instructions for use when handing over the device to third parties.

- Do not use the pulse oximeter:
 - if the device or the finger you are using is damp.
 - if you are allergic to rubber products.
 - while taking a blood pressure measurement on the same arm.
 - during an MRI or CT scan.
 - on fingers that have nail polish, are dirty, have other coatings on them, or have false nails applied.
 - on large fingers that do not fit into the device easily (fingertip width greater than 0.43 in / 11 mm).
 - on fingers with anatomical changes, edemas, scars or burns.
 - on fingers that are too small, as with small children for example (width approx. < 0.27 in / 7 mm).
 - on people who are not steady at the site of application (e.g. trembling).
 - near flammable or explosive gas mixtures.
- Extended use may cause pain for people with circulatory disorders. Do not use the pulse oximeter for longer than two hours on one finger.
- The pulse oximeter displays an instantaneous measurement but cannot be used for continuous monitoring.
- Measurements are for your information only and are not substitute for medical examination.
- The pulse oximeter is not for continuous monitoring and does not have alarm functions. Never use this unit in situations where alarms are required.
- Do not look directly inside the housing during the measurement. The red light and the invisible infrared light in the pulse oximeter are harmful to your eyes.
- Do not wrap the wrist strap around your neck to avoid an accident.
- This device is not intended for use by people (including children) with restricted physical, sensory or mental skills or a lack of experience and/or a lack of knowledge, unless they are supervised by a person who has responsibility for their safety or they receive instructions from this person on how to use the device. Children should be supervised around the device to ensure they do not play with it.
- Neither of the displays for the pulse wave and pulse bar allows the strength of the pulse or circulation to be evaluated at the measurement site. Rather, they are exclusively used to display the current visual signal variation at the measurement site and do not enable diagnostics for the pulse.
- In cases of carbon monoxide poisoning, the pulse oximeter will display a measurement value that is too high.

⚠ CAUTION

- Non-observance of the following instructions can lead to incorrect or failed measurements.
- Ensure that the finger nail on the finger to be measured is short enough that the fingertip covers the sensor element in the housing.
- Keep your hand, finger and body steady during the measurement.
- To avoid incorrect results, there should not be any strong light sources (e.g. fluorescent lamps or direct sunlight) in the immediate vicinity of the pulse oximeter.

NOTICE:

- Check the pulse oximeter regularly before each use to ensure that there is no visible damage to the device and the batteries are still sufficiently charged. In case of doubt, do not use the device and contact customer services or an authorized retailer.
- Do not use any additional parts that are not recommended by the manufacturer.
- Under no circumstances should you open or repair the device yourself. Failure to comply will result in voiding of the warranty. For repairs, please contact customer services or an authorized retailer.
- Protect the pulse oximeter from dust, shocks, moisture, extreme temperatures and explosive materials.

Battery Handling Safety Precautions

- Use only the size and type of batteries specified.
- Be sure to follow the correct polarity when installing the batteries. Reversed batteries may cause damage to the device.
- Do not mix different types of batteries together (e.g. Alkaline and Carbon-zinc or rechargeable batteries) or old batteries with fresh ones. Always replace batteries as a simultaneous set.
- If the batteries in the device are depleted or the device will not be used for a long period of time, remove the batteries to prevent damage or injury from possible battery leakage.
- Do not try to recharge batteries not intended to be recharged; they can overheat and rupture (follow battery manufacturer's directions).
- Do not dispose of batteries in fire, batteries may explode or leak.
- Clean the battery contacts and also those of the device prior to battery installation.
- Remove discharged batteries from the product and dispose/recycle in compliance with all applicable laws.
- Keep batteries away from children and pets. Batteries may be harmful if swallowed. Should a child or pet swallow a battery, seek medical assistance immediately.

Notes on electromagnetic compatibility

- The device is suitable for use in all environments listed in these instructions for use, including domestic environments.
- The use of the device may be limited in the presence of electromagnetic disturbances. This could result in issues such as error messages or the failure of the display/device.
- Avoid using this device directly next to other devices or stacked on top of other devices, as this could lead to faulty operation. If, however, it is necessary to use the device in the manner stated, this device as well as the other devices must be monitored to ensure they are working properly.
- The use of accessories other than those specified or provided by the manufacturer of this device can lead to an increase in electromagnetic emissions or a decrease in the device's electromagnetic immunity; this can result in faulty operation.
- Failure to comply with the above can impair the performance of the device.


2. INTENDED USE

Only use the Beurer PO 13 kids pulse oximeter on humans to measure the arterial oxygen saturation (SpO2) and the pulse rate. The pulse oximeter is suitable for private use (at home). This device is not intended for continuous monitoring.

3. PACKAGE CONTENTS

1 x PO 13 Kids pulse oximeter
1 x Wrist strap
1 x Belt pouch
1 x Small screwdriver
2 x 1.5V AAA alkaline batteries
1 x Instruction manual


4. PARTS AND CONTROLS

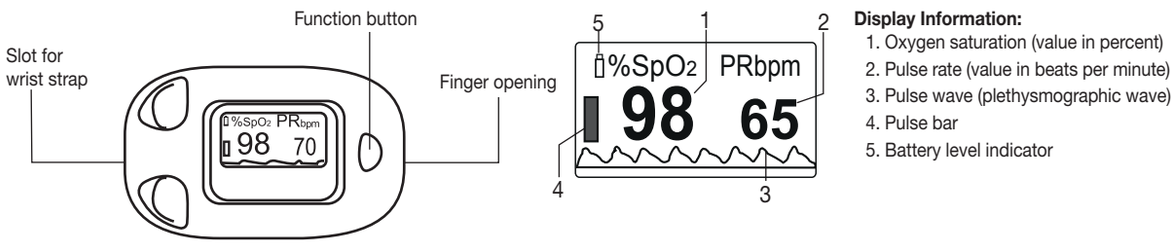Display Information:
1. Oxygen saturation (value in percent)
2. Pulse rate (value in beats per minute)
3. Pulse wave (plethysmographic wave)
4. Pulse bar
5. Battery level indicator

5. BEFORE FIRST USE


5.1 Inserting the batteries

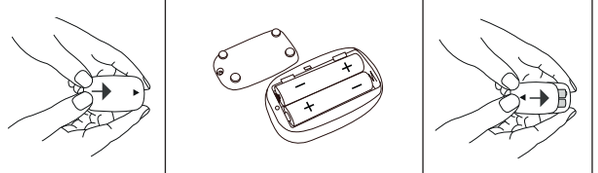1. Slide the battery compartment lid open.
2. Place two "AAA" alkaline batteries (included) into the battery compartment observing the correct polarity as indicated.
3. Reattach the battery compartment lid and tighten the screw.


5.2 Attaching the wrist strap

To transport the pulse oximeter more easily attach the wrist strap as follows:

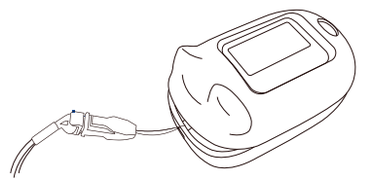1. Insert the narrow end of the wrist strap through the holes as shown.
2. Route the other end through the loop at the narrow end and tighten.

6. TO USE

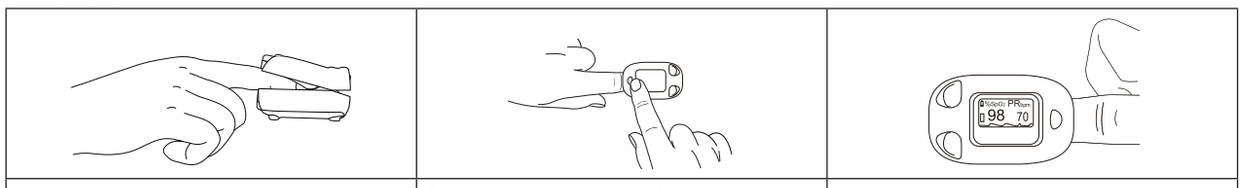1. Insert one finger into the opening as shown and hold it steady.
2. Press the function button. Do not move during the measurement.
3. Your measurement values will appear on the screen after a few seconds.

ⓘ Note

When you remove your finger from the pulse oximeter, the device will switch off after approx. five seconds.

Features

The pulse oximeter has the following features:
- **Power**: When the pulse oximeter is off, press the function button to switch it on.
- **Display**: The display format (vertical, horizontal etc.) changes automatically if you change the direction of the pulse oximeter.
- **Brightness:** To select your desired display brightness, press and hold the function button during operation until the brightness level changes.

7. CARE, MAINTENANCE AND DISPOSAL

- Clean the device with a slightly moistened cloth.
- Do not use detergents or solvents.
- Never immerse the device in water or other liquids.
- When storing the device, first remove the batteries. Make sure no heavy objects are placed on top of it while is it stored. Store in a dry location away from direct sunlight.

Disposal

Observe the local regulations for material disposal. Dispose of the device in accordance with local regulations. If you have any questions, please contact the local authorities responsible for waste disposal.

8. TROUBLESHOOTING GUIDE

Problem	Possible cause	Solution
Measurements do not display.	Batteries are depleted.	Replace the batteries.
	Batteries not inserted correctly.	Reinsert the batteries. If after reinserting the batteries correctly there are still no measurement values displayed, contact customer service.
Measurements are erratic.	Insufficient circulation in the measurement finger.	Observe the Important Safety Notes section.
	The finger is not properly inserted.	Please insert the finger fully into the device until you hit the rubber stop.
	Finger, hand or body is moving.	Keep your finger, hand and body still during the measurement.
	Cardiac arrhythmia.	Seek medical advice.

9. TECHNICAL SPECIFICATIONS

Type	PO 13
Measurement method	Non-invasive measurement of arterial oxygen saturation of hemoglobin and pulse rate in finger
Measurement range	SpO_2 0–100%, pulse 30–250 beats per minute
Accuracy	SpO_2 70–100%, ± 2%, pulse rate 30–250, ± 2 bpm
Dimensions	2.32 (L) x 1.45 (W) x 1.38 (H) in (59 (L) × 37 (W) × 35 (H) mm)
Weight	1.3 oz (36 g) without batteries
Sensor	Red light (wavelength 660 nm); infrared (wavelength 905 nm); silicon receiver diode
Operating conditions	50°–104° F (10°–40° C), ≤ 75% relative humidity, 700–1060 hPa ambient pressure
Storage conditions	-40°–140° F (-40°–60° C), ≤ 95% relative humidity, 500–1060 hPa ambient pressure
Power supply	2 x 1.5V AAA alkaline batteries
Battery life	Approx. one year for three 60-second measurements per day
Classification	IP22, application part, type BF

The serial number can be found on the unit itself or in the battery compartment.

10. FCC COMPLIANCE INFORMATION

Kids Pulse Oximeter PO 13

Responsible Party – U.S. Contact Information:
Beurer North America LP, 1 Oakwood Boulevard, Suite 255, Hollywood, FL 33020, United States
1-800-536-0366, info@beurer.com

FCC Compliance Statement

This device complies with Part 15 of the FCC Rules.
Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.
Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.
NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:
- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

11. WARRANTY

Limited Lifetime Warranty For Original Purchaser

Your Beurer Kids Pulse Oximeter, model PO 13, is warranted to be free from defects in materials and workmanship for the life of the product under normal conditions of intended use and service. This warranty extends only to the original retail purchaser and does not extend to retailers or subsequent owners.
We will, at our option, repair or replace the Beurer Kids Pulse Oximeter, model PO 13, without additional charge, for any part or parts covered by these written warranties. No refunds will be given. Repair or replacement is our only responsibility and your only remedy under this written warranty. If replacement parts for defective materials are not available, Beurer reserves the right to make product substitutions in lieu of repair or replacement.
For warranty service contact our customer service department at **1-800-536-0366 or at info@beurer.com** to provide a description of the problem. If the problem is deemed to be within the scope of the limited lifetime warranty, you will be asked to mail the product at your costs in its original package with proof of purchase, your name, address and phone number. If the problem is not deemed to be within the scope of the limited lifetime warranty, we will provide a quotation for repair respectively replacement and return shipping fee.
This warranty does not cover damage caused by misuse or abuse; accident; the attachment of unauthorized accessory; alteration to the product; improper installation; misapplication; lack of reasonable care with respect to the product; unauthorized repairs or modifications; improper use of electrical/power supply; normal wear; loss of power; dropped product; malfunction or damage of an operating part as a result of failure to comply with instructions for use or to provide manufacturer's recommended maintenance; transit damage; theft; neglect; vandalism; or environmental conditions; loss of use during the period the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions whatsoever that are beyond the control of Beurer. This warranty is void if the product is ever used in a commercial or business environment. The maximum liability of Beurer under this warranty is limited to the purchase price actually paid by the customer for the product covered by the warranty, as confirmed by proof of purchase, regardless of the amount of any other direct or indirect damage suffered by the customer.
This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. A product that requires modifications or adaptation to enable it to operate in any other country than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under this warranty.
THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO THE DURATION OF THIS APPLICABLE WRITTEN WARRANTY. Some states do not allow limitations on how long an implied warranty lasts, so the above limitations may not apply to you.
IN NO EVENT SHALL BEURER BE LIABLE FOR ANY SPECIAL, INCIDENTAL, INDIRECT OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS, IMPLIED OR ANY OTHER THEORY OF LIABILITY, WHATSOEVER. Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.
Beurer does not authorize anyone, including, but not limited to, retailers, the subsequent consumer purchaser of the product from a retailer or remote purchaser, to obligate Beurer in any way beyond the terms set forth herein.
This warranty does not extend to the purchase of opened, used, repaired, repackaged and/or resealed products, including but not limited to sale of such products on Internet auction sites and/or products by surplus or bulk resellers. Any and all warranties or guarantees shall immediately cease and terminate in connection with any products or parts thereof which are repaired, replaced, altered, or modified, without the prior explicitly written consent of Beurer.
This warranty gives you specific legal rights, and you may also have other rights which may vary from state to state.

Distributed by:
Beurer North America LP, 1 Oakwood Boulevard, Suite 255, Hollywood, FL 33020, USA
Assembled in China

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