beurer

HA 69 Pair

ENGLISH

In The Ear Sound Amplifiers HA 69

Instructions for use..... 3

ESPAÑOL

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Questions or Comments? Call our US-based customer service toll free at 1-800-536-0366.

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READ THIS MANUAL COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT

Keep this manual in a safe location for future reference

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ENGLISH

Contents

1. IMPORTANT SAFETY INSTRUCTIONS 5
2. Package Contents 12
3. Introduction 13
4. Parts and Controls..... 14
5. To Use 15
6. Care, Maintenance and Disposal 17
7. Troubleshooting Guide 19
8. Technical Specification..... 20
9. FCC Compliance Information..... 22
10. Warranty 24

SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE

If the instruction manual is damaged or if you no longer have the instruction manual in your possession, please contact customer service. See warranty for service contact.

Dear Valued Customer
Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for application in the areas of heating, weight management, blood pressure, body temperature, pulse, gentle therapy, massage, and beauty.

Register your product today to receive further benefits:
- added protection for your purchase
- health and wellness tips
- new product information
- continual outstanding customer service and support.

Visit www.registerbeurer.com or contact customer service at 1-800-536-0366.

Please read these instructions for use carefully and keep them for later use. Make them accessible to other users and observe the information they contain.

With kind regards,
Your Beurer team

1. IMPORTANT SAFETY INSTRUCTIONS

Signs and symbols

Whenever used, the following signs identify safety and property damage messages and designate a level of hazard or seriousness.

READ THIS ENTIRE MANUAL, THE SAFETY SECTION AND ALL INSTRUCTIONS AND WARNINGS COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT. FOLLOW ALL SAFETY INSTRUCTIONS AND WARNINGS TO AVOID HAZARDOUS SITUATIONS AND TO MAKE CORRECT USE OF THIS PRODUCT.

	This is the safety alert symbol. It alerts you to potential personal injury hazards. Obey all safety messages that follow this symbol to avoid possible injury or death.
 **WARNING**	Indicates a hazardous situation which, if not avoided, could result in death or serious injury.

⚠ **CAUTION**	Indicates a hazardous situation which, if not avoided, may result in minor or moderate injury.
NOTICE	Addresses practices not related to personal injury, such as product and/or property damage.
IP21	Protected against solid foreign objects 12.5 mm in diameter and larger, and against vertically falling drops of water

⚠ WARNING – To reduce the risk of fire, electric shock, or serious personal injury:

- Have your hearing checked by a doctor or qualified hearing care professional before using the sound amplifiers.
- In the event of an existing or earlier ear problem, you should consult your doctor before using these sound amplifiers.

- NOT permissible for use by patients suffering from an ear inflammation of the middle ear or ear canal or infection.
- If you feel unwell or experience pain when wearing these sound amplifiers, stop using them immediately and contact your doctor.
- Stop using the sound amplifiers and seek medical attention if the following occurs:
 - Rash or eczema on areas of skin that contact the sound amplifiers. (The sound amplifier is made from a material that is gentle on the skin. In rare cases, however, allergic reactions may occur that are similar to eczema.)
 - Use of the sound amplifiers results in discharge from the ear and/or necessitates ear treatment.
 - Sudden, significant decrease in sound ability.
- If parts of the sound amplifier become stuck in the ear, seek medical attention immediately. Do not attempt to remove the parts yourself.
- At first, set the volume to a low setting and increase it gradually once you have become accustomed to the sound amplifiers.
- Do NOT set the volume control too high, as this may cause sound damage.

- Do NOT lend your sound amplifiers to others or use sound amplifiers belonging to others, as there is a risk of sound damage and infection.
- Sound amplifiers are NOT suitable for use by children less than 36 months of age.
- Keep the sound amplifiers, the accessories and the batteries away from children, as there is a risk of injury and choking.
- This device is NOT intended for use by people (including children) with restricted physical, sensory or mental skills or a lack of experience and/or a lack of knowledge, unless they are supervised by a person who has responsibility for their safety or they receive instructions from this person on how to use the device. Supervise children around the device to ensure they do not play with it.
- Do NOT use the sound amplifiers in potentially explosive and hazardous environments.
- Do NOT use the devices in oxygen-rich environments and over an extended period.
- If you have any doubts, do NOT use a damaged sound amplifier; contact Customer Service.

- Switch the devices off immediately if they are faulty or not working properly.
- Do not attempt to open and/or repair the sound amplifiers yourself. Repairs may only be performed by Customer Service. Failure to comply will invalidate the warranty.

⚠ CAUTION – To reduce the risk of personal injury or product/ property damage:

- Do NOT wear the sound amplifiers in the vicinity of strong electromagnetic fields or X-rays.
- Do NOT wear the sound amplifiers while bathing, showering or swimming.
- Protect the devices from impact, humidity, dirt, extreme temperature fluctuations, direct sunlight and high humidity.
- High temperatures may develop in closed vehicles during summer. You should therefore never leave your sound amplifiers inside a vehicle.
- If the sound amplifiers fall into water, dry them with a cloth and then have it inspected by your retailer.

- Never attempt to dry the sound amplifiers in a microwave, a tumble dryer or with a hair dryer, as doing so may damage the sound amplifier.
- Do NOT use the sound amplifiers together with a mobile phone. (The radio waves of the mobile phone may cause background noise or a loss of volume.)
- Do NOT wear the sound amplifiers while applying cosmetics (e.g. make-up, hairspray, perfume or sunscreen).
- Take care when operating the controller. Do NOT press too firmly.
- Do not insert any sharp objects into the microphone opening as this may damage the microphone.
- Do NOT touch the magnets of the sound amplifiers as this may damage the devices.
- Do NOT touch the inside of the devices. There are many sensitive components that have been precisely configured. Under no circumstances should you poke inside the sound amplifiers with a needle or similar object.

2. Package Contents

Check that the exterior of the delivery cardboard packaging is intact and make sure that all components are present. Before use, ensure that the device and accessories show no visible damage and all packaging material has been removed. If you have any doubts, do not use the device and contact your retailer or the specified Customer Service address.

2 x Sound Amplifiers HA 69
6 x Exchangeable ear tips
1 x Cleaning Brush
1 x Cleaning Stick
1 x Charging Case
1 x USB Charging cable
1 x Instruction Manual

3. Introduction

This device amplifies ambient sounds and transmits them directly into your ear. It cannot prevent or alleviate hearing loss. It can be used to hear sounds louder and more clearly, such as nature sounds, listening to lectures with a distant speaker, and listening to soft sounds (e.g., distant conversations).

4. Parts and Controls

Sound amplifier, accessories and replacement parts

1. Volume Control
2. Microphone
3. Earpiece and Ear Tip
4. Charging Contacts
5. Charging Slot
6. Charge Indicator LEDs
7. Power LED
8. Micro USB Charging Port (on rear)
9. Cleaning Brush
10. Cleaning Stick
11. USB Charging cable


5. To Use

1. Select an appropriately sized set of Ear Tips and attach them to the Earpieces as shown below.

2. The Sound Amplifiers are right and left specific.There is a L or R marking on the sound amplifier. This will tell you if it is a left or a right sound amplifier; if you are unsure which is which, place them back into the Storage Case to check, since they fit into the case only one way.
3. Place the sound amplifiers in your ears so that the ear tip fits into the ear canal snugly; do not force it.
4. Turn the Sound Amplifier on by pressing the Volume Button at least three seconds. The Sound Amplifier will beep twice and the LED will illuminate in white for two seconds.

5. Press the Volume Button to increase the volume level; the Sound Amplifier will emit an increasing loud beep as you repeatedly press it, indicating the volume setting. When the volume level is at maximum, the Sound Amplifier will beep three times. Pressing the Volume Button again will return the volume level to its lowest setting.

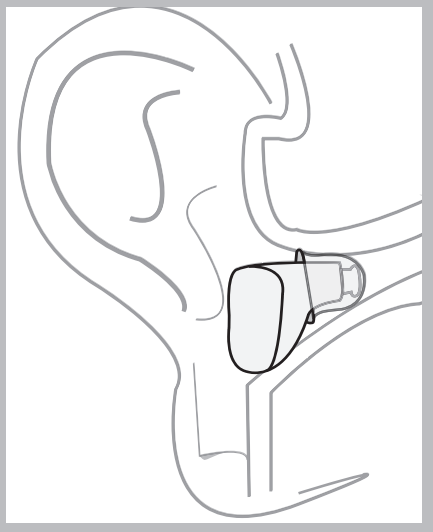6. If the Sound Amplifier emits a whistling sound, ensure the earpiece is properly fitted in your ear. If necessary, use a different ear tip. If sound levels are too low even at the higher volume settings, the batteries may need to be charged.
7. Turn the Sound Amplifier off by pressing the Volume Button at least three seconds. The Sound Amplifier will beep twice and the LED will illuminate in blue for two seconds.

Charging the Batteries

1. With the Sound Amplifiers inserted correctly into the Storage Case, connect the supplied Charging Cable to the Micro USB Charging Port on the rear of the Storage Case and to a 5V USB charging source, such as a computer or AC USB Charger (not included). The Charge Indicator LEDs and Power LED will illuminate in blue. If the Charge Indicator LEDs blink red or do not illuminate, the sound amplifiers are inserted incorrectly; check and reinsert.
2. When each sound amplifier is charged, its corresponding Charge Indicator LED will change to white.

6. Care, Maintenance and Disposal

NOTICE: Never let water enter the device.

1. Clean each sound amplifier with a soft, dry cloth after each use. If extra cleaning is needed, moisten the cloth with a mild detergent. Do not use abrasive or corrosive cleaners.
2. Let each sound amplifier dry completely before reuse.


3. The supplied cleaning brush and stick are suitable for removing earwax from the ear cushions and earpiece hole; see below for details.

Changing the Earwax Guard

Earwax is one of the main causes of problems with sound amplifiers. During use, earwax can get into your sound amplifiers and partially or completely block the passage of sound.

A cleaning brush and cleaning stick are included with this product to enable easy cleaning of the earpiece hole and ear tips.

1. Remove the ear tip.
2. Clean any residual ear wax using the cleaning brush and cleaning stick.
3. Reattach the ear tip after cleaning it.
4. Clean the brush and stick if needed, using a moistened cloth with ethyl alcohol.

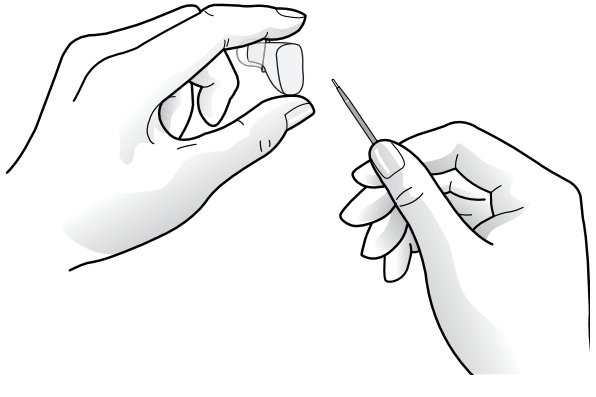Disposal

Observe the local regulations for material disposal. Dispose of the device in accordance with local regulations. If you have any questions, please contact the local authorities responsible for waste disposal.

7. Troubleshooting Guide

Problem	Possible causes
No sound	• Battery needs charging • Microphone is clogged/dirty • Sound amplifier is damaged
Sound is very weak/faint	• Volume setting on the sound amplifier is too low • Battery is weak • Earplug and/or Microphone is dirty or damaged • Ear Tip is dirty and needs to be cleaned or replaced

Sound is distorted	• Ear Tip not attached correctly • Volume too high; lower the volume • Switch the device to the other ear, although using the device in the specified ear is recommended whenever possible. • Make sure the microphone is not obscured by your ear or hair, as this may increase acoustic feedback
Device is not charged	• Device is not inserted correctly in the charging dock. • The device was not shut off before charging. Please keep in mind that it won‘t properly charge if not fully powered off.

8. Technical Specification

Max Sound Pressure Level	≤120 dB (+3 dB)
Max Amplification	≤42 dB (+3 dB)

Frequency Range	300 - 3400 Hz
Total Harmonic Wave Distortion	≤ 10 % + 3 %
Current	≤ 3 mA
Input Noise	≤ 30 dB + 3 dB
Working Voltage	3.7 V
Input Voltage	5 V
Operating Time	≥ 12 hours
Operating Temperature Range	36° - 104° F (2° - 40° C)
Storage Temperature Range	14° - 131° F (-10° - 55° C)

9. FCC Compliance Information

In the Ear Sound Amplifiers HA 69

Responsible Party – U.S. Contact Information:
Beurer North America LP
1 Oakwood Boulevard, Suite 255
Hollywood, FL 33020
United States
1-800-536-0366
info@beurer.com

FCC Compliance Statement

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

1. This device may not cause harmful interference, and
2. This device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications to the product not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

– Reorient or relocate the receiving antenna.
– Increase the separation between the equipment and receiver.

– Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
– Consult the dealer or an experienced radio/TV technician for help.

10. Warranty

Limited Lifetime Warranty For Original Purchaser

Your Beurer In the Ear Sound Amplifiers Amplifier, model HA 69, is warranted to be free from defects in materials and workmanship for the life of the product under normal conditions of intended use and service. This warranty extends only to the original retail purchaser and does not extend to retailers or subsequent owners.

We will, at our option, repair or replace the Beurer In the Ear Sound Amplifiers Amplifier, model HA 69, without additional charge, for any part or parts covered by these written warranties. No refunds will be given. Repair or replacement is our only responsibility and your only remedy under this written warranty. If replacement parts for defective materials are not available, Beurer reserves the right to make product substitutions in lieu of repair or replacement.

For warranty service contact our customer service department at 1-800-536-0366 or at info@beurer.com to provide a description of the problem. If the problem is deemed to be within the scope of the limited lifetime warranty, you will be asked to mail the product at your costs in its original package with proof of purchase, your name, address and phone number. If the problem is not deemed to be within the scope of the limited lifetime warranty, we will provide a quotation for repair respectively replacement and return shipping fee.

This warranty does not cover damage caused by misuse or abuse; accident; the attachment of unauthorized accessory; alteration to the product; improper installation; misapplication; lack of reasonable care with respect to the product; unauthorized repairs or modifications; improper use of electrical/power supply; old worn batteries; normal wear; loss of power; dropped product; malfunction or damage of an operating part as a result of failure to comply with instructions for use or to provide manufacturer's recommended maintenance; transit damage; theft; neglect; vandalism; or environmental conditions; loss of use during the period the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions

whatsoever that are beyond the control of Beurer. This warranty is void if the product is ever used in a commercial or business environment. The maximum liability of Beurer under this warranty is limited to the purchase price actually paid by the customer for the product covered by the warranty, as confirmed by proof of purchase, regardless of the amount of any other direct or indirect damage suffered by the customer.

This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. A product that requires modifications or adaptation to enable it to operate in any other country than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under this warranty.

THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO THE DURATION OF THIS

APPLICABLE WRITTEN WARRANTY. Some states do not allow limitations on how long an implied warranty lasts, so the above limitations may not apply to you.

IN NO EVENT SHALL BEURER BE LIABLE FOR ANY SPECIAL, INCIDENTAL, INDIRECT OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS, IMPLIED OR ANY OTHER THEORY OF LIABILITY, WHATSOEVER. Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.
Beurer does not authorize anyone, including, but not limited to, retailers, the subsequent consumer purchaser of the product from a retailer or remote purchaser, to obligate Beurer in any way beyond the terms set forth herein. This warranty does not extend to the purchase of opened, used, repaired, repackaged and/or resealed products, including but not limited to sale of such products on Internet auction sites and/or products by surplus or bulk resellers. Any and all warranties or guarantees shall immediately cease and terminate in connection with any products or parts thereof which are

repaired, replaced, altered, or modified, without the prior explicitly written consent of Beurer.

This warranty gives you specific legal rights, and you may also have other rights which may vary from state to state.
For more information regarding our product line in the USA, please visit:
www.beurer.com
Distributed by:
Beurer North America LP
1 Oakwood Boulevard, Suite 255
Hollywood, FL 33020
USA
www.beurer.com

Questions or comments? Call our US-based customer service toll free at 1-800-536-0366.

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