beurer

BM 31

EN Blood pressure monitor
Instructions for use 3
ES Tensiómetro
Manual de instrucciones 20
FR Tensiomètre
Mode d'emploi 40

Distributed by/Distribuido por/
Distribué par :
Beurer North America LP
1 Oakwood Boulevard,
Suite 255
Hollywood, FL 33020, USA
www.shop-beurer.com

For Canada only/Solo para Canadá/
Pour le Canada seulement:
Manufacturer/Fabricante/
Fabricant :
Beurer GmbH
Söflinger Str. 218
89077 Ulm Germany

Questions or Comments? Call our US-based customer service toll free at 1-800-536-0366.

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READ THIS MANUAL COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT

Keep this manual in a safe location for future reference

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ENGLISH

Contents

1. Important safety notes 4
2. Intended Use 9
3. Package Contents 9
4. Parts and Controls 9
5. Preparing for Measurement 10
6. Measuring Blood Pressure 11
7. Evaluating Results 13
8. Displaying and Deleting Measurements 14
9. Care and Maintenance 15
10. Troubleshooting Guide 15
11. Specifications 15
12. Power Supply Unit 16
13. FCC Compliance Information 17
14. Warranty 17

SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE
If the instruction manual is damaged or if you no longer have the instruction manual in your possession, please contact customer service. See warranty for service contact.

Dear Valued Customer,
Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for application in the areas of heating, weight management, blood pressure, body temperature, pulse, gentle therapy, massage, and beauty.

Register your product today to receive further benefits:
- added protection for your purchase
- health and wellness tips
- new product information
- continual outstanding customer service and support.

Visit www.registerbeurer.com or contact customer service at **1-800-536-0366**.

Please read these instructions for use carefully and keep them for later use.

Make them accessible to other users and observe the information they contain.

With kind regards,
Your Beurer team

1. Important safety notes

Signs and symbols

Whenever used, the following signs identify safety and property damage messages and designate a level of hazard or seriousness.

READ THIS ENTIRE MANUAL, THE SAFETY SECTION AND ALL INSTRUCTIONS AND WARNINGS COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT. FOLLOW ALL SAFETY INSTRUCTIONS AND WARNINGS TO AVOID HAZARDOUS SITUATIONS AND TO MAKE CORRECT USE OF THIS PRODUCT.

Definitions	
	This is the safety alert symbol. It alerts you to potential personal injury hazards. Obey all safety messages that follow this symbol to avoid possible injury or death.
WARNING	Indicates a potentially hazardous situation, which, if not avoided, could result in death or serious injury.
CAUTION	Indicates a potentially hazardous situation, which, if not avoided, may result in minor/ moderate injury.

NOTICE	Addresses practices not related to personal injury, such as product and/or property damage.
	**Product information** Note on important information
	**Observe the instructions** Read the instructions before starting work and/or operating devices or machines
	**Isolation of applied parts Type BF** Galvanically isolated applied part (F stands for “floating”); meets the requirements for leakage currents for type B
	**Direct current** The device is suitable for use with direct current only
	Permissible storage temperature and humidity
	Permissible operating temperature and humidity

IP21	**IP class** Device protected against foreign objects ≥ 12.5 mm and against water dripping at an angle
SN	Serial number

WARNING – To reduce the risk of fire, electric shock, or serious personal injury:

- For consistency, always measure your blood pressure at the same time each day.
- Before each measurement, relax for about five minutes.
- When performing multiple measurements on an individual, wait one minute between each measurement.
- Do not take a measurement within 30 minutes of eating, drinking, smoking or exercising.
- Repeat the measurement if you are unsure of the results.
- Measurements are for your information only – they are no substitute for a medical examination. Discuss the measurements with your doctor and never base any medical decisions on them (e.g. medicines and their ad- ministration).
- Using the blood pressure monitor outside your home environment or whilst on the move (e.g. whilst travelling in a car, ambulance or helicopter, or whilst undertaking physical activity such as playing sport) can influence the measurement accuracy and cause incorrect measurements.
- Do not use the blood pressure monitor on newborns or patients with pre-eclampsia. Consult your doctor before using the blood pressure monitor during pregnancy.
- Cardiovascular diseases may lead to incorrect measurements or have a detrimental effect on measurement accuracy. The same also applies to very low blood pressure, diabetes, circulatory disorders and arrhythmias as well as chills or shaking.
- This device is not intended for use by people (including children) with restricted physical, sensory or mental skills or a lack of experience and/or a lack of knowledge, unless they are supervised by a person who is responsible for their safety or are instructed by such a person in how to use the device. Supervise children around the device to ensure they do not play with it.
- The blood pressure monitor must not be used in connection with a high-frequency surgical unit.
- Only use the device on people who have the specified upper arm measurement for the device.

- Please note that when inflating, the functions of the limb in question may be impaired.
- During the blood pressure measurement, the blood circulation must not be stopped for an unnecessarily long time. If the device malfunctions, remove the cuff from the arm.
- Avoid any mechanical restriction, compression or bending of the cuff line.
- Do not allow sustained pressure in the cuff or frequent measurements. The resulting restriction of the blood flow may cause injury.
- Make sure that the cuff is not placed on an arm in which the arteries or veins are undergoing medical treatment, e.g. intravascular access or therapy, or an arteriovenous (AV) shunt.
- Do not use the cuff on people who have undergone a mastectomy.
- Do not place the cuff over wounds as this may cause further injury.
- Place the cuff on your upper arm only. Do not place the cuff on other parts of the body.
- You can either use the blood pressure monitor with batteries or with a power adapter. Please note that data transfer and data storage is only possible when your blood pressure monitor is supplied with power. As soon as the batteries are empty or the power adapter is disconnected from the power supply, the blood pressure monitor loses the date and time.
- To conserve the batteries, the blood pressure monitor switches off automatically if no buttons are pressed for one minute.
- The device is only intended for the purpose described in these instructions for use. The manufacturer is not liable for damage resulting from improper or careless use.

⚠ **CAUTION – To reduce the risk of personal injury or product/ property damage:**

- The device should not be used when your arm has been wounded/injuried or when a catheter has been inserted.Such use may result in injury.
- Remove any kind of upper arm jewelry or the like before taking a measurement. This could cause bruises.
- Do not place the arm cuff over heavy clothing (e.g. a jacket or sweater sleeve) as the blood pressure monitor will not be able to take a proper measurement and there is an elevated danger of acquiring hematoma or skin marks during the course of the measurement.
- In case the cuff does not stop inflating, interrupt the measurement by pressing the start/stop button and open the cuff at once.
- Do not disassemble the unit, it may result in injuries.

- When applying the arm cuff, make sure there are no wrinkles in the cuff. This could cause bruises.
- Do not exert any kind of pressure on the hose during measurement, e.g. laying your arms or any other object on the hose. This could cause incorrect measurements.
- Blood pressure measurements can lead to temporary marks on the skin at the site of the cuff placement. This is especially the case in high repetition rates, in hypertonic patients and in patients with weak pulses. In rare cases a mark may persist for a couple of days. Please contact your physician about these specific risks of cuff pressure in your specific case.
- For hygienic reasons the cuff is intended only for the use by one person.
- Prolonged overinflation (cuff pressure exceeds 300 mmHg or maintained above 15mmHg for longer than 3 minutes) of the bladder may cause ecchymoma of your arm.

NOTICE – To reduce the risk of product/property damage:

- The blood pressure monitor is made up of precision electronic components. Accuracy of readings and the instrument's service life depend on careful handling. Protect the unit against hard knocks (e.g. dropping the unit), moisture, water, dirt, dust, chemicals, extreme hot or cold temperatures, major temperature fluctuations, direct exposure to sunlight and heat sources which are too close (e.g. stoves, heating radiators). This may damage the unit. The device must be stored in the specified ambient conditions. Please see section "9. Care and Maintenance" and section "11. Specifications" for details.
- Do not use any aggressive solvents, cleaning agents, detergents or any other strong chemicals to clean the device.
- Never immerse and or spill water or any other liquid onto the monitor or any components, otherwise liquid will enter it and cause damage.
- Do not wash the unit in a washing machine, dish washer or tumble dryer. Observe the local regulations for material disposal. Dispose of the device, components and optional accessories according to applicable local regulations. Unlawful disposal may cause environmental pollution.
- Never attempt to repair, open and/or disassemble the unit (including arm cuff and optional accessory) or adjust it yourself. This may damage the unit and impair the functions. If you need to have the unit repaired, please contact our customer service. Please see warranty for service contact. Before submitting any complaint, first check the batteries and replace them if necessary.

- Changes or modifications to the device will nullify the user warranty.
- Do not drop or insert any object into any opening or hose. This may damage the unit.
- Do not press the buttons with excessive force or with pointed objects.

Battery Handling Safety Precautions

- Use only the size and type of batteries specified.
- Be sure to follow the correct polarity when installing the batteries. Reversed batteries may cause damage to the device.
- Do not mix different types of batteries together (e.g. Alkaline and Carbon-zinc or rechargeable batteries) or old batteries with fresh ones. Always replace batteries as a simultaneous set.
- If the batteries in the device are depleted or the device will not be used for a long period of time, remove the batteries to prevent damage or injury from possible battery leakage.
- If your skin or eyes come into contact with battery fluid, rinse the affected areas with water and seek medical assistance.
- Do not try to recharge batteries not intended to be recharged; they can overheat and rupture (follow battery manufacturer's directions.)
- Do not dispose of batteries in fire, batteries may explode or leak.
- Clean the battery contacts and also those of the device prior to battery installation.
- Remove discharged batteries from the product and dispose/recycle in compliance with all applicable laws.
- Keep batteries away from children and pets. Batteries may be harmful if swallowed. Should a child or pet swallow a battery, seek medical assistance immediately.
- Batteries can contain toxins that are harmful to the environment. Always dispose of batteries in accordance with applicable local regulations.

Notes on electromagnetic compatibility

- The device is suitable for use in all environments listed in these instructions for use, including domestic environments.
- The use of the device may be limited in the presence of electromagnetic disturbances. This could result in issues such as error messages or the failure of the display/device.
- Avoid using this device directly next to other devices or stacked on top of other devices, as this could lead to faulty operation. If, however, it is necessary to use the device in the manner stated, this device as well as

the other devices must be monitored to ensure they are working properly.

- The use of accessories other than those specified or provided by the manufacturer of this device can lead to an increase in electromagnetic emissions or a decrease in the device's electromagnetic immunity; this can result in faulty operation.
- Failure to comply with the above can impair the performance of the device.


2. Intended Use

This blood pressure monitor performs non-invasive measurement and monitoring of arterial blood pressure and heart rate on adult patients.

3. Package Contents


1 x Blood Pressure Monitor BM 31
1 x Upper Arm Cuff 8.7 in – 16.5 in (22 – 42 cm)
4 x 1.5 V LR03 AAA batteries
1 x Storage pouch
1 x Quick start guide
1 x Instruction manual

4. Parts and Controls

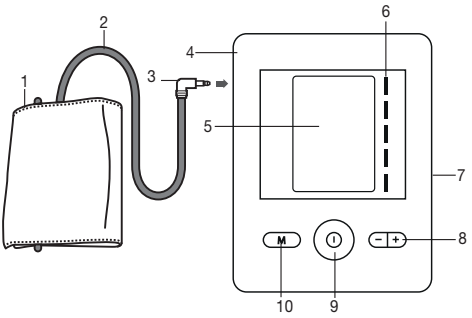1. Arm Cuff
2. Air Tube
3. Air Tube Connector
4. Air Tube Port
5. LCD Display
6. Risk Classification Scale
7. AC Adapter Jack (Adapter Not Included)
8. +/- Buttons
9. Start/Stop (⏻) Button
10. Memory (**M**) Button


LCD Display Information:

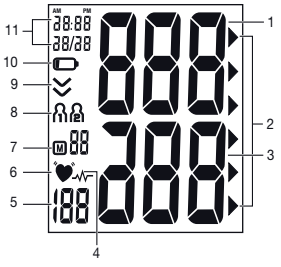1. Systolic Pressure
2. Risk Classification Indicator
3. Diastolic Pressure
4. Irregular Heartbeat Icon
5. Pulse Rate
6. Pulse Icon
7. Memory Space Number/Average (A)/ Morning Reading (AM) /Evening Reading (PM)
8. User Number
9. Deflation Icon
10. Low Battery Icon
11. Time and Date

5. Preparing for Measurement

Inserting/changing batteries:

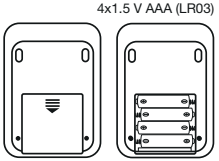1. Remove the Battery Compartment Cover from the rear of the device.
2. Insert four "AAA" alkaline batteries (included) according to the polarity markings inside the compartment.
 NOTE: Do not use rechargeable batteries.
3. Reattach the Battery Compartment Cover.

All LCD Display elements will briefly appear and the year value will flash. Set the date and time as described below. Any saved measurements will be retained.

Operating with an AC Power Adapter:

You can also operate this device with an AC Power Adapter (not included). You MUST remove the batteries before powering this unit with an AC Adapter. Only use an AC Power Adapter with specifications according to section "Power supply unit".

1. Insert the round power plug from the AC Power Adapter into the AC Adapter Jack on the right hand side of the blood pressure monitor. Plug the AC Power Adapter into a standard electrical outlet.
2. When finished using the device, unplug the AC Power Adapter from the electrical outlet and disconnect it from the blood pressure monitor. When you unplug the AC Power Adapter, date and time settings will be lost but saved measurements will remain.

Setting the hour format, date, and time:

If date and time are not set, stored measurements will not have a valid timestamp associated with them.

NOTE: Press and hold the + or - Buttons to advance values rapidly.


1. Press and hold the Memory **M** Button for about three seconds so the hour format is flashing.
2. The hour format (12h or 24h) flashes on the LCD Display. Set the hour format with the +/- Buttons and press the **M** Button to confirm.
3. The year flashes on the LCD Display. Set the year with the +/- Buttons and press the **M** Button to confirm.
4. The month flashes on the LCD Display. Set the month with the +/- Buttons and press the **M** Button to confirm.
5. The day flashes on the LCD Display. Set the day with the +/- Buttons and press the **M** Button to confirm.
6. The hour flashes on the LCD Display. Set the hour with the +/- Buttons and press the **M** Button to confirm.
7. The minutes flash on the LCD Display. Set the minutes with the +/- Buttons and press the **M** Button to confirm and exit.

6. Measuring Blood Pressure


The blood pressure monitor needs to be at room temperature when measuring. Measurements can be performed on either arm.

Positioning the Cuff


1. Place the cuff round your bare left upper arm 1–1½ in (2-3 cm) above the bend of the elbow and above the artery, with the Air Tube in line with the center of the palm. Ensure blood circulation is not restricted by tight clothing or other objects.

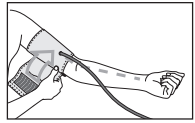2. Fasten the cuff, making sure it is not too tight around your arm. You should be able to fit two fingers under the cuff after it is in place. The cuff is suitable for you if the index mark (▼) is within the **OK** range on the cuff when it is in place.

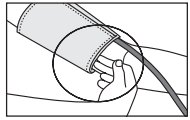
3. Insert the Air Tube Connector into the Air Tube Port. Blood pressure can vary between the right and left arm, which may mean that the measured blood pressure values are different. Always perform the measurement on the same arm. If the values between the two arms are significantly different, please consult your doctor to determine which arm should be used for the measurement.

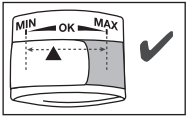CAUTION: Only use genuine Beurer cuffs with this device. The included cuff is suitable for an arm circumference of 8.7- 16.5 inches (22 to 42 cm). For arm circumferences outside this range, contact customer service.

Correct Posture

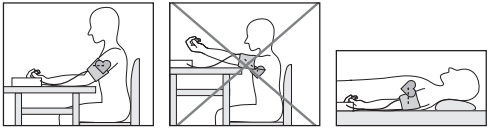- Rest for at least five minutes before each measurement.
- Perform the measurement either sitting or lying down. Always make sure that the cuff is level with your heart.
- To perform a measurement, sit comfortably with your arms and back supported. Do not cross your legs and place your feet flat on the floor.
- Keep still during the measurement and do not talk.

Performing a Measurement:

1. To select a desired user memory, press the **+/-** Buttons until the desired user memory icon appears in the LCD Display. Wait a few seconds until the selected user memory icon stops flashing and the lighted ring around the **START/STOP** Button ⏻ shuts off. Each of the two user memories can store 60 sets of measurements. If you do not select a user memory, the measurement will be stored in the most recently-used user memory.
2. With the cuff attached properly and in a correct posture, press the **START/STOP** Button ⏻. All LCD Display elements will briefly display and measurement will begin after about three seconds. To cancel a measurement at any time, press the **START/STOP** Button ⏻.
3. As a measurement concludes, air pressure in the Cuff will release. The pulse symbol (♥) will display, along with systolic and diastolic pressure.
4. E_ (where the blank space is a number 1-6) will appear if the measurement was not performed properly. Refer to the "Troubleshooting Guide" section for details about specific error message numbers and then repeat the measurement.

5. Press the **START/STOP** Button ⏻ to shut off the blood pressure monitor. The measurement is now stored in the selected user memory. The device will also shut off automatically after about one minute.
6. Wait at least five minutes before taking another measurement.

7. Evaluating Results

Irregular Heartbeat:

This device can detect irregular heartbeat patterns during measurement and will indicate them by displaying the Irregular Heartbeat Icon .
This could be an indicator for arrhythmia, a condition where the heart rhythm is abnormal. These symptoms (omitted or premature heartbeats, slow or excessively fast heart rate) may be caused, among other things, by heart disease, age, physical predisposition, excessive use of stimulants, stress, or lack of sleep. Only your doctor can determine if you do have arrhythmia.
Repeat the measurement if the Irregular Heartbeat Icon displays after a measurement. Rest for at least five minutes between measurements and do not talk or move during measurement. If the icon appears often, contact your doctor, since self-diagnosis and treatment based solely on the test results may be dangerous. It is vital to follow your doctor‘s instructions.

Risk Classification

Measurements can be classified and assessed by the table below. However, these standard values are only a general guideline, since individual blood pressure varies in different people and different age groups, etc.

Hypertension Category	Systole (in mmHg)	Diastole (in mmHg)	Action
Grade 3: Severe hypertension	≥180	≥110	Seek medical advice
Grade 2: Moderate hypertension	160–179	100–109	Seek medical advice
Grade 1: Mild hypertension	140–159	90–99	Have your doctor check your blood pressure regularly
High normal	130–139	85–89	Have your doctor check your blood pressure regularly
Normal	120–129	80–84	Continue self-monitoring
Optimal	<120	<80	Continue self-monitoring

Source: WHO, 1999 (World Health Organization)

It is important to consult your doctor regularly for advice. Your doctor will tell you your individual values for normal blood pressure as well as the value above which your blood pressure is classified as dangerous.

The classification on the display and the scale on the unit show which category the recorded blood pressure values fall into. If the systolic and diastolic values fall into two different categories, the Risk Indicator on the device always shows the higher category.

8. Displaying and Deleting Measurements

User Memory:

Every successful measurement is stored along with its date and time. If there are more than 60 measurements in a particular user memory bank, the oldest measurement is deleted.
Press the **+/-** Buttons to select a user memory and confirm by pressing the **M** Button. The most recent measurement values will display.

Average Values:

- Press the **M** Button, after about three seconds of displaying the most recent measurement, "A" flashes on the display and the average of all saved measurements in the selected user memory is displayed.

- Press the **M** Button again. "AM" flashes on the display. The average value of the morning measurements for the last seven days is displayed (morning: 5:00 AM – 9:00 AM).

- Press the **M** Button a third time. "PM" flashes on the display. The average value of the evening measurements for the last seven days is displayed (evening: 6:00 PM – 10:00 PM).

Individual Measured Values:

- Press the **M** Button again and the last individual measurement is displayed. Continue in this fashion to review older measurements.
- To shut the device off, press the Start/Stop Button.

You can exit the menu at any time by pressing the **START/ STOP** Button ⏻

Clearing Memory:

- To clear one of the user memories, you must first select it.
- Start individual measurement access as described above and select a user memory. After waiting for three seconds for the average values to display, press and hold the both the "+" and "-" Buttons at the same time for five seconds.

All values in the selected user memory will be deleted. Individual measurements cannot be deleted.

9. Care and Maintenance

- Clean the device and cuff with a slightly moistened cloth.
- Do not use detergents or solvents.
- Never immerse the device in water or other liquids.
- When storing the device, make sure no heavy objects are placed on top of it. Remove the batteries. Ensure the Air Tube does not have any sharp kinks.

10. Troubleshooting Guide

In case of errors, "E_" will appear in the LCD Display, where "_" is a number 1-6. Error messages may appear if:

- Pulse is not detected correctly and could not be measured (E1)
- No Measurement could be taken (E2)
- The cuff is fastened too tightly or loosely (E3)
- A system or device error occurs (E4)
- Air pressure is higher than 300 mmHg (E5)
- A system error occurred (E6). If this error reoccurs, contact Customer Service.
- Batteries are almost depleted .

In such cases, repeat the measurement. Do not move or speak. If necessary, replace the batteries.

11. Specifications

Model	BM 31
Measurement method	Oscillometric, non-invasive blood pressure measurement on the upper arm
Measurement range	Cuff pressure 0-300 mmHg, Systolic 50-280 mmHg, Diastolic 30-200 mmHg, Pulse 40-199 beats per minute
Display accuracy	Pressure ±3 mmHg, Pulse ±5%
Dimensions	3.9 in L x 5.5 in W x 2.1in H (98 mm L x 140 mm W x 53 mm H)
Weight	11.6 oz (329 g) with cuff, without batteries
Cuff size	8.7 - 16.5 (22 - 42 cm)
Operating Range	50° F - 104° F (10° C - 40° C), 10 - 85% relative humidity (non-condensing)
Storage Range	-4° F - 131° F (-20° C - 55° C), 10 - 90% relative humidity
Power supply	4 x 1.5V ⎓ AAA batteries
Battery life	Approx. 500 measurements

Expected Service life	500hours
Classification	Internal supply, IP21, no AP or APG, continous operation, application part type BF

Technical information is subject to change without notification to allow for updates. The serial number is located on the device or in the battery compartment.

- This device complies with European Standard EN 60601-1-2 (in accordance with CISPR 11, IEC 61000-4-2, IEC 61000-4-3, IEC 61000-4-4, IEC 61000-4-5, IEC 61000-4-6, IEC 61000-4-8, IEC 61000-4-11) and is subject to particular precautions with regard to electromagnetic compatibility. The device is suitable for use in all environments listed in these instructions for use, including domestic environments. Please note that portable and mobile HF communication systems may interfere with this device.
- The device complies with the EU Medical Devices Directive 93/42/EEC, the German Medical Devices Act (Medizinproduktgesetz) and the standards EN 1060-3 (Non-invasive sphygmomanometers – Part 3: Supplementary requirements for electro-mechanical blood pressure measuring systems) and IEC 80601-2-30 (Medical electrical equipment – Part 2 – 30: Particular requirements for the basic safety and essential performance of automated non-invasive sphygmomanometers).
- The accuracy of this blood pressure monitor has been carefully checked and developed with regard to a long useful life. If the device is used for commercial medical purposes, it must be regularly tested for accuracy by appropriate means. Precise instructions for checking accuracy may be requested from the service address.

12. Power Supply Unit

Input	100–240V, 50–60 Hz, 0.5 A max.
Output	6V DC, 600 mA, only in connection with Beurer blood pressure monitor.
Protection	This device is double insulated and protected against short circuit and overload by a primary thermal fuse Make sure to take the batteries out of the compartment before using the mains adapter.
	Polarity of the the DC voltage connection
	Double insulated/equipment class 2

Enclosures and Protective Covers	Equipment enclosed to protect against contact with live parts, and with parts which can become live (finger, pin, hook test). The operator shall not contact the patient and the output plug of AC mains part simultaneously.

13. FCC Compliance Information

Upper Arm Blood Pressure Monitor BM 31

Responsible Party – U.S. Contact Information
Beurer North America LP,
1 Oakwood Boulevard, Suite 255
Hollywood, FL 33020, United States
1-800-536-0366
info@beurer.com

FCC Compliance Statement

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

14. Warranty

Limited Lifetime Warranty For Original Purchaser

Your Beurer Upper Arm Blood Pressure Monitor, model BM 31, is warranted to be free from defects in materials

and workmanship for the life of the product under normal conditions of intended use and service. This warranty extends only to the original retail purchaser and does not extend to retailers or subsequent owners.

We will, at our option, repair or replace the Beurer Upper Arm Blood Pressure Monitor, model BM 31, without additional charge, for any part or parts covered by these written warranties. No refunds will be given. Repair or replacement is our only responsibility and your only remedy under this written warranty. If replacement parts for defective materials are not available, Beurer reserves the right to make product substitutions in lieu of repair or replacement.

For warranty service contact our customer service department at 1-800-536-0366 or at info@beurer.com to provide a description of the problem. If the problem is deemed to be within the scope of the limited lifetime warranty, you will be asked to mail the product at your costs in its original package with proof of purchase, your name, address and phone number. If the problem is not deemed to be within the scope of the limited lifetime warranty, we will provide a quotation for repair respectively replacement and return shipping fee.

This warranty does not cover damage caused by misuse or abuse; accident; the attachment of unauthorized accessory; alteration to the product; improper installation; misapplication; lack of reasonable care with respect to the product; unauthorized repairs or modifications; improper use of electrical/power supply; old worn batteries; normal wear; loss of power; dropped product; malfunction or damage of an operating part as a result of failure to comply with instructions for use or to provide manufacturer's recommended maintenance; transit damage; theft; neglect; vandalism; or environmental conditions; loss of use during the period the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions whatsoever that are beyond the control of Beurer. This warranty is void if the product is ever used in a commercial or business environment. The maximum liability of Beurer under this warranty is limited to the purchase price actually paid by the customer for the product covered by the warranty, as confirmed by proof of purchase, regardless of the amount of any other direct or indirect damage suffered by the customer.

This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. A product that requires modifications or adaptation to enable it to operate in any other country than the country for which it was designed, manufactured,

approved and/or authorized, or repair of products damaged by these modifications is not covered under this warranty. **THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO THE DURATION OF THIS APPLICABLE WRITTEN WARRANTY. Some states do not allow limitations on how long an implied warranty lasts, so the above limitations may not apply to you.**

IN NO EVENT SHALL BEURER BE LIABLE FOR ANY SPECIAL, INCIDENTAL, INDIRECT OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS, IMPLIED OR ANY OTHER THEORY OF LIABILITY, WHATSOEVER. Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.

Beurer does not authorize anyone, including, but not limited to, retailers, the subsequent consumer purchaser of the product from a retailer or remote purchaser, to obligate Beurer in any way beyond the terms set forth herein.

This warranty does not extend to the purchase of opened, used, repaired, repackaged and/or resealed products, including but not limited to sale of such products on Internet auction sites and/or products by surplus or bulk resellers. Any and all warranties or guarantees shall immediately cease and terminate in connection with any products or parts thereof which are repaired, replaced, altered, or modified, without the prior explicitly written consent of Beurer.

This warranty gives you specific legal rights, and you may also have other rights which may vary from state to state. For more information regarding our product line in the USA, please visit: www.shop-beurer.com

For Canada only:

Manufacturer:
Beurer GmbH, Söflinger Str. 218
89077 Ulm, Germany

Distributed by:
Beurer North America LP. 1 Oakwood Boulevard,
Suite 255, Hollywood, FL 33020, USA
www.shop-beurer.com

Questions or comments? Call our US-based customer service toll free at 1-800-536-0366.

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