READ THIS MANUAL COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT
Keep this manual in a safe location for future reference

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ENGLISH

SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE

If the instruction manual is damaged or if you no longer have the instruction manual in your possession, please contact Beurer customer service. Please see warranty for service contact.

Dear Valued Customer,
Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for application in the areas of heating, weight management, blood pressure, body temperature, pulse, gentle therapy, massage, and beauty.

Register your product today to receive further benefits:
- added protection for your purchase
- health and wellness tips
- new product information
- continual outstanding customer service and support.

Visit www.registerbeurer.com or contact customer service at 1-800-536-0366.

Please read these instructions for use carefully and keep them for later use.
Make them accessible to other users and observe the information they contain.

With best regards,
Your Beurer Team

IMPORTANT SAFETY NOTES

Signs and symbols
Whenever used, the following signs identify safety and property damage messages and designate a level of hazard or seriousness.

READ AND UNDERSTAND THIS ENTIRE MANUAL INCLUDING THE SAFETY SECTION AND ALL INSTRUCTIONS AND WARNINGS COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT. FOLLOW ALL SAFETY INSTRUCTIONS AND WARNINGS TO AVOID HAZARDOUS SITUATIONS AND TO MAKE CORRECT USE OF THIS PRODUCT.

⚠	This is the safety alert symbol. It alerts you to potential personal injury hazards. Obey all safety messages that follow this symbol to avoid possible injury or death.
⚠ **WARNING**	WARNING indicates a hazardous situation which, if not avoided, could result in death or serious injury.
⚠ **CAUTION**	CAUTION indicates a hazardous situation which, if not avoided, may result in minor or moderate injury.
NOTICE	NOTICE addresses practices not related to personal injury, such as product and/or property damage.

IMPORTANT SAFEGUARDS

⚠ Signs and symbols

The following symbols are used in these instructions for use, on the packaging and on the type plate for the device and accessories:

	Observe the instructions for use
	Application part, type BF
IP 64	Dust-proof and protected against splashes from all directions
LOT	Lot code
Storage / Transport	Permissible storage and transport temperature and humidity
Operating	Permissible operating temperature and humidity

READ ALL INSTRUCTIONS BEFORE USE.

⚠ **WARNING** – To reduce the risk of serious personal injury:
- Keep the product out of the reach of children and pets. This product is not intended for use by children.

⚠ **CAUTION** – To reduce the risk of personal injury or product/property damage:
- This product is intended solely for household, non-commercial, non-industrial use in measuring basal body temperature (BBT) as an aid in ovulation prediction to aid in conception **(not to be used for contraception).**
- The thermometer is intended for the measurement of human body temperature only.
- The thermometer is only designed for the measuring areas on the human body stated in these instructions.
- This thermometer features a flexible measuring tip and therefore provides more comfort and safety during measurement. Do not bend the tip more than 90°.
- Do not leave the product unattended while it is in use, especially when children are present. The product should only be used by an adult.
- Check the thermometer frequently for signs of wear and tear or damage. If detected, discontinue use and contact customer service.
- Do not walk, run, or talk while measuring a temperature.
- The thermometer contains delicate electronic components. Do not bend or drop the product or expose it to strong shocks, high temperatures or sunlight.
- Avoid exposing the product to dust, high humidity, or extreme temperatures.
- Do not attempt to disassemble the product, except to replace the battery.
- Replace the battery when the Low Battery symbol shows on the display, or inaccurate readings could result.
- Electromagnetic interference: Avoid strong electrical or electromagnetic fields in the direct vicinity of the device (e.g. mobile telephones, microwave ovens) while it is in operation, as inaccurate measurements may result. To prevent such interference, use the unit at a sufficient distance from such devices or turn the devices off.
- Influencing factors such as alcohol consumption or unusual body signs, fever, disrupted sleep, your mood or extreme or unusual physical activity may distort your basal temperature and make it unusable for evaluation. Consult your doctor if in doubt.
- Consult your doctor before use if you have any exsiting conditions that might affect use of this product.
- If you are unsure about the temperature results, consult your doctor.

NOTICE – To reduce the risk of product/property damage:
- Always observe the minimum duration of measurement until the acoustic signal.
- The thermometer conducts a self-test when switched on. There is no need to check the measurement accuracy.
- Do not allow the product to be bitten or bent.
- Clean the product after each use following the instructions provided under Care, Maintenance and Disposal.
- The thermometer is only intended to measure basal temperature. It is not intended to evaluate values or detect ovulation, or to diagnose fever or other contagious disease. Present the results to your physician or family planning counselor for interpreting results.
- Do not place the product in a dishwasher, washing machine or dryer.
- There is no guarantee that this product can assist you in getting pregnant.

⚠ Notes on electromagnetic compatibility

- The device is suitable for use in all environments listed in these instructions for use, including domestic environments.
- The use of the device may be limited in the presence of electromagnetic disturbances. This could result in issues such as error messages or the failure of the display/device.
- Avoid using this device directly next to other devices or stacked on top of other devices, as this could lead to faulty operation. If, however, it is necessary to use the device in the manner stated, this device as well as the other devices must be monitored to ensure they are working properly.
- The use of accessories other than those specified or provided by the manufacturer of this device can lead to an increase in electromagnetic emissions or a decrease in the device's electromagnetic immunity; this can result in faulty operation.
- Failure to comply with the above can impair the performance of the device.

Battery Handling Safety Precautions

- Use only the size and type of batteries specified.
- Be sure to follow the correct polarity when installing the batteries. Reversed batteries may cause damage to the device.
- If the batteries in the device are depleted or the device will not be used for a long period of time, remove the batteries to prevent damage or injury from possible battery leakage.
- Do not try to recharge batteries not intended to be recharged; they can overheat and rupture (follow battery manufacturer's directions).
- Do not dispose of batteries in fire, batteries may explode or leak.
- Clean the battery contacts and also those of the device prior to battery installation.
- Remove discharged batteries from the product and dispose/recycle in compliance with all applicable laws.
- Keep batteries away from children and pets. Batteries may be harmful if swallowed. Should a child or pet swallow a battery, seek medical assistance immediately
- If your skin or eyes come into contact with battery fluid, rinse the affected areas with water and seek medical assistance.


Intended Use

The basal thermometer is a special digital thermometer designed to measure the basal temperature in women of childbearing age as an aid in ovulation prediction to aid in conception (not to be used for contraception). The thermometer has 30 memory spaces for measurements. The thermometer is only intended to measure basal temperature. It is not intended to evaluate values or detect ovulation, or to diagnose fever or other contagious disease.
Transferring the measurements to the Ovy app enables fertile and non-fertile days in the cycle to be predicted, as well as periods. It can be used for oral, rectal or vaginal measurement.

Package Contents

1 x Basal thermometer OT 30
1 x 3V battery CR 2032 (included)
1 x Protective cover
1 x Instruction manual

Parts and Controls

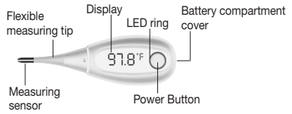To Use

NOTE: A Basal Thermometer is used to chart basal body temperature (BBT), the basis of a fertility chart. Basal temperature is the lowest normal body temperature after waking in the morning (before getting up, talking, moving around, eating, drinking, or smoking). Take your temperature at the same time every day, as soon as you wake up and before you get out of bed. You should have slept for at least four hours before measurement. If you have slept for less than four hours, enter the amount of sleep you have had in the Ovy app or select "Disrupted sleep." Do not eat anything before the measurement and avoid physical exertion.
Clean the thermometer before and after each use with a soft cloth and isopropyl alcohol diluted with water. Wipe with a clean, damp cloth and let dry. **NOTICE:** Do not soak the product or sterilize by boiling or steaming it.

1. A battery is included with your thermometer and is already in place; you will need to remove the battery isolation strip from the battery compartment cover. To do so, unscrew the battery compartment cover using a small coin. Pull out the battery isolation strip and reattach the cover.
2. Remove the protective cap from the measuring tip.
3. Press the On/Off Button; the thermometer will beep. The thermometer will carry out a self-test for about 2 seconds. During the selftest all segments of the LCD Display will activate briefly. 
4. The product will then briefly display the reference reading of 97.70 °F (36.5 °C) and then "Lo" with a flashing "°F", which indicates the thermometer is ready for use. **NOTICE:** During measurement, the display shuts off to save the battery and the LED ring flashes slowly. Measuring ends when a stable temperature is reached. A beep will sound and the measured temperature displays.
5. Promptly begin temperature measuring. Check section below for different measuring areas. Note that it may take up to three minutes for the temperature to stabilize and record a reliable measurement, therefore a measuring time of 3 minutes is recommended. Please always wait until the end of the measurement before you remove the thermometer from the place of measurement.
6. Press the ON/OFF Button to turn off the thermometer. **NOTE:** The thermometer switches off automatically around 30 seconds after measurement.

Tips for correct use:
- Influencing factors, such as alcohol consumption or unusual body signs, fever, disrupted sleep, your mood or extreme or unusual physical activity may distort your basal temperature and make it unusable for evaluation.
- Do not eat and avoid physical exertion before measurement.
- You should have slept at least 4 hours before the measurement. An unusual lack of sleep may distort your basal temperature. In that case, enter the amount of sleep in the Ovy app or your basal temperature chart or mark the measurement as "disrupted sleep".
- Take your temperature as soon as you wake up and before you get out of bed.
- The most accurate reading is taken at the same time every morning. Your resting (basal) temperature will measure low and stable during the first weeks of your monthly cycle. Recording your basal temperature on a regular basis will document a temperature fluctuation pattern in your monthly cycle.
- Do not change the thermometer or measurement area during your cycle. If you change the thermometer or the measurement area, mark it as a disruption. Reliable evaluation is no longer guaranteed for the current cycle.
- The thermometer does not detect ovulation or evaluate any measured values.

Measuring Areas

You can take oral, rectal, or vaginal measurements, but always use the same place once you begin. We recommend taking the measurement under your tongue with your mouth closed.

⚠ **CAUTION** Do not talk, eat, drink, bathe, exercise, or smoke before a temperature measurement.

Oral Temperature Method
- Place the Probe Tip under the tongue as near as possible to a heat pocket, to the left or right of the root of your tongue, and keep the mouth closed (see Fig. 1). The Probe Tip must have good contact with tissue.
- Close your mouth and breathe calmly through your nose during measurement.

Fig. 1

Rectal Temperature Method
- Apply a water-soluble lubricant to the Probe Tip. Do not use petroleum jelly.
- Gently insert the Probe Tip into the rectum approximately ½ inch (1 cm). Normal rectal temperature is 97.90 °F (36.6 °C).

Vaginal Temperature Method
- Gently insert the Probe Tip into the vagina approximately ½ inch (1 cm).

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NOTE: Do not use a thermometer that has been used for rectal temperature measuring for, vaginal or oral temperature measuring.

Memory

The thermometer has 30 memory spaces for measurements. The memory menu cannot be opened during data transfer.

1. To access the memory menu, press the Power Button to switch the thermometer on. After the self-test, briefly press the Power Button again. The thermometer now sends measurements to your smartphone. During data transfer, the LED ring on the thermometer lights up blue. Once the data transfer has been successfully completed, the LED ring flashes twice.
2. The measurement is displayed during data transfer. After data transfer, the memory space number first appears in the display and then the relevant measurement.
3. To scroll among saved measurements, press the Power Button repeatedly. To exit memory mode and shut the thermometer off, press and hold the Power Button for three seconds.

Connecting the thermometer with a smartphone

The thermometer can communicate with the Ovy app to document and keep track of your cycles.
Before using the thermometer for the first time, you will need to pair it with the Ovy mobile app in order to synchronize the time. For more details about the app, refer to the "Ovy app" section of this manual. Before you can transfer measurements to your smartphone, the thermometer must be paired with your phone as follows:

1. Download the Ovy app from the Apple App Store (iOS) or from Google Play (Android™).
2. Activate Bluetooth® on your smartphone and then open the Ovy app.
3. Press the Power Button on the thermometer to power it on.
4. Press the Power Button again to access stored measurements.
5. Follow the instructions in the Ovy app. Once the Ovy app pairs with the thermometer, a six-number code appears in the display of the thermometer.
6. Enter the code in your smartphone to complete the pairing process.

NOTE: You can only set the measurement units (°F/°C) and turn beeps on or off via the Ovy app.

The Ovy App

Ovy is designed for use specifically with Beurer basal thermometers by woman who want to conceive. It is NOT to be used for contraception. The app is available for download in the Apple App Store and Google Play.

How does the app work?
The Ovy app works out your cycle automatically, without pen and paper. This includes your period, ovulation and your fertile and non-fertile windows. To do this, take your temperature when you wake up in the morning and enter the value (with two decimal places) in the Ovy app. To make the prediction even more accurate, add other body signals each day. You should ideally start recording this on the first day of your period; i.e. the first day of your menstrual cycle. However, you can also start recording on any day in your cycle. It is important that you record when potential external temperature influences such as a fever, consumption of alcohol or a lack of sleep have occurred. You can also enter your basal temperature, body signals and influencing factors at a later time, but the temperature measurement should be taken on the relevant day. The Ovy app is based on a self-learning algorithm. This means that the more often you enter your basal temperature, the more accurately Ovy will be able to predict your fertile days and your next period.

Do you have any Questions about the Ovy app?

The Ovy team is on hand to help! You are welcome to e-mail us at hello@ovyapp.com. The Ovy app also has a contact function in your profile, which you can also use to send questions and feedback directly. The Ovy website and the information sections in the app already offer some answers and diverse background information about your cycle.

Website: www.ovyapp.com | Email: hello@ovyapp.com
Facebook: Ovy | Instagram: @ovyapp | Twitter: @ovyapp

Care, Maintenance and Disposal

Clean the basal thermometer with a soft, clean cloth before and after each use.

· For stubborn stains, wipe with a damp cloth and a mild detergent. Then dry the basal thermometer with a soft, dry cloth.
· You can also disinfect the basal thermometer with 75% ethanol or isopropyl alcohol diluted with water. When disinfecting, do NOT immerse the measurement sensor (the thermometer's tip) in alcohol or hot water (≥ 122° F/50° C).
· To prevent product damage, do not use benzene, gasoline or other strong thinners or solvents.
· Do NOT boil the thermometer to sterilize.

To store, replace the protective cap and place the thermometer in a cool, dry place away from direct sunlight and temperature extremes.

Disposal
Observe the local regulations for material disposal. Dispose of the device in accordance with local regulations. If you have any questions, please contact the local authorities responsible for waste disposal.

Battery Replacement

⚠ **CAUTION**
Do not replace the battery in the presence of small children or pets, or during conditions of extreme heat.

1. Replace the battery when the low battery indicator (▭) appears in the LCD Display.
2. Unscrew the battery compartment cover using a coin.
3. Remove the button cell battery from the battery compartment and replace with a new CR2032 battery with the + symbol facing upwards.
4. Reattach battery compartment cover.

Troubleshooting

Lo °F	The measured temperature is below 89.60 °F (32 °C) and therefore outside the measurement range.
Hi °F	The measured temperature is above 109.39 °F (43 °C) and therefore outside the measurement range.
Err	Electronic fault in the device. If this occurs repeatedly, please contact customer service.
▭	The battery is flat. Replace the battery.

Technical Specifications

Type	OT 30
Measurement range	89.60 - 111.18 °F (32.00 - 43.99 °C)
Measurement accuracy	± 0.09 °F (0.05 °C) between 95.00 – 100.00 °F (35 - 38 °C) ±0.18 °F/0.1 °C outside of this measurement and temperature range
Display	6-digit OLED
Storage and transportation conditions	Temperature: -4 -131 °F (-20 -55 °C) Relative humidity: 15 % to 95 % Ambient pressure: 700 to 1060 hPa
Operating conditions	Temperature: 41-104 °F (5-40 °C) Relative humidity: 15 % to 95 % 700 to 1060 hPa ambient pressure
Power	3V CR2032 battery (included)
Data transfer:	2402 - 2480 MHz Max. transmission power: 2.8 dBM
System requirements:	iOS ≥ 10.0, AndroidTM ≥ 5.0, Bluetooth® ≥ 4.0

List of compatible devices:

Technical information is subject to change without notification.
The serial number can be found on the unit itself or in the battery compartment.

This device meets the requirements of the EU Directive 93/42/EEC concerning medical devices, as well as those of the Medizinproduktegesetz (German Medical Devices Act), the European Standard EN 12470-3: Clinical thermometers – Part 3: Performance of compact electrical thermometers (non-predictive and predictive) with maximum device and the European Standard EN 60601-1-2 (In accordance with CISPR 11, IEC 61000-4-2, IEC 61000-4-3 and IEC 61000-4-8), and is subject to particular precautions with regard to electromagnetic compatibility. Please note that portable and mobile HF communication systems may interfere with this device.

FCC Compliance Information

Basal Thermometer OT 30
Responsible Party – U.S. Contact Information:
Beurer North America LP
1 Oakwood Boulevard, Suite 255
Hollywood, FL 33020
United States
1-800-536-0366
info@beurer.com

FCC Compliance Statement
This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:
1. This device may not cause harmful interference, and
2. This device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interferencein a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. Ifthis equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:
— Reorient or relocate the receiving antenna.
— Increase the separation between the equipment and receiver.
— Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
— Consult the dealer or an experienced radio/TV technician for help.

Warranty

Limited Lifetime Warranty For Original Purchaser
Your Beurer Basal Thermometer, model OT 30, is warranted to be free from defects in materials and workmanship for the life of the product under normal conditions of intended use and service. This warranty extends only to the original retail purchaser and does not extend to retailers or subsequent owners. We will, at our option, repair or replace the Beurer Basal Thermometer, model OT 20, without additional charge, for any part or parts covered by these written warranties. No refunds will be given. Repair or replacement is our only responsibility and your only remedy under this written warranty. If replacement parts for defective materials are not available, Beurer reserves the right to make product substitutions in lieu of repair or replacement.
For warranty service contact our customer service department at 1-800-536-0366 or at info@beurer.com to provide a description of the problem. If the problem is deemed to be within the scope of the limited lifetime warranty, you will be asked to mail the product at your costs in its original package with proof of purchase, your name, address and phone number. If the problem is not deemed to be within the scope of the limited lifetime warranty, we will provide a quotation for repair respectively replacement and return shipping fee.
This warranty does not cover damage caused by misuse or abuse; accident; the attachment of unauthorized accessory; alteration to the product; improper installation; misapplication; lack of reasonable care with respect to the product; unauthorized repairs or modifications; improper use of electrical/power supply; old worn batteries; normal wear; loss of power; dropped product; malfunction or damage of an operating part as a result of failure to comply with instructions for use or to provide manufacturer's recommended maintenance; transit damage; theft; neglect; vandalism; or environmental conditions; loss of use during the period the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions whatsoever that are beyond the control of Beurer. This warranty is void if the product is ever used in a commercial or business environment. The maximum liability of Beurer under this warranty is limited to the purchase price actually paid by the customer for the product covered by the warranty, as confirmed by proof of purchase, regardless of the amount of any other direct or indirect damage suffered by the customer.
This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. A product that requires modifications or adaptation to enable it to operate in any other country than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under this warranty.

THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO THE DURATION OF THIS APPLICABLE WRITTEN WARRANTY. Some states do not allow limitations on how long an implied warranty lasts, so the above limitations may not apply to you.

IN NO EVENT SHALL BEURER BE LIABLE FOR ANY SPECIAL, INCIDENTAL, INDIRECT OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS, IMPLIED OR ANY OTHER THEORY OF LIABILITY, WHATSOEVER. Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you. Beurer does not authorize anyone, including, but not limited to, retailers, the subsequent consumer purchaser of the product from a retailer or remote purchaser, to obligate Beurer in any way beyond the terms set forth herein.
This warranty does not extend to the purchase of opened, used, repaired, repackaged and/or resealed products, including but not limited to sale of such products on Internet auction sites and/or products by surplus or bulk resellers. Any and all warranties or guarantees shall immediately cease and terminate in connection with any products or parts thereof which are repaired, replaced, altered, or modified, without the prior explicitly written consent of Beurer.

This warranty gives you specific legal rights, and you may also have other rights which may vary from state to state. For more information regarding our product line in the USA, please visit: www.beurer.com

Distributed by:
Beurer North America LP
1 Oakwood Blvd., Suite 255
Hollywood, FL 33020, USA

Made in China

Questions or comments? Call our US-based customer service toll free at 1-800-536-0366.

Subject to error and change

Apple and the Apple logo are trademarks of Apple Inc., registered in the U.S. and other countries.
App Store is a service mark of Apple Inc., registered in the U.S. and other countries.
Google Play and the Google Play logo are trademarks of Google LLC. Android is a trademark of Google LLC.

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