beurer

FT 65

ENGLISH Multi-Function Thermometer
Instructions for use2

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READ THIS MANUAL COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT
Keep this manual in a safe location for future reference

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Questions or Comments?
Call our US-based customer service toll free at 1-800-536-0366.

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ENGLISH

Contents

1. Important safety notes 3
2. Getting to Know your Instrument 6
3. Package Contents 8
4. Parts and Controls 8
5. To Set Up 9
6. To Use 9
7. Care and Maintenance 11
8. Troubleshooting Guide 12
9. Specifications 12
10. FCC Compliance Information 13
11. Warranty 14

SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE

If the instruction manual is damaged or if you no longer have the instruction manual in your possession, please contact customer service. See warranty for service contact.

Dear Valued Customer,

Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for application in the areas of heating, weight management, blood pressure, body temperature, pulse, gentle therapy, massage, and beauty.

Register your product today to receive further benefits:

- added protection for your purchase
- health and wellness tips
- new product information
- continual outstanding customer service and support.


Visit www.registerbeurer.com or contact customer service at 1-800-536-0366.
Please read these instructions for use carefully and keep them for later use.
Make them accessible to other users and observe the information they contain.


With kind regards, Your Beurer Team

1. Important safety notes

Signs and symbols
Whenever used, the following signs identify safety and property damage messages and designate a level of hazard or seriousness.

READ THIS ENTIRE MANUAL INCLUDING THE SAFETY SECTION AND ALL INSTRUCTIONS AND WARNINGS COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT. FOLLOW ALL SAFETY INSTRUCTIONS AND WARNINGS TO AVOID HAZARDOUS SITUATIONS AND TO MAKE CORRECT USE OF THIS PRODUCT.

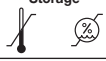	This is the safety alert symbol. It alerts you to potential personal injury hazards. Obey all safety messages that follow this symbol to avoid possible injury or death.
 **CAUTION**	CAUTION – Indicates a hazardous situation, which, if not avoided, may result in minor/moderate injury.
NOTICE	NOTICE – Addresses practices not related to personal injury, such as product and/or property damage.
	**Note** Note on important information.
	Observe the instructions for use
	Application part, type BF
Storage 	Permissible storage temperature and humidity

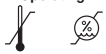Operating 	Permissible operating temperature and humidity
SN	Serial number

CAUTION – To reduce the risk of personal injury or product/property damage:

- This device is designed solely for non-commercial, household use in non-invasively measuring temperatures via an infrared sensor of the auditory canals of adults and children, as well as detecting surface temperatures of objects and ambient temperatures of surrounding areas. Do not use for any other purpose.
- Do not disassemble this unit.
- This thermometer is a sensitive electronic device. Treat it carefully and do not expose it to any sharp impacts.
- Do not expose the thermometer to direct sunlight.
- Place thermometer in the room where measurements will be made for at least 30 minutes before use.
- The thermometer is NOT waterproof; keep away from water and other liquids.
- Clean the measuring tip after each use with a soft cloth moistened with disinfectant.
- Before each use, ensure the sensor is undamaged; otherwise, contact customer service.
- This device is not intended for commercial or clinical use.
- Different measurement results are possible between the left and right ears. In order to record temperature changes, always measure in the same ear for the same person.
- If you have been lying on one ear for some time, if your ears have been covered or exposed to heat or cold, you have just been swimming or bathing, measurements will be inaccurate; wait about 20 minutes before measuring.
- Do not use on children under six months of age. In infants under six months, the ear passage is still very narrow and the temperature of the eardrum can often not be recorded, so any results obtained may be inaccurate.

- The measurement must not be taken in an ear affected by inflammatory diseases (e.g. discharging pus or secretion), after possible ear injuries (e.g. eardrum damage) or after any ear operation while the ear is healing.
- Do not use the thermometer on multiple people if any of them have an ear infection.
- Read these instructions for use carefully and keep them for later use.

NOTICE:

- Taking your temperature provides information about your current state of well-being. If you are uncertain about interpreting the results or if the result is abnormal (e.g. fever), consult your doctor. This also applies in the case of slight temperature changes if there are other symptoms of illness such as agitation, severe sweating, flushed skin, fast pulse rate, tendency to collapse, etc.
- Temperatures measured with different thermometers cannot be compared. You should tell your doctor or keep in mind if you are measuring for yourself what type of thermometer you use and in which part of the body.

Notes on electromagnetic compatibility

- The device is suitable for use in all environments listed in these instructions for use, including domestic environments.
- The use of the device may be limited in the presence of electromagnetic disturbances. This could result in issues such as error messages or the failure of the display/device.
- Avoid using this device directly next to other devices or stacked on top of other devices, as this could lead to faulty operation. If, however, it is necessary to use the device in the manner stated, this device as well as the other devices must be monitored to ensure they are working properly.
- The use of accessories other than those specified or provided by the manufacturer of this device can lead to an increase in electromagnetic emissions or a decrease in the device's electromagnetic immunity; this can result in faulty operation.
- Failure to comply with the above can impair the performance of the device.

Battery Handling Safety Precautions

- Use only the size and type of batteries specified.
- Be sure to follow the correct polarity when installing the batteries. Reversed batteries may cause damage to the device.
- Do not mix different types of batteries together (e.g., alkaline and carbon-zinc or rechargeable batteries) or old batteries with fresh ones. Always replace batteries as a simultaneous set.
- If the batteries in the device are depleted or the device will not be used for a long period of time, remove the batteries to prevent damage or injury from possible battery leakage.
- Do not try to recharge batteries not intended to be recharged; they can overheat and rupture (follow battery manufacturer's directions.)
- Do not dispose of batteries in fire, batteries may explode or leak.
- Clean the battery contacts and also those of the device prior to battery installation.
- Remove discharged batteries from the product and dispose/recycle in compliance with all applicable laws.
- Should battery fluid leak and come into contact with your eyes or skin, immediately rinse with plenty of clean water. Please call 911 immediately.
- Batteries should not be taken apart, thrown in the fire or short circuited.
- If a battery has leaked, put on protective gloves and clean the battery compartment with a dry cloth.
- Keep batteries away from children and pets. Batteries may be harmful if swallowed. Should a child or pet swallow a battery, seek medical assistance immediately.

2. Getting to Know your Instrument

This infrared thermometer measures:

- Temperature in the human ear
- Temperature of the human forehead
- Surface temperature of objects

Additional features include:

- Stores last 10 measurements for plotting changes in temperature
- Stores date and time with each measurement

- Fever alarm for temperatures over 99.5 °F (37.5° C)
- Switchable between °F and °C

Different types of thermometers measure body temperature differently. The resulting temperature varies depending on the part of the body where the measurement is taken. Temperatures in various parts of the body can vary up to two degrees Fahrenheit.

Measurement Location	**Normal Temperature Range**
Forehead temperature	96.4 °F to 99.7 °F (35.8 °C to 37.6 °C)
Ear temperature	96.8 °F to 100.0 °F (36.0 °C to 37.8 °C)
Oral temperature	96.8 °F to 99.3 °F (36.0 °C to 37.4 °C)

NOTE: Temperatures measured with different thermometers cannot be compared. You should tell your doctor or keep in mind if you are measuring for yourself what type of thermometer you use and in which part of the body.

Temperature is also influenced by:
- Individual metabolism
- Age (body temperature is higher in babies and small children and falls with increasing age. Greater temperature fluctuations occur faster and more often in children)
- Clothing
- Outside temperature
- The time of day (body temperature is lower in the morning and increases slightly throughout the day).
- Physical and, to a lesser extent, mental activity.


NOTE: Taking your temperature provides information about your current state of well-being. If you are uncertain about interpreting the results or if the result is abnormal (e.g. fever), consult your doctor. This also applies in the case of slight temperature changes if there are other symptoms of illness such as agitation, severe sweating, flushed skin, fast pulse rate, tendency to collapse, etc.


3. Package Contents

- 1 x Multi-Function thermometer
- 2 x 1.5V AAA batteries (LR03)
- 1 x Instruction manual

4. Parts and Controls


Overview

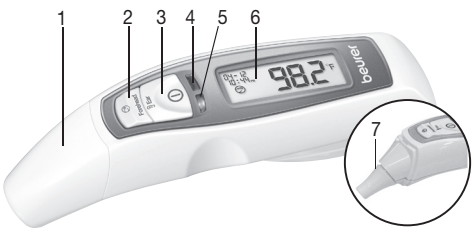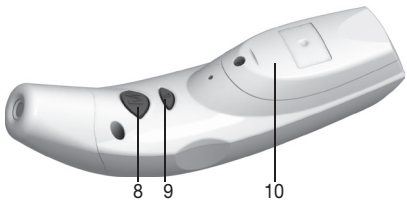1. Sensor Cover
2. Forehead Thermometer Mode Button
3. Ear thermometer mode /Power Button
4. Red LED - High Temperature Indicator
5. Green LED - Normal Temperature Range Indicator
6. LCD Display
7. Temperature Sensor (Cover Removed)
8. Memory Button (M)
9. SET Button (S)
10. Battery Compartment Cover

11. Date
12. Time
13. Measurement mode (Ear/Forehead/Object)
14. Low Battery Indicator
15. Temperature Measurement
16. Memory storage space number
17. Measurement in Progress Indicator
18. Unit of measurement (°F or °C)

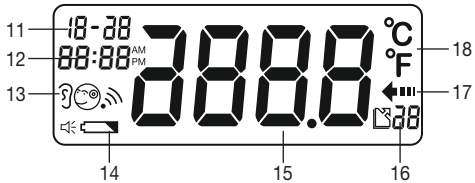5. To Set Up

Before first use, pull out the clear battery isolation strip to activate the batteries.


Setting the time and date:

1. With the device switched on, press and hold the SET Button for about three seconds; the Celsius/Fahrenheit scale indicator will appear on the LCD Display. Press the SET Button again to select Celsius (°C) or Fahrenheit (°F). Press the M button to save and continue.
2. The time format option flashes on the display.Use the SET button to select 12 or 24-hour time format and confirm with the M button. Press the M button to save and continue. The year will appear on the LCD Display. Press the SET button to set the year and then press the M Button Ⓜ to save and continue.
3. Continue in this fashion to set month, day, hours, and minutes, in turn.
4. Press the M button Ⓜ to save and exit; the device will shut off.


6. To Use

NOTICE: Before each use, make sure the thermometer, and especially the sensor, are undamaged. Note that the thermometer should be in the same room where it is to be used for at least 30 minutes in order to acclimatize it.

Ear Temperature Measurement

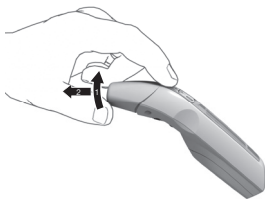1. Press and hold the Power Button ⏻ for one second to switch the thermometer on. Following a successful self test, the thermometer will beep twice, indicating it is ready for use. Remove the Sensor Cover by gently pushing it up (1) and then pulling it off (2).
2. Make sure the sensor tip and ear canal are clean. Since the ear canal is slightly curved, you will need to pull the outer ear slightly up and backwards as shown before measurement. This enables the sensor tip to be pointed directly at the eardrum for a more accurate measurement.
3. Insert the sensor tip carefully; press and hold the Power/Ear Button for one second, then release. The thermometer will beep when the measurement concludes and will display the ear temperature.
4. If the measured value is below 100.4 °F/38 °C, the green Normal Temperature Range LED will illuminate for about three seconds. If the measured value is equal to or higher than 100.4 °F/38 °C, the red High Temperature Indicator LED will illuminate.

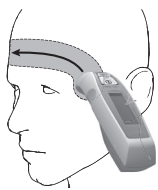Forehead Temperature Measurement

1. Press the Power Button ⏻ to switch the thermometer on. With the Sensor Cover in place, position the sensor against one of the temples on the person to be measured.
2. Press and hold the Forehead Button. While holding the Forehead Button down, move the thermometer across the person's forehead to the opposite temple, then release.
3. The thermometer will beep twice and display the measurement.
4. If the measured value is below 100.4 °F/38 °C, the green Normal Temperature Range LED will illuminate for about three seconds. If the measured value is equal to or higher than 100.4 °F/38 °C, the red High Temperature Indicator LED will illuminate.

Surface temperatures:

1. Switch the thermometer on and then press and hold the M Ⓜ and Power ⏻ Buttons for about three seconds; the symbol will appear in the LCD Display indicating the thermometer is in Object Temperature Mode.
2. Hold the sensor 1-1 ½ inches (2-3 cm) from the object to be measured and press and hold the Power Button ⏻. After about a second, release the Power Button ⏻; the thermometer will beep once and display the measurement.
3. Shut the thermometer off to exit Object Temperature Mode.

Viewing Stored Measurements

This device stores 10 measurements, along with the date/time and measurement mode for each. With the device switched on, press the M Button Ⓜ to display the last measurement. Press the M Button Ⓜ repeatedly to cycle through the older measurements.

7. Care and Maintenance

- The sensor is the most sensitive part of the thermometer. Work carefully when cleaning the sensor to avoid product damage. Use a clean cloth or cotton swab moistened with disinfectant, alcohol, or warm water.
- To clean the outer surface of the thermometer, use a soft cloth slightly moistened with a mild soapy solution.
- Do not use any harsh cleaning products.
- Before prolonged storage, remove the batteries.
- Do not use or store the device in excessively high or low temperature or humidity, in direct sunlight, near an electrical current or in dusty locations.
- Always store the thermometer with the sensor cover on.
- Never immerse the device in water.
- To protect the sensor tip, always replace sensor cover after use.
- To clean the device, use a soft, damp cloth and a mild detergent. Ensure water does not get into the device. If the device does get wet, remove the battery immediately and let dry completely before reuse.

Battery Replacement
Remove the battery compartment cover and replace the depleted batteries with two 1.5V AAA (LR03) batteries, observing the correct polarity as indicated inside the compartment. Replace the battery compartment cover.

8. Troubleshooting Guide

Error message	Problem	Solution
Err	Operating temperature is outside the acceptable range of 59 °F – 95 °F (15 °C – 35 °C)	Only use thermometer within measurement range.
HI	Ear/Forehead temperature is above 109.4 °F (43 °C) or surface temperature above 212 °F (100 °C).	Only measure temperatures within the thermometer‘s range.
Lo	Ear/Forehead temperature below 93.2 °F (34 °C) or surface temperature below 32 °F (0 °C).	Only measure temperatures within the thermometer‘s range.

9. Specifications

Temperature Measurement Range	Forehead: 93.2 °F – 109.4 °F (34.0 °C – 43 °C) Object: 32 °F – 212 °F (0 °C – 100 °C)
Measurement Accuracy	Forehead/Ear temperature: ±0.4 °F (0.2 °C) from 95 °F – 107.6 °F (35 °C – 42 °C); Object temperature: ±2.7 °F (1.5 °C)
Operating conditions	59 °F – 95 °F (15 °C – 35 °C), Humidity < 95% (noncondensing)
Storage conditions	-13 °F – 131 °F (-25 °C – 55 °C), Humidity < 85% (noncondensing)

Memory Slots	10
Display resolution	0.1 °
Power	2 x 1.5V AAA (LR03)
Dimensions	6.3 x 1.7 x 1.4 in (161 x 43 x 36 mm)
Weight	Approx. 3.2 oz (90 g) with batteries
Automatic shutoff	60 seconds

These specifications are subject to change without notice for purpose of improvement.
The serial number can be found on the unit itself or in the battery compartment.
This device complies with the EU Directive 93/42/EEC concerning medical products, the Medizinproduktegesetz (German Medical Devices Act), the ASTM (American Society for Testing and Materials) E 1965-98 and the European Standard EN 12470-5: Clinical thermometers – Part 5: Requirements for infrared ear thermometers (with maximum device) and the European Standard EN 60601-1-2 (in accordance with CISPR 11, IEC 61000-4-2, IEC 61000-4-3, IEC 61000-4-8), and is subject to particular precautions with regard to electromagnetic compatibility.
Please note that portable and mobile HF communication systems may interfere with this device.

10. FCC Compliance Information

Multi-Function Thermometer FT 65

Responsible Party U.S. Contact Information:
Beurer North America LP
1 Oakwood Boulevard, Suite 255
Hollywood, FL 33020
United States
1-800-536-0366
info@beurer.com

FCC Compliance Statement
This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:
(1) this device may not cause harmful interference, and
(2) this device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.
NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:
— Reorient or relocate the receiving antenna.
— Increase the separation between the equipment and receiver.
— Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
— Consult the dealer or an experienced radio/TV technician for help.

11. Warranty

Limited Lifetime Warranty For Original Purchaser
Your Beurer Multi-Function Thermometer, model FT 65, is warranted to be free from defects in materials and workmanship for the life of the product under normal conditions of intended use and service. This warranty extends only to the original retail purchaser and does not extend to retailers or subsequent owners.
We will, at our option, repair or replace the Beurer Multi-Function Thermometer, model FT 65, without additional charge, for any part or parts covered by these written warranties. No refunds will be given. Repair or replacement is our only responsibility and your only remedy under this written warranty. If replacement parts for defective materials are not available, Beurer reserves the right to make product substitutions in lieu of repair or replacement.

For warranty service contact our customer service department at 1-800-536-0366 or at info@beurer.com to provide a description of the problem. If the problem is deemed to be within the scope of the limited lifetime warranty, you will be asked to mail the product at your costs in its original package with proof of purchase, your name, address and phone number.
If the problem is not deemed to be within the scope of the limited lifetime warranty, we will provide a quotation for repair respectively replacement and return shipping fee.
This warranty does not cover damage caused by misuse or abuse; accident; the attachment of unauthorized accessory; alteration to the product; improper installation; misapplication; lack of reasonable care with respect to the product; unauthorized repairs or modifications; improper use of electrical/power supply; old worn batteries; normal wear; loss of power; dropped product; malfunction or damage of an operating part as a result of failure to comply with instructions for use or to provide manufacturer's recommended maintenance; transit damage; theft; neglect; vandalism; or environmental conditions; loss of use during the period the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions whatsoever that are beyond the control of Beurer. This warranty is void if the product is ever used in a commercial or business environment. The maximum liability of Beurer under this warranty is limited to the purchase price actually paid by the customer for the product covered by the warranty, as confirmed by proof of purchase, regardless of the amount of any other direct or indirect damage suffered by the customer.
This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. A product that requires modifications or adaptation to enable it to operate in any other country than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under this warranty.
THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO THE DURATION OF THIS APPLICABLE WRITTEN WARRANTY. Some states do not allow limitations on how long an implied warranty lasts, so the above limitations may not apply to you.
IN NO EVENT SHALL BEURER BE LIABLE FOR ANY SPECIAL, INCIDENTAL, INDIRECT OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS, IMPLIED OR ANY OTHER THEORY OF LIABILITY, WHATSOEVER. Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.

Beurer does not authorize anyone, including, but not limited to, retailers, the subsequent consumer purchaser of the product from a retailer or remote purchaser, to obligate Beurer in any way beyond the terms set forth herein.
This warranty does not extend to the purchase of opened, used, repaired, repackaged and/or resealed products, including but not limited to sale of such products on Internet auction sites and/or products by surplus or bulk resellers. Any and all warranties or guarantees shall immediately cease and terminate in connection with any products or parts thereof which are repaired, replaced, altered, or modified, without the prior explicitly written consent of Beurer.
This warranty gives you specific legal rights, and you may also have other rights which may vary from state to state.
For more information regarding our product line in the USA, please visit: www.beurer.com

Distributed by:
Beurer North America LP
1 Oakwood Boulevard, Suite 255
Hollywood, FL 33020, USA
www.beurer.com

Questions or comments? Call our US-based customer service toll free at 1-800-536-0366.

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Contracted Manufacturer: AViTA (WuJiang) Co. LTD